



Cadenza **CX-C18**  
 Unison **UK-C8**  
 Cadenza **CW-C18**  
 Unison **UK-C18**  
 Cadenza **CD-C18**

75 x 4.6 mm

Application

**Ofloxacin (fluoroquinolone)**

オフロキサシン(キノロン系抗菌薬)

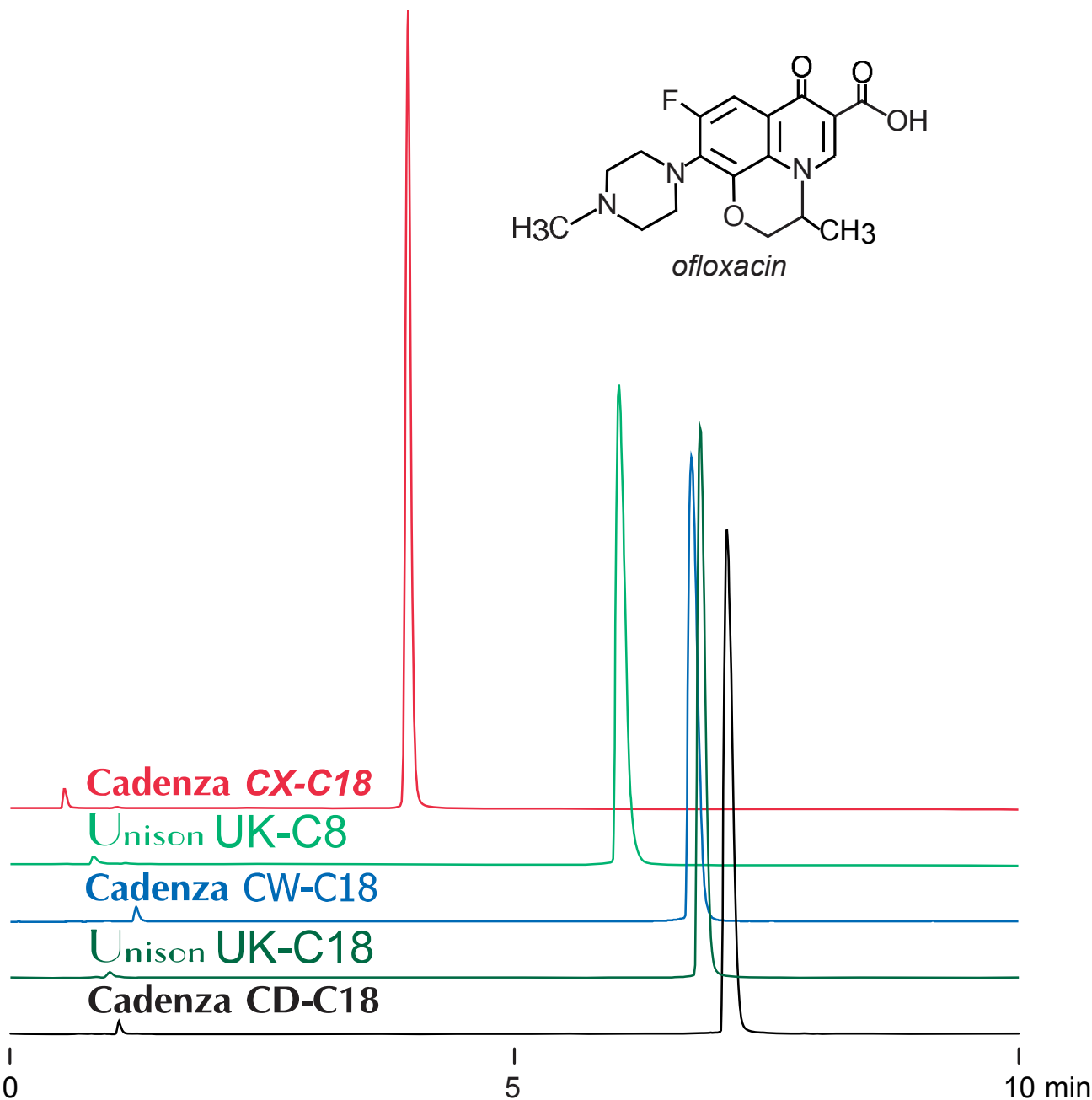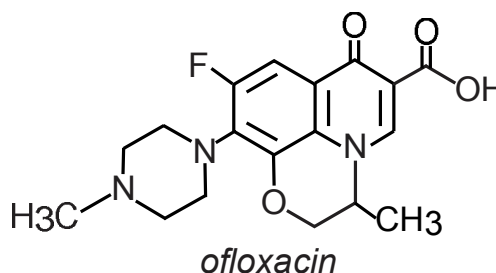

75 x 4.6 mm  
 A: water / formic acid = 100 / 0.1  
 B: acetonitrile / formic acid = 100 / 0.1  
 0-30%B (0-10min)  
 1.0 mL/min (7MPa), 37°C, 290nm  
 10 μL (0.1mg/mL, 1%HCCOH)