



Cadenza CX-C18
Cadenza CD-C18
Cadenza CW-C18
Unison UK-C18

75 x 4.6 mm

Application

Lomefloxacin (fluoroquinolone)

ロメフロキサシン(キノロン系抗菌薬)

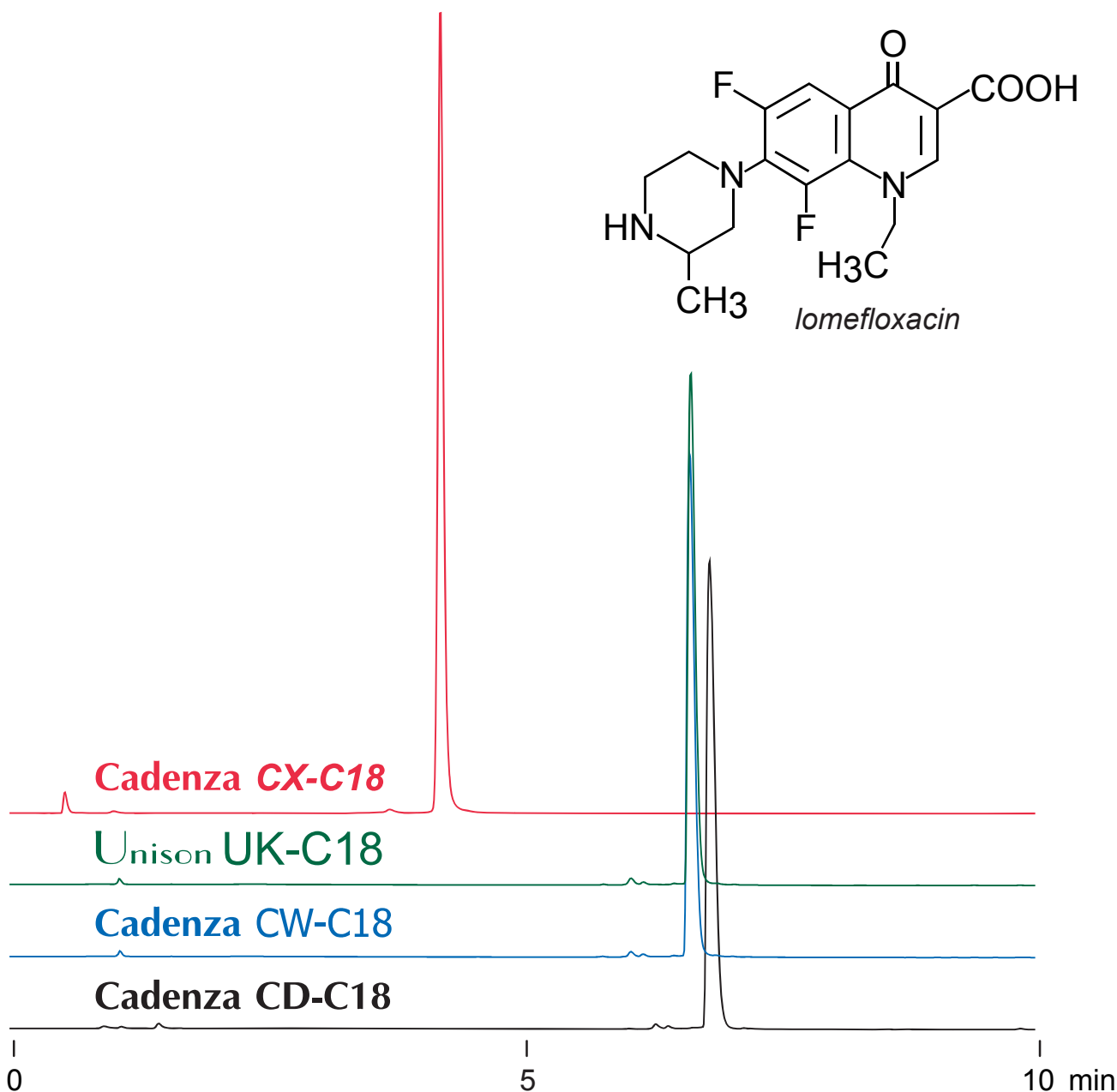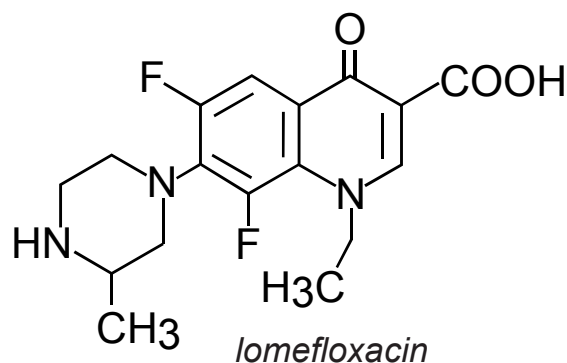

75 x 4.6 mm

A: water / formic acid= 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-30%B (0-10min)

1.0 mL/min (7MPa), 37°C, 290nm

10 µL (0.1mg/mL, 1%HCCOH)