



Cadenza **CX-C18**  
 Unison **UK-C8**  
 Cadenza **CW-C18**  
 Unison **UK-C18**  
 Cadenza **CD-C18**

75 x 4.6 mm

Application

**Oxolinic acid (fluoroquinolone)**

オキソリン酸(キノロン系抗菌薬)

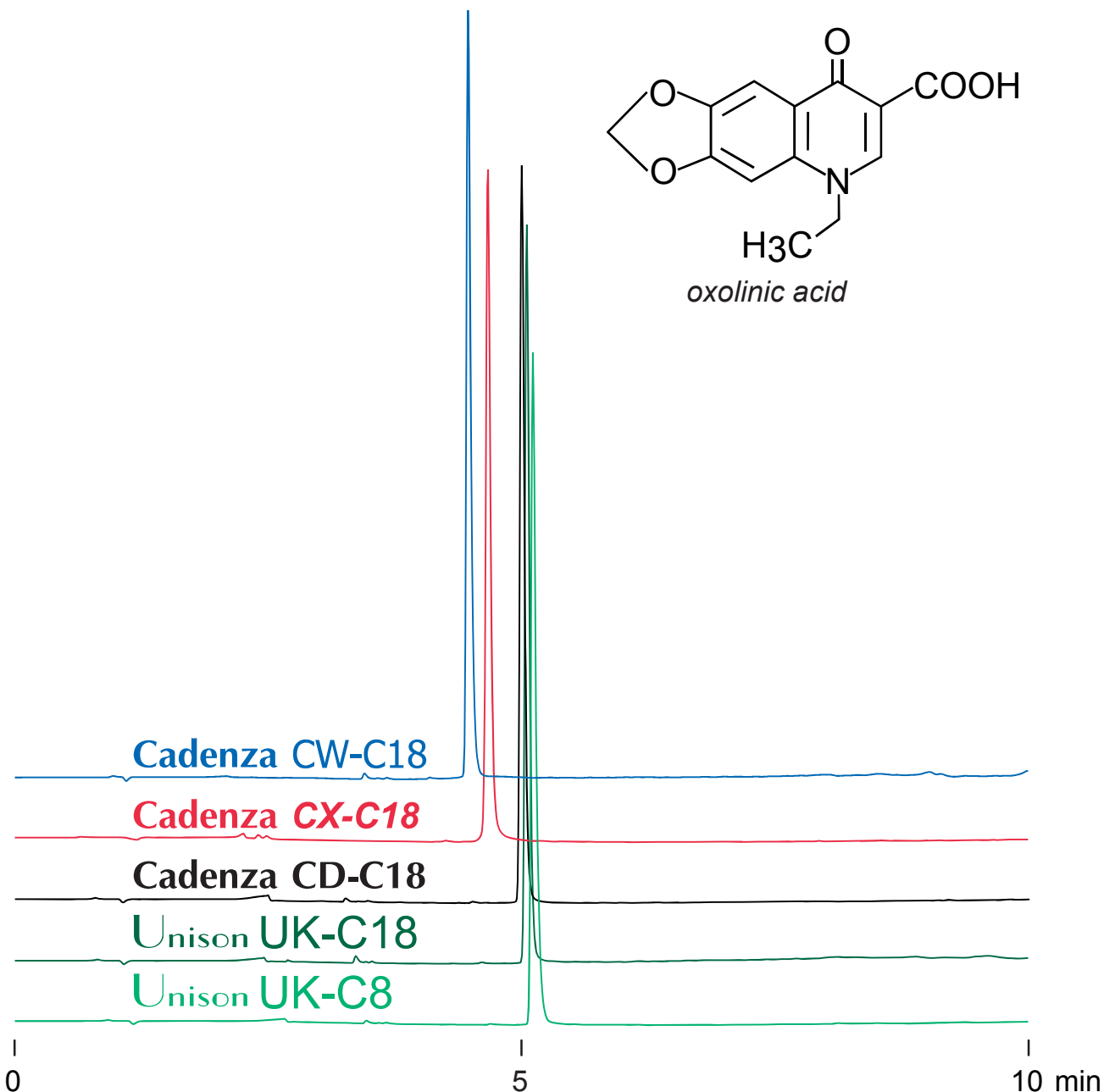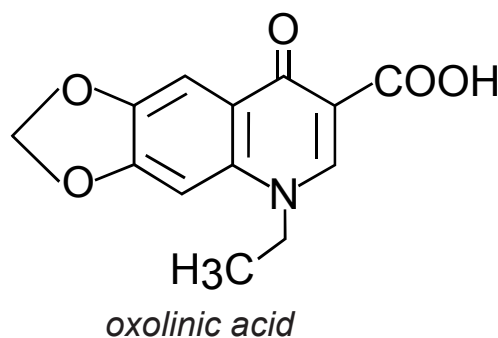

75 x 4.6 mm

A: water / formic acid= 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-10min)

1.0 mL/min (7MPa), 37°C, 290nm

10 µL (0.1mg/mL, 50mM NH<sub>4</sub>OH)