


Cadenza CD-C18

250 x 4.6 mm

Application

**Galactosylceramide from bovine brain**  
 スフィンゴ糖脂質 (ガラクトセレブロシド, 牛脳由来)

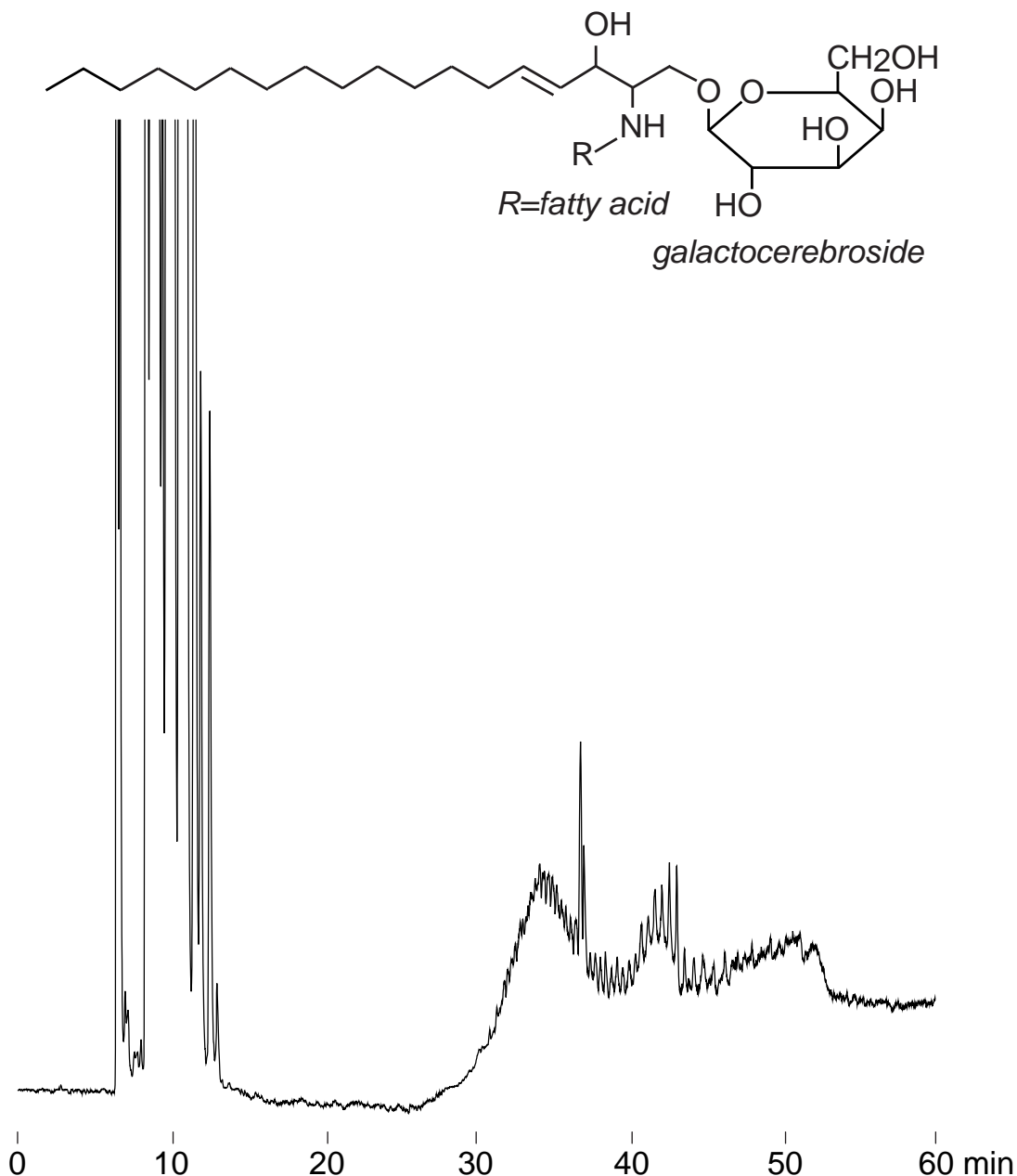

Cadenza CD-C18, 250 x 4.6 mm  
 A: acetonitrile / formic acid = 100 / 0.3  
 B: tetrahydrofuran  
 25-27%B (0-12min), 27-100%B (12-50min), 100%B (50-60min)  
 1mL/min (8MPa), 37°C, ELSD (spray chamber 50°C, drift chamber 100°C)  
 5µL (4mg/mL, chloroform / methanol = 2 / 1)