

Unison UK-Phenyl

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

Application

Bromhexine (Expectorant drug)

ブロムヘキシン (去痰薬)

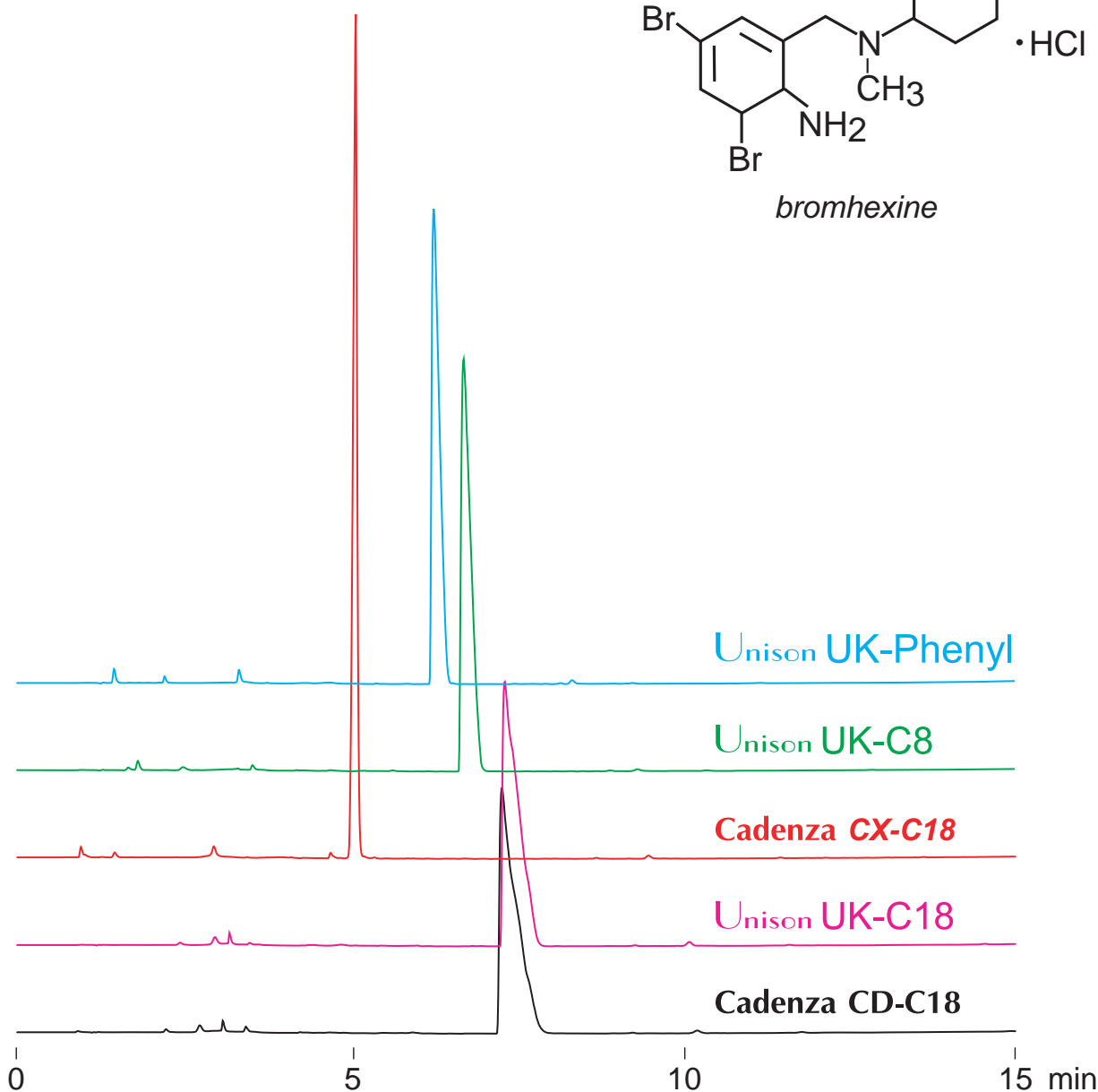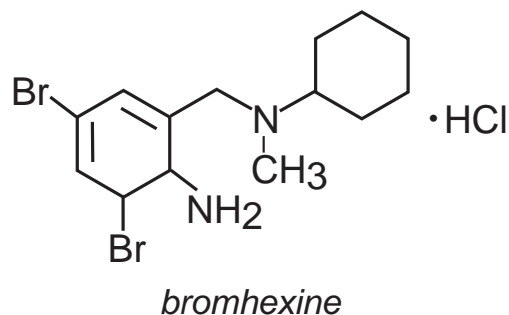

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-15min)

1.0mL/min (6MPa), 37°C, 260nm

10µL (1mg/mL, water / formic acid = 100 / 0.1)