



Cadenza CD-C18

100 x 4.6 mm

Application

*Nisoldipine (Antihypertensive, Calcium channel blocker)*

ニソルジピン ( 血圧降下薬 , カルシウム拮抗薬 )

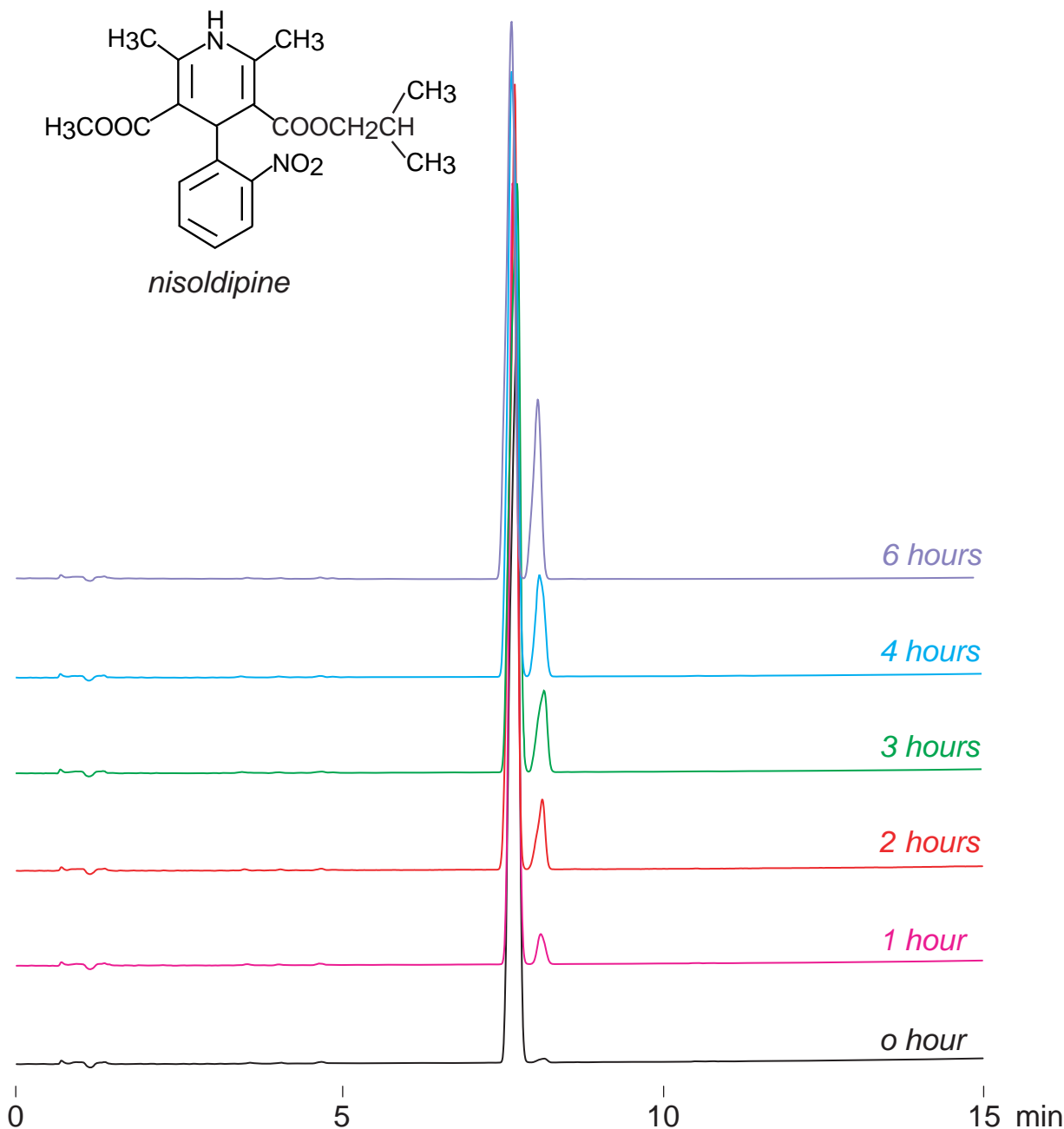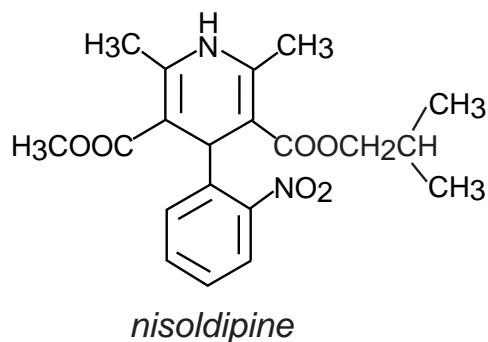

Cadenza CD-C18, 100 x 4.6 mm  
 A: water / formic acid = 100 / 0.1  
 B: acetonitrile / formic acid = 100 / 0.1  
 40-100%B (0-15min)  
 1.0mL/min (8MPa), 37°C, 260nm  
 10µL (1mg/mL, methanol)