

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Loperamide (Antidiarrheal, Opioid receptor agonist)

ロペラミド (止瀉剤 , オピオイド受容体作動薬)

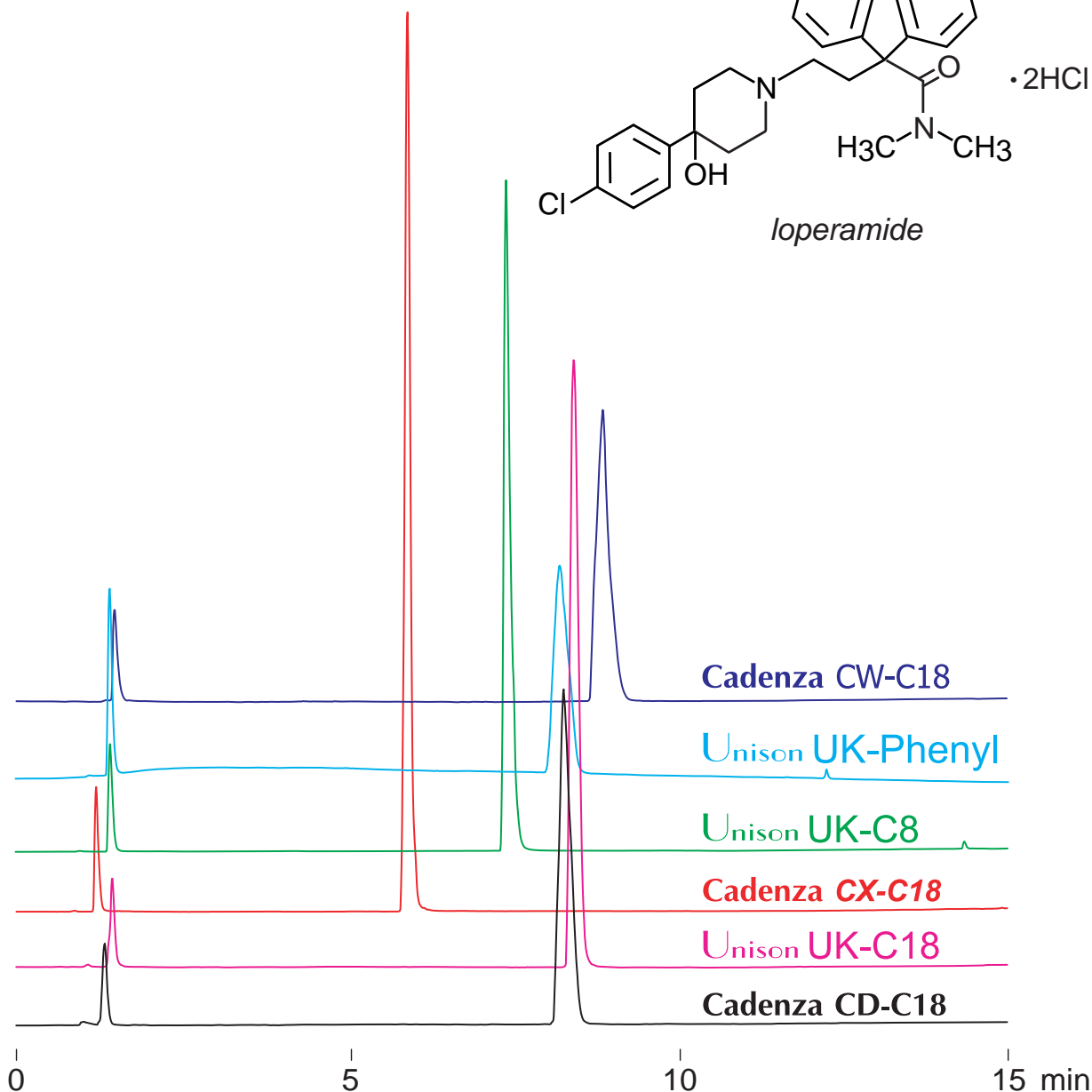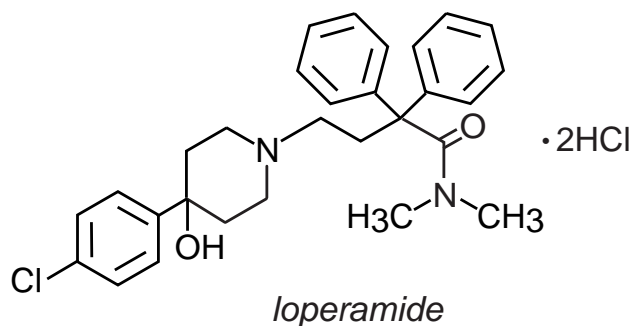

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C

ELSD (spray chamber 50°C, drift chamber 100°C)

10µL (0.5mg/mL, ethanol)