



Unison UK-Phenyl Cadenza CX-C18

Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm HPLC COLUMN APPLICATION

Zaltoprofen (Anti-inflammatory, Cyclooxygenase inhibitor)

ザルトプロフェン (鎮痛薬 , シクロオキシゲナーゼ 2 阻害薬)

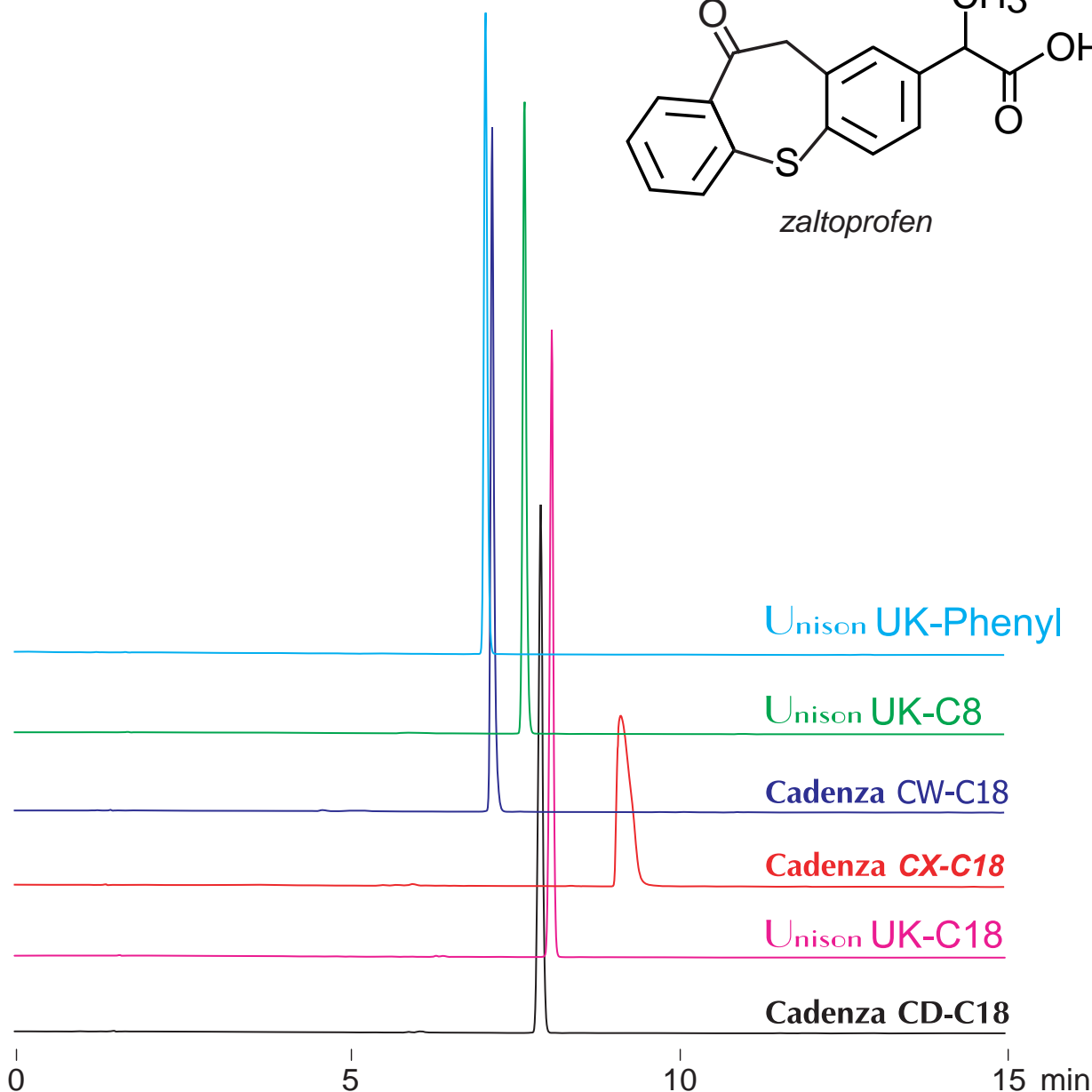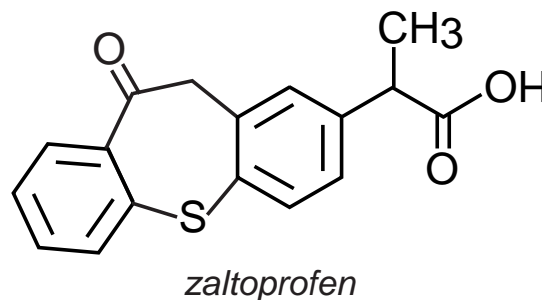

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate= 90 / 10

10-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate)