

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Pioglitazone (Antidiabetic)

ピオグリタゾン (糖尿病治療薬)

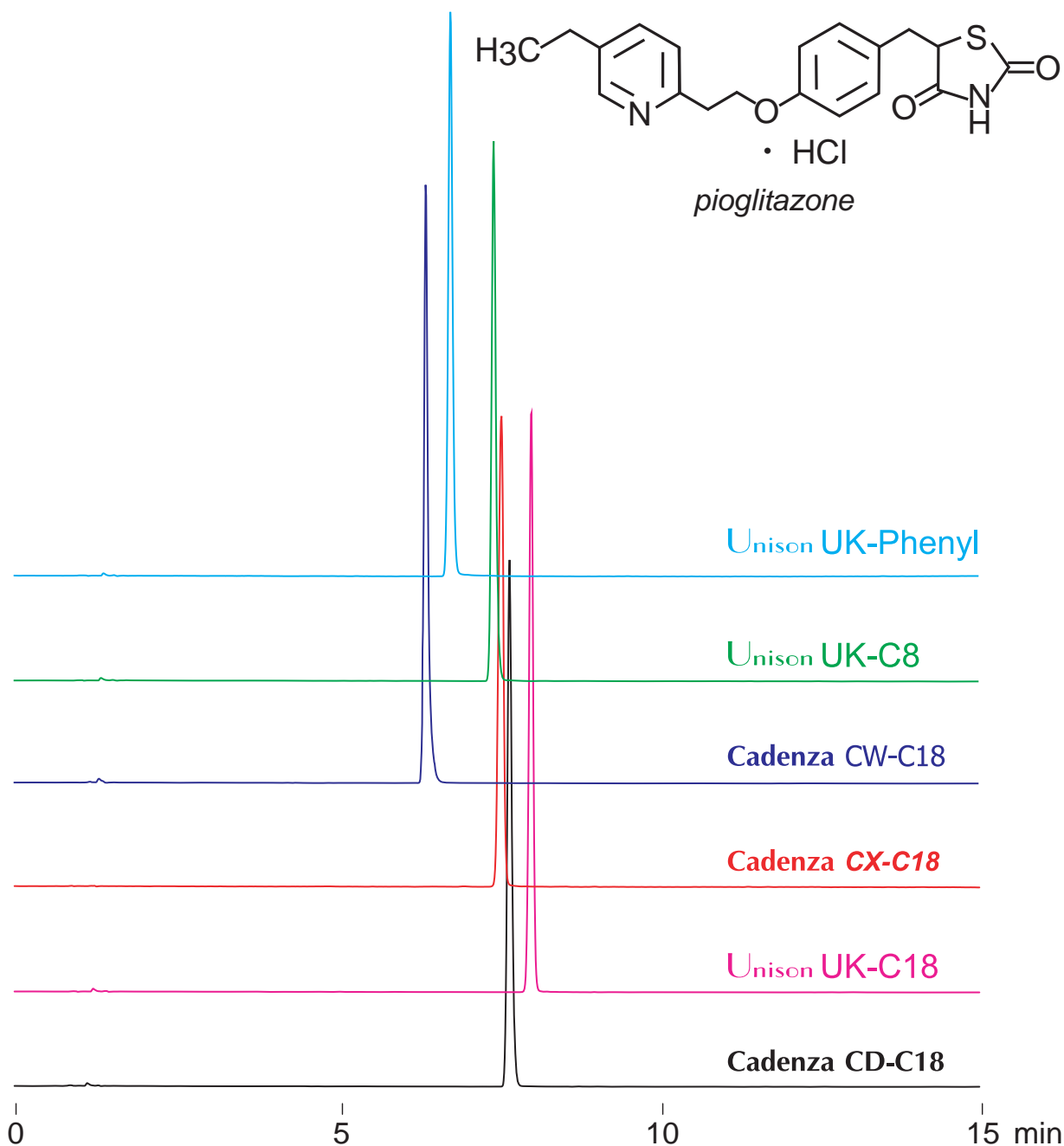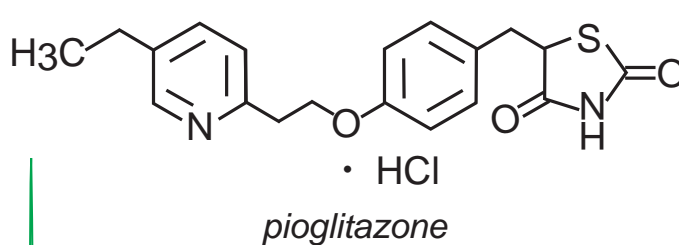

100 x 4.6 mm
 A: 100mM ammonium formate
 B: acetonitrile / 100mM ammonium formate= 90 / 10
 30-100%B (0-15min)
 1.0mL/min (8MPa), 37°C, 260nm
 10μL (0.5mg/mL, water/ formic acid = 100 / 1)