

Unison UK-Phenyl Cadenza CX-C18


Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Mosapride (Prokinetic)

モサプリド (消化管運動改善薬)

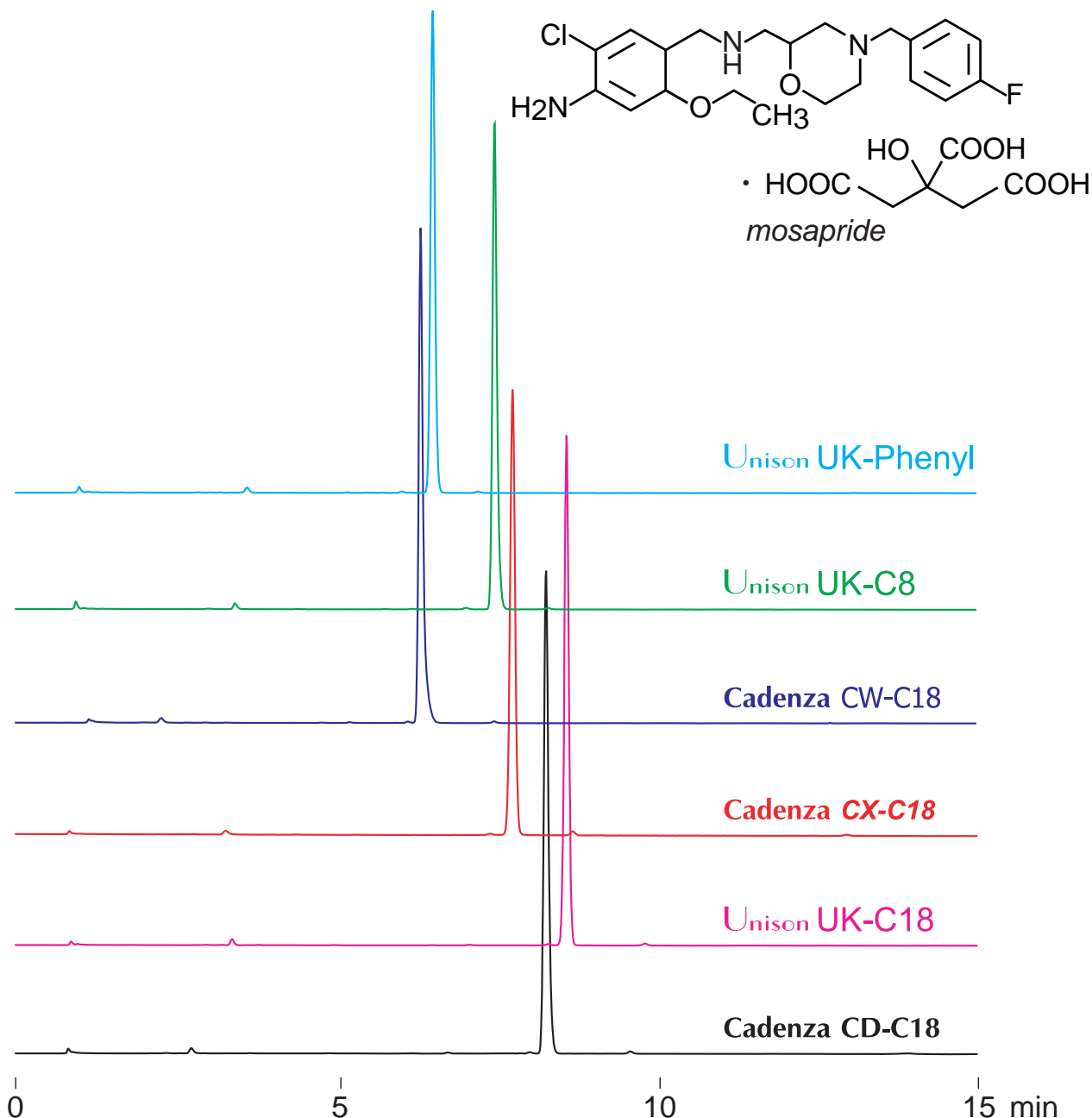

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

35-100%B (0-15min)

1.0 mL/min (8MPa), 37°C, 260nm

10µL (0.5mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)