

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

## Amlodipine ( Antidepressant )

アムロジピン ( 血圧降下薬 )

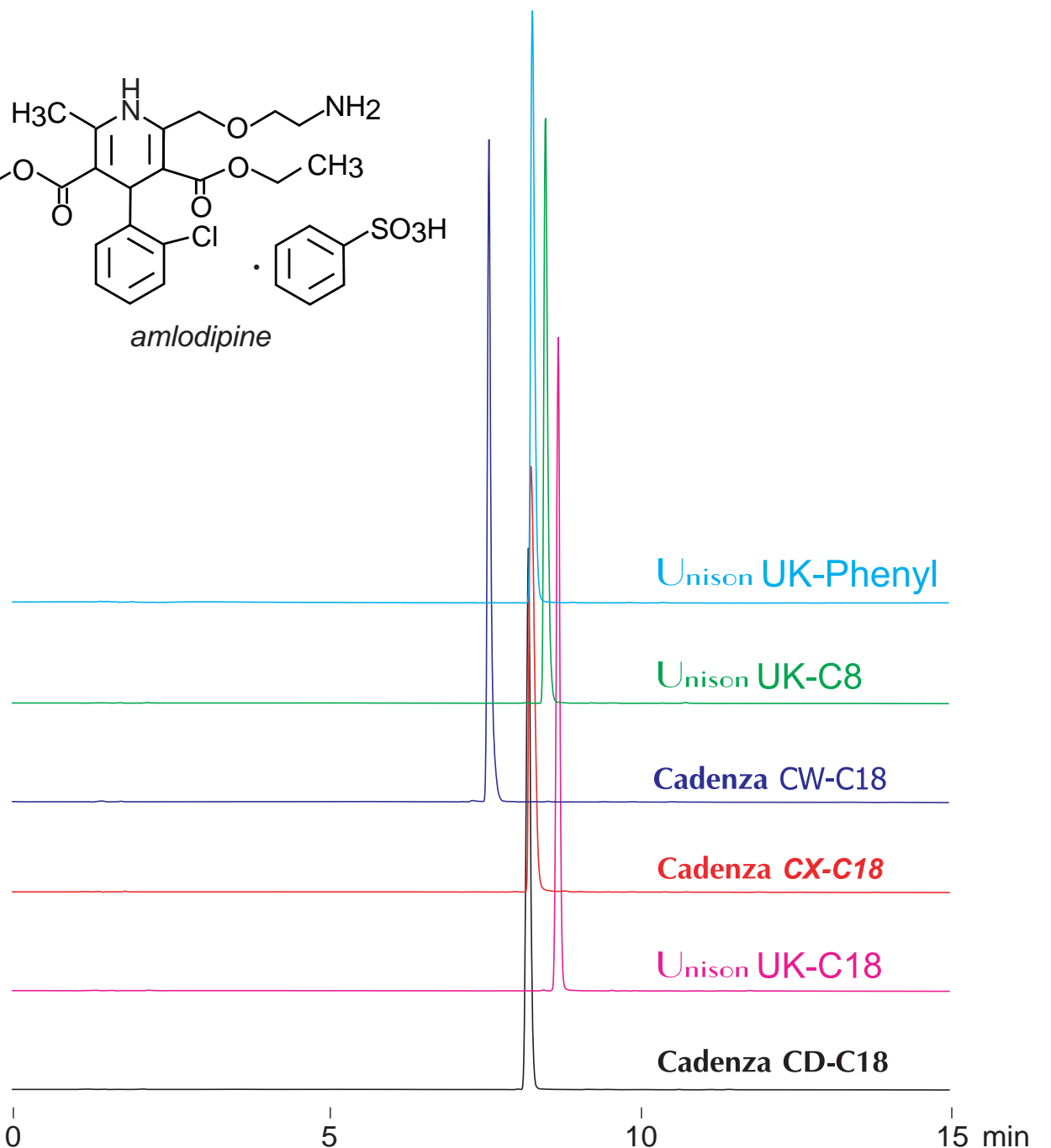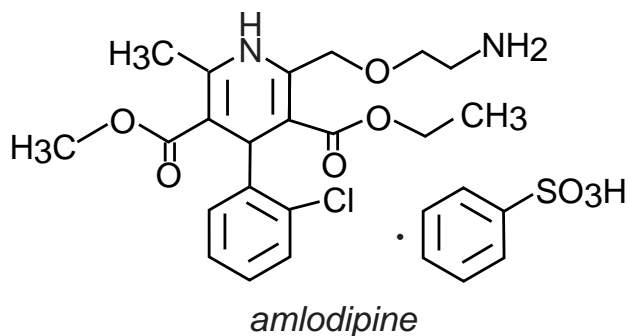

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate= 90 / 10

10-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.5mg/mL, methanol / formic acid = 100 / 1)