

Unison UK-Phenyl Cadenza CX-C18



Unison UK-C8 Cadenza CW-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Rebamipide (Gastrointestinal agent)

レバミピド (消化性潰瘍薬)

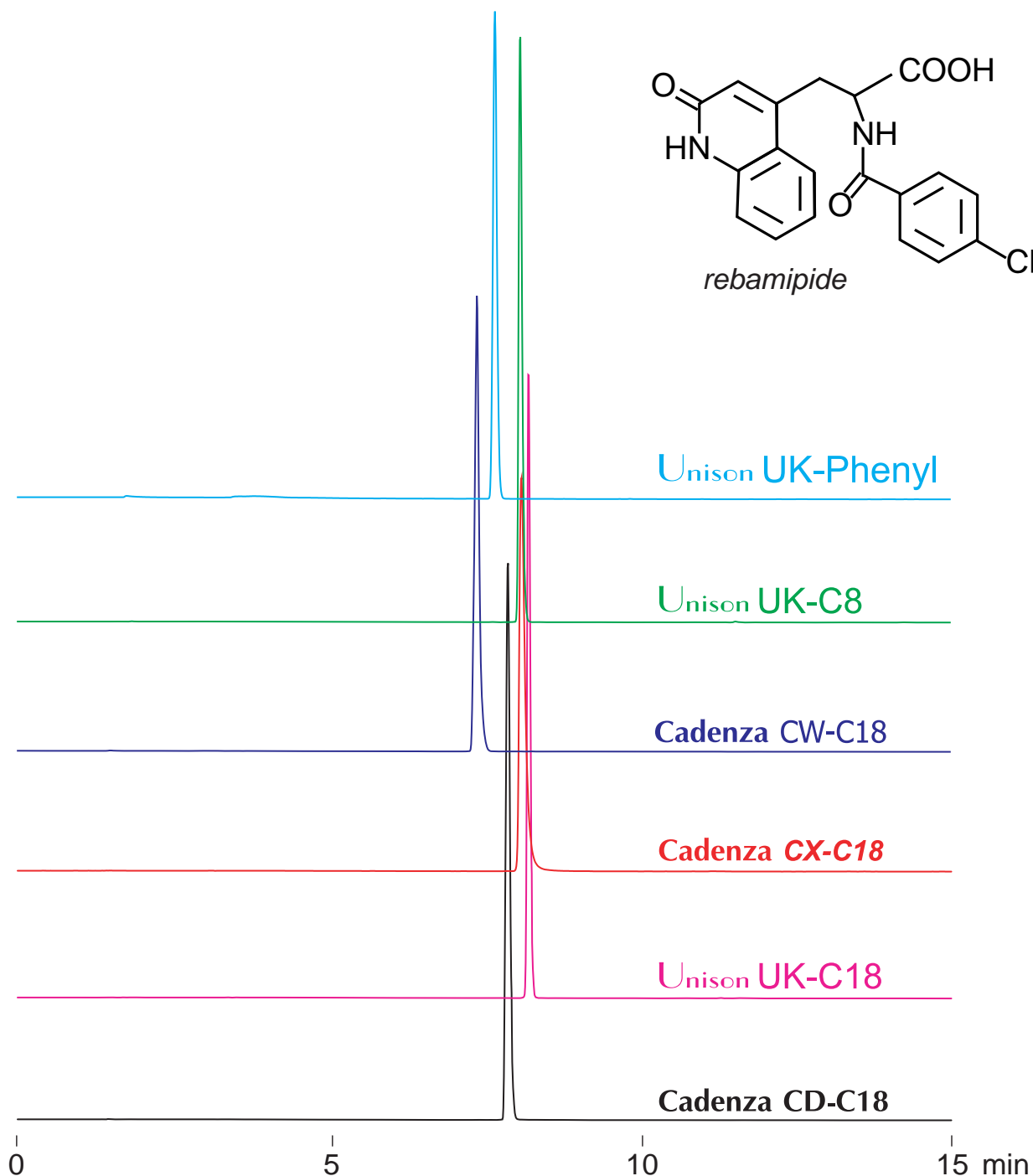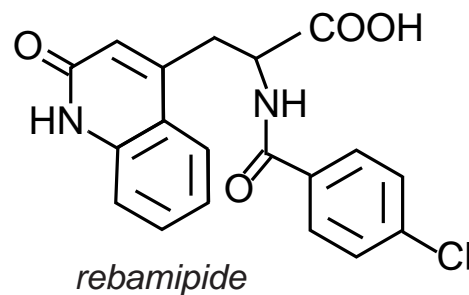

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

0-70%B (0-15min)

1.0mL/min (10MPa), 37°C, 260nm

10µL (0.5mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)