

Unison UK-C1

Unison UK-Phenyl Cadenza CX-C18

Unison UK-C8 Cadenza CW-C18


Unison UK-C18 Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

### Cetylpyridinium ( Disinfectant )

塩化セチルピリジニウム ( 口腔殺菌消毒薬 )

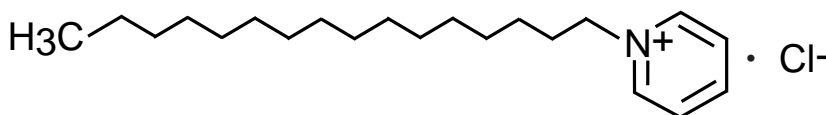

cetylpyridinium chloride

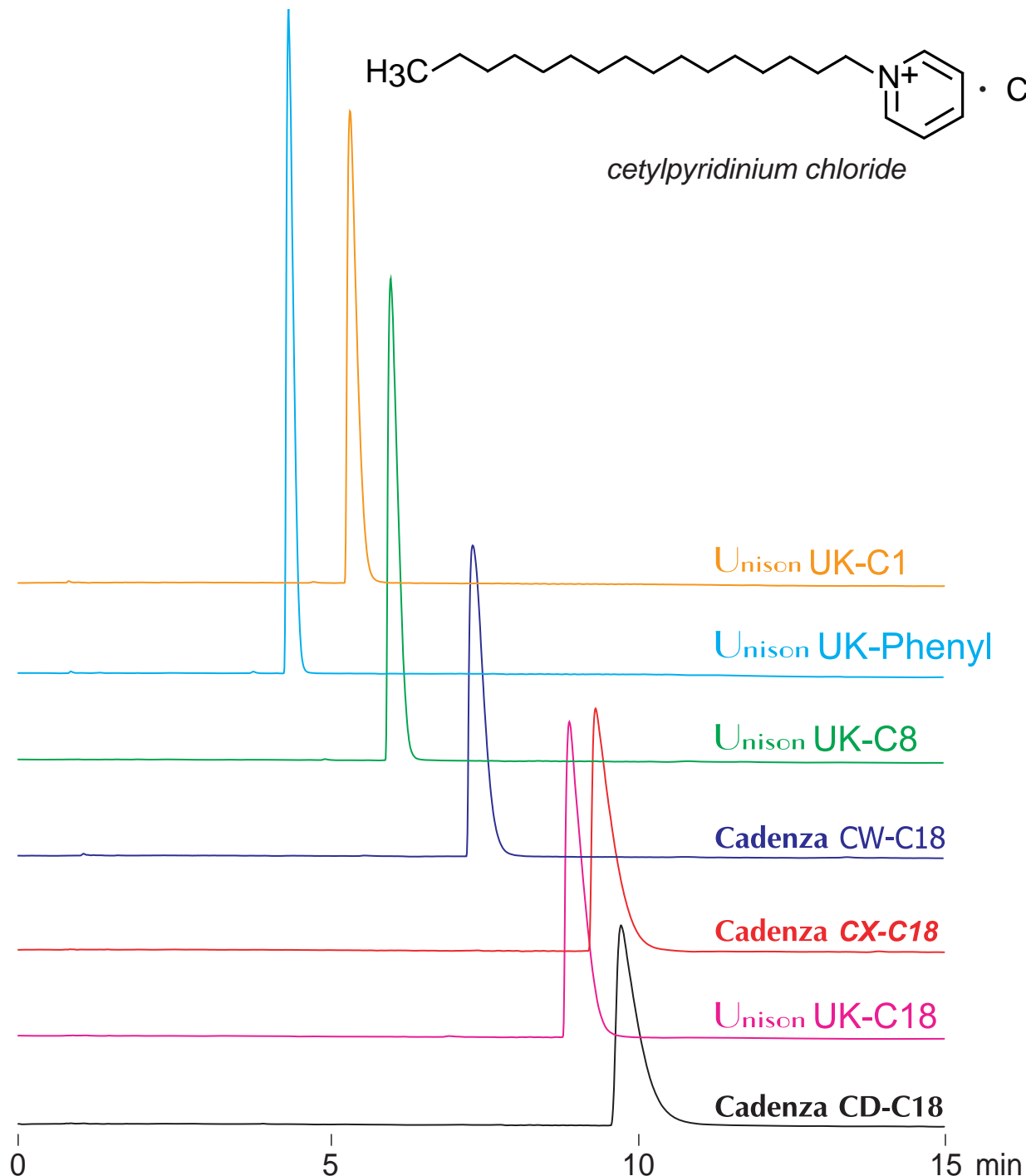

100 x 4.6 mm

A: acetonitrile / 100mM ammonium formate= 50 / 50

B: acetonitrile / 100mM ammonium formate= 90 / 10

0-100%B (0-10min), 100%B (10-15min)

1.0mL/min (7MPa), 37°C, 260nm

10µL (0.5mg/mL, 100mM ammonium formate)