

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

## Drospirenone ( Dysmenorrhea medicine )

ドロスピレノン ( 月経困難症治療薬 )

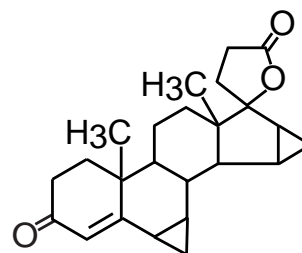

drospirenone

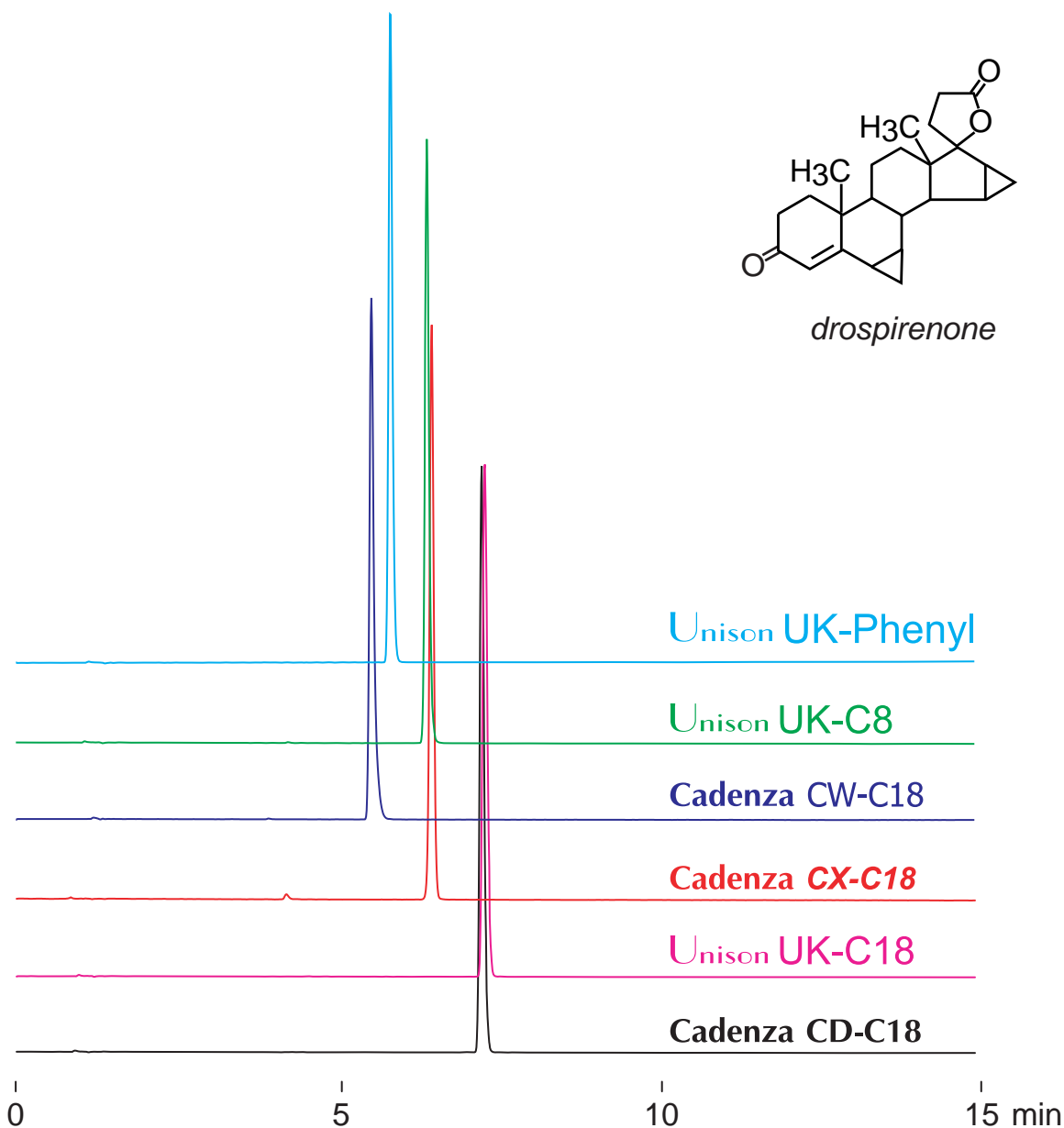

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

30-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 260nm

10µL (0.2mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)