

Unison UK-C8

Unison UK-Phenyl

Cadenza GL-C18

Cadenza CW-C18

Cadenza CD-C18



Unison UK-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Glycyrrhizic acid (Anti-inflammatory)

グリチルリチン酸 (抗炎症薬)

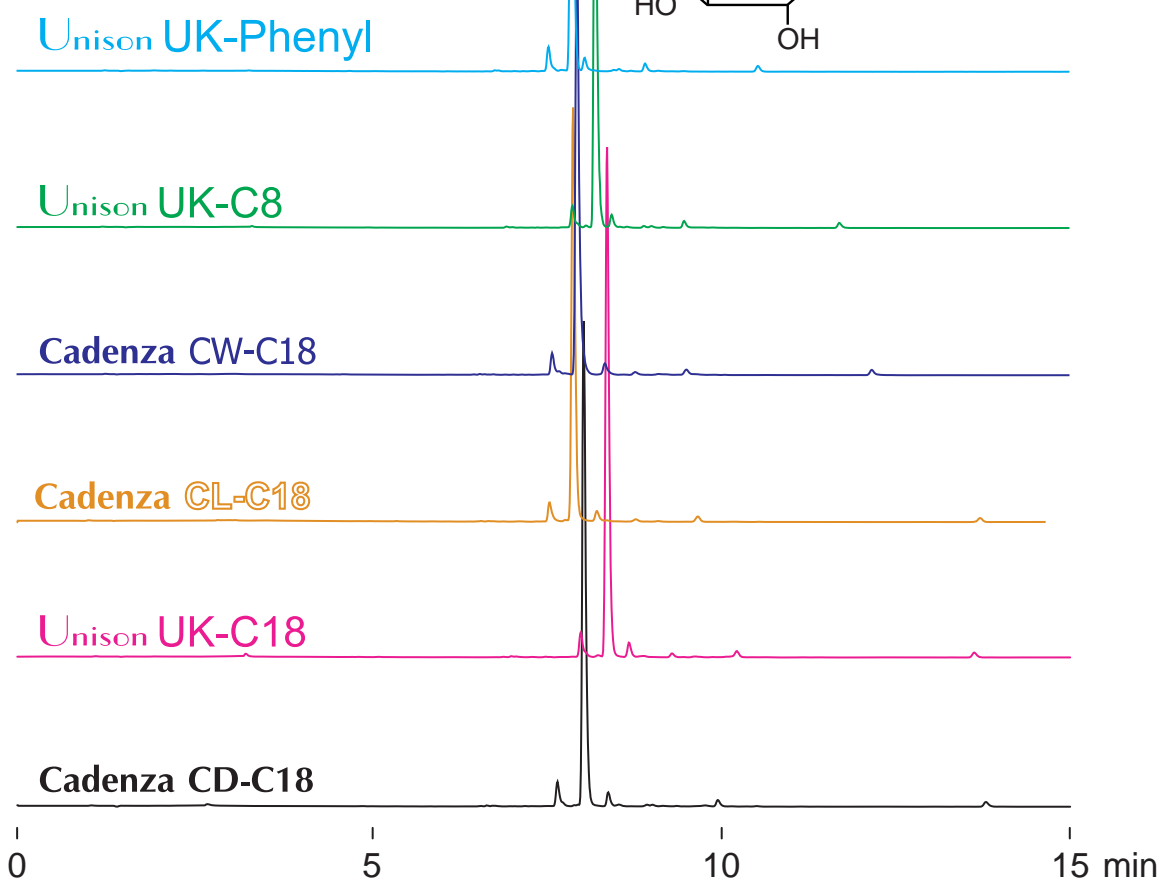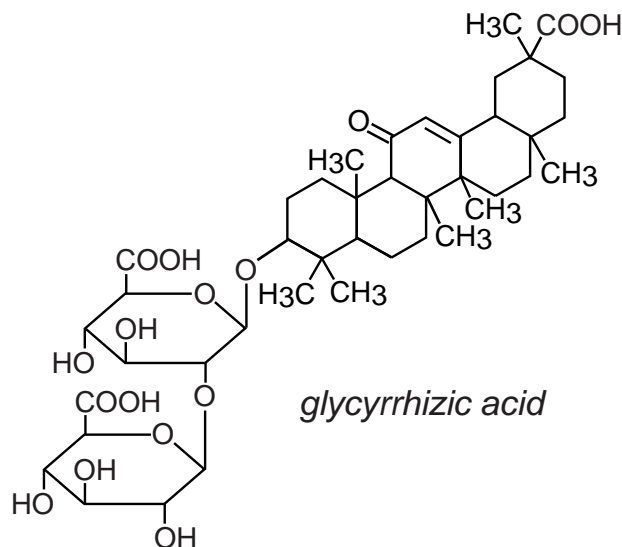

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-100%B (0-15min),

1.0mL/min (9MPa), 37°C, 260nm

10μL (1mg/mL, 100mM ammonium formate)