


Unison UK-Amino

100 x 4.6 mm

HPLC COLUMN APPLICATION

Tranexamic acid ( Anti-inflammatory )

トラネキサム酸 ( 抗炎症薬 )

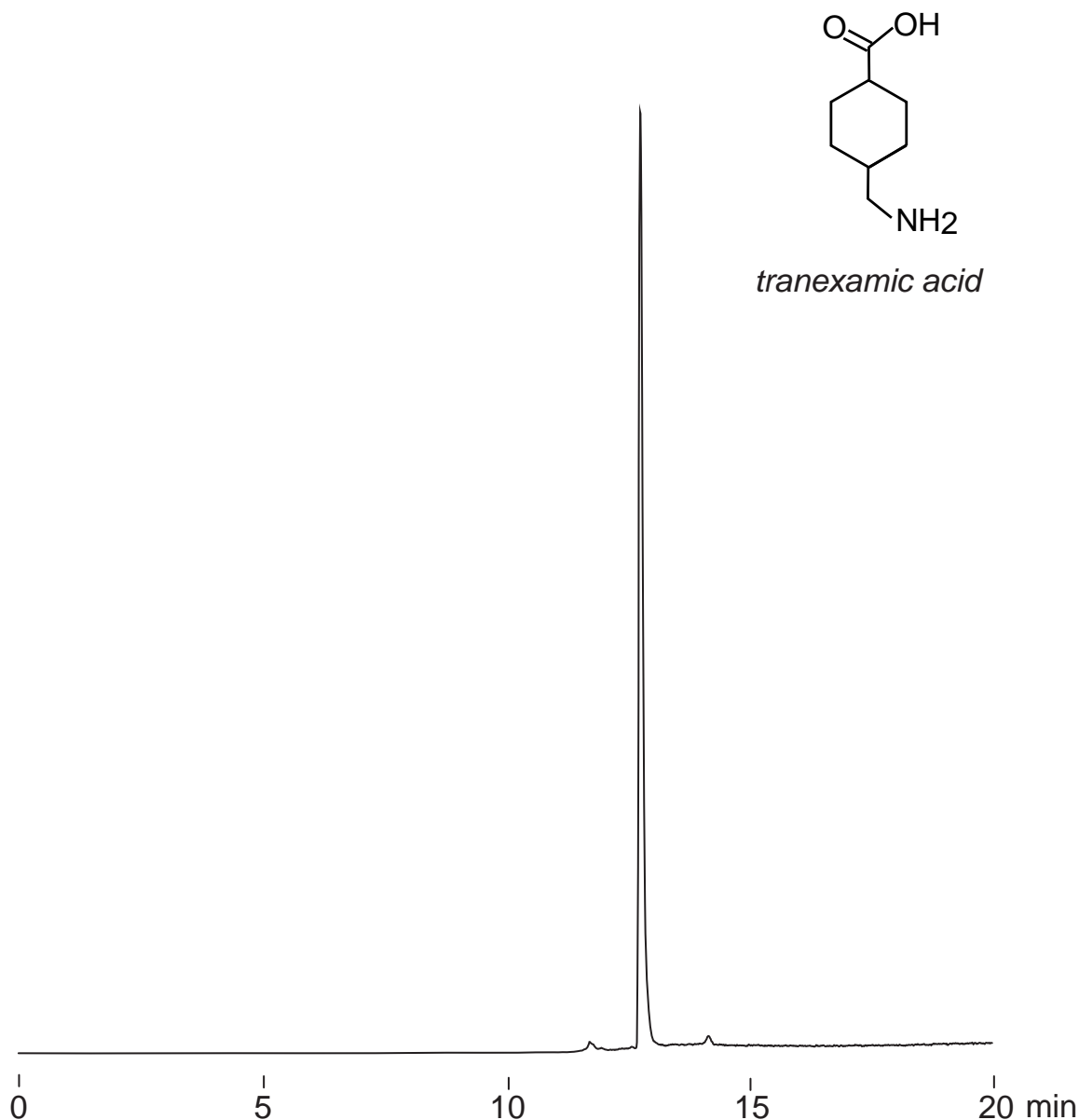

Unison UK-Amino, 250 x 4.6 mm

A: acetonitrile

B: 100mM ammonium formate

0%B (0-2min), 0-100%B (2-15min), 100%B (15-20min)

1.0mL/min (7MPa), 37°C

ELSD (spray chamber 30°C, drift chamber 70°C)

10µL (1mg/mL, 100mM ammonium formate)