

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Alprazolam (Antianxiety)

アルプラゾラム (抗不安薬)

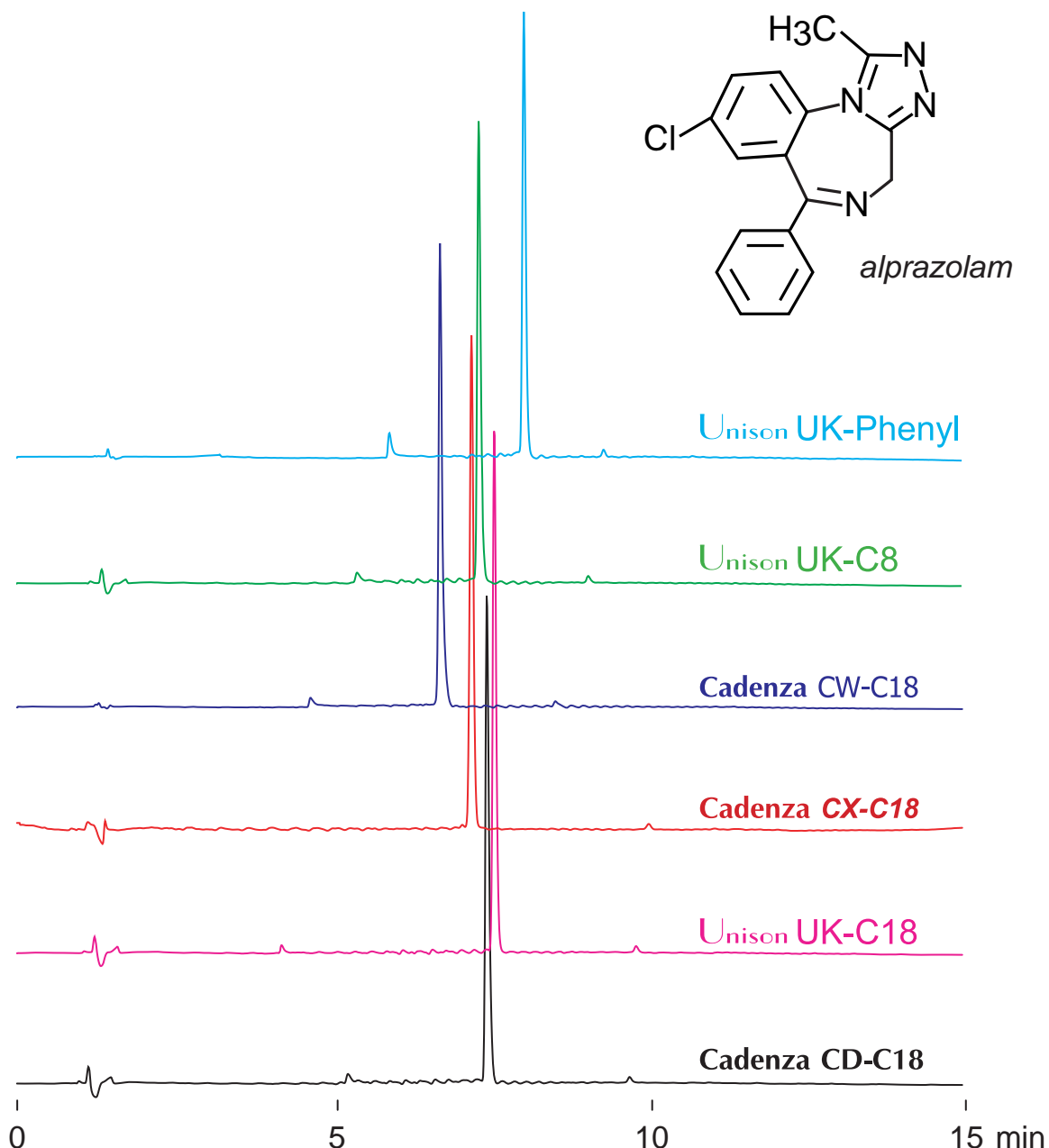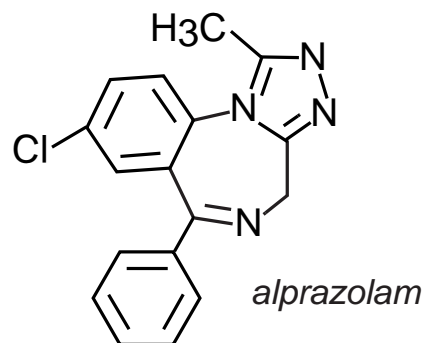

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

10-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 250nm

10µL (0.5mg/mL, methanol / 100mM ammonium formate = 50 / 50)