

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8  
Cadenza CX-C18



Unison UK-C18  
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Alprazolam ( Antianxiety )

アルプラゾラム ( 抗不安薬 )

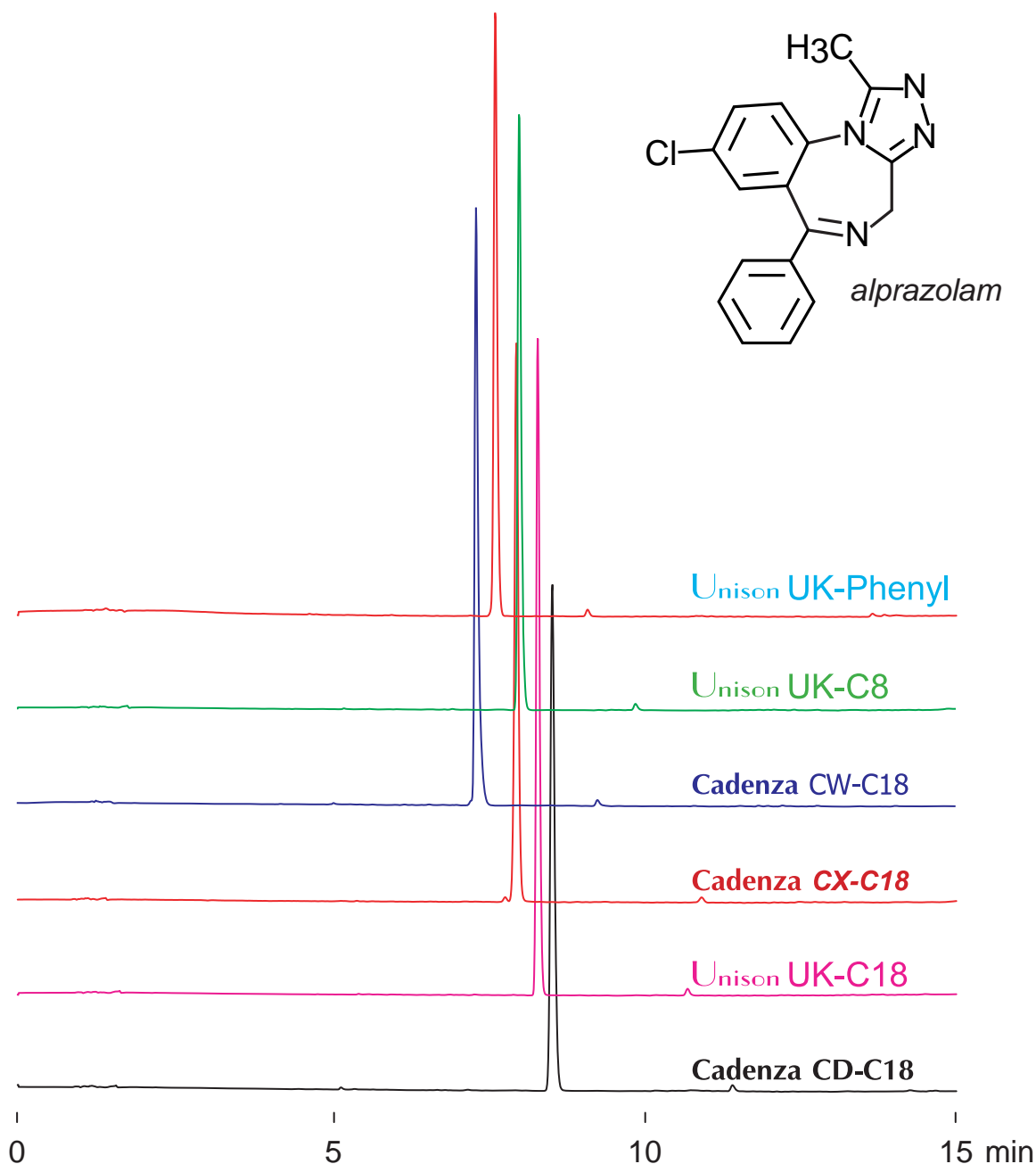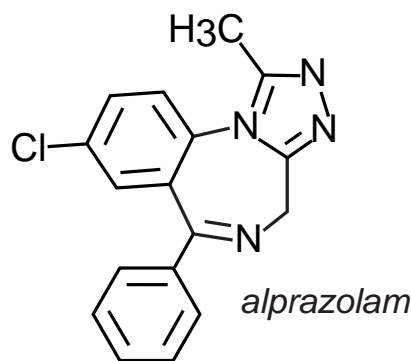

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

10-100%B (0-15min)

1.0mL/min (8MPa), 37°C, 250nm

10µL (0.5mg/mL, methanol / 100mM ammonium formate = 50 / 50)