

Cadenza CW-C18

Unison UK-C18 Unison UK-Phenyl



Cadenza CD-C18 Unison UK-C8

100 x 4.6 mm

HPLC COLUMN APPLICATION

Picosulfate (Laxative)

ピコスルファート (緩下薬)

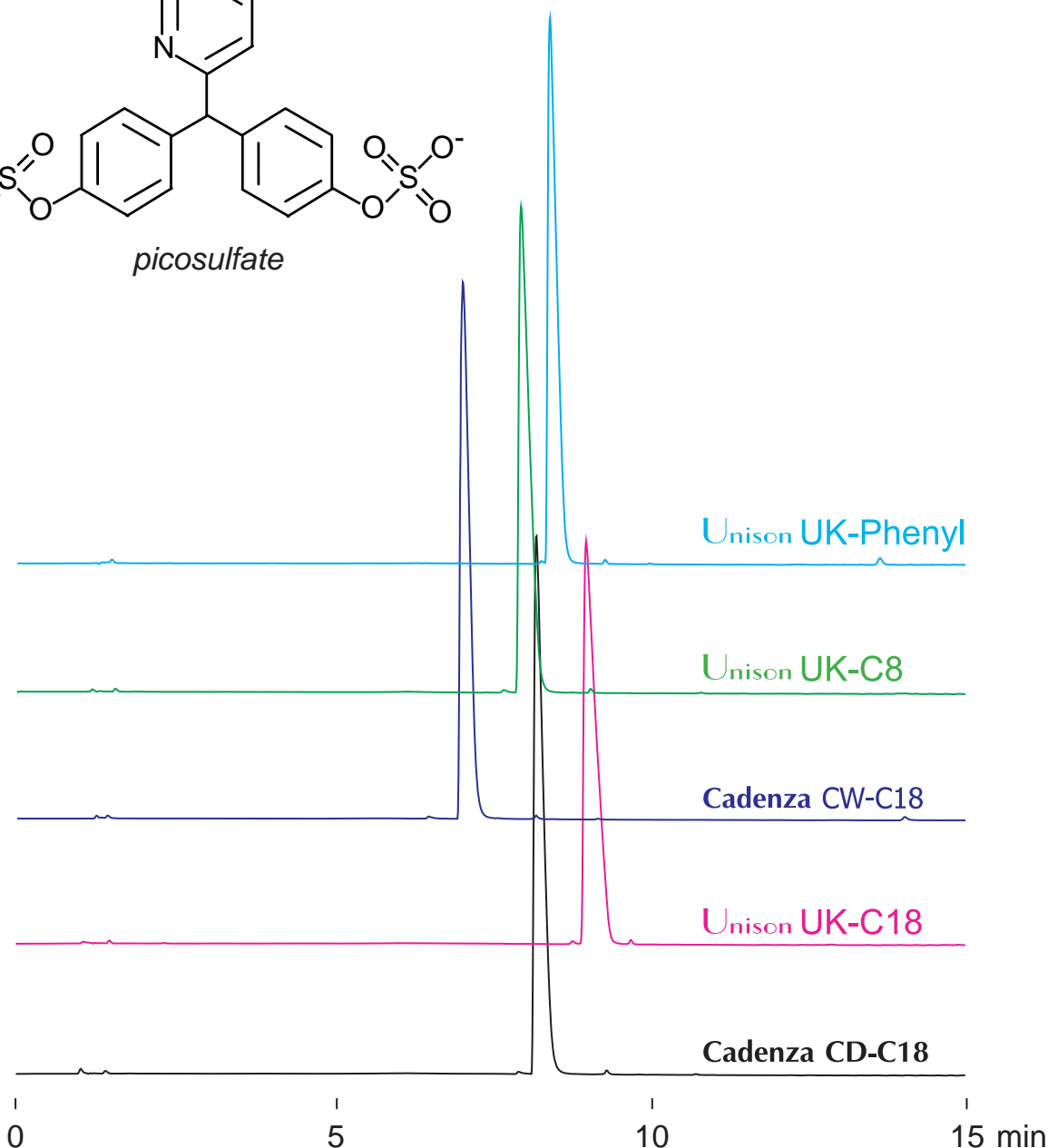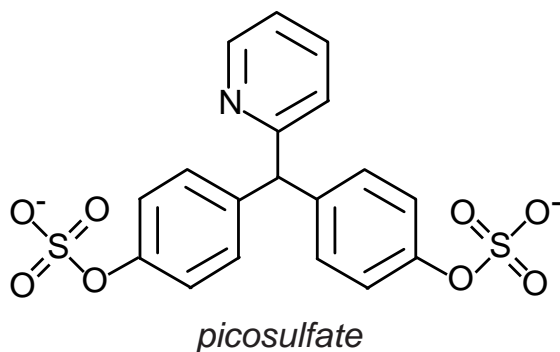

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0%B (0-3min), 50%B (3-15min)

1.0mL/min (8MPa), 37°C, 260nm

10μL (1mg/mL, water / formic acid = 100 / 0.1)