

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18


Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Azoxybenzene ( Liquid-crystal materials )

アゾキシベンゼン ( 液晶材料 )

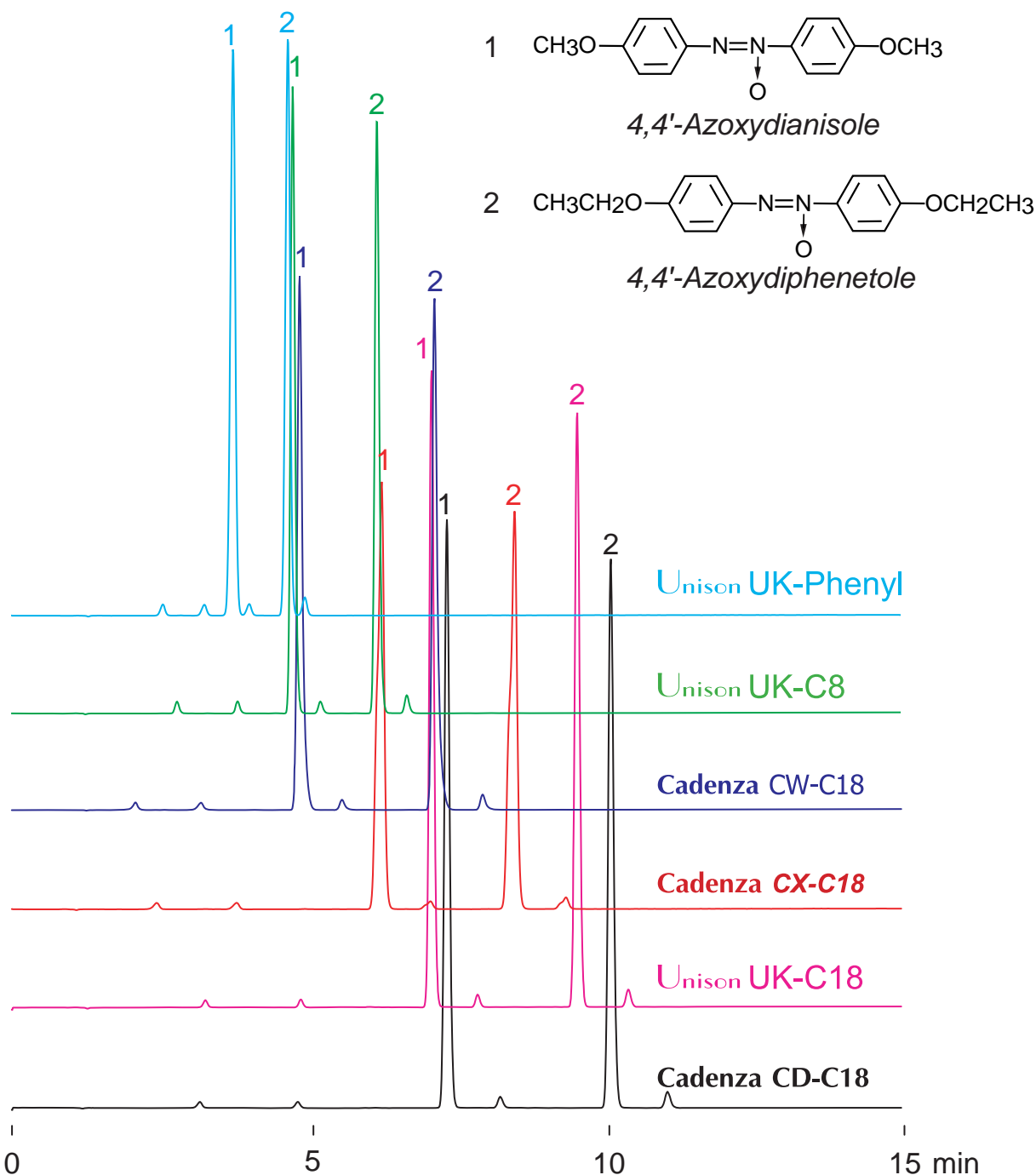

100 x 4.6 mm

A: acetonitrile / water / acetic acid = 50 / 50 / 0.5

B: acetonitrile / acetic acid = 100 / 0.5

0-100%B (0-15min)

1.0 mL/min (7MPa), 37°C, 360nm

10µL (0.1mg/mL, acetonitrile)