

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Unison UK-C18

Cadenza CX-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Fexofenadine (Allergic rhinitis, Opioid receptor agonist)

フェキソフェナジン (アレルギー鼻炎薬 , ヒスタミン受容体拮抗薬)

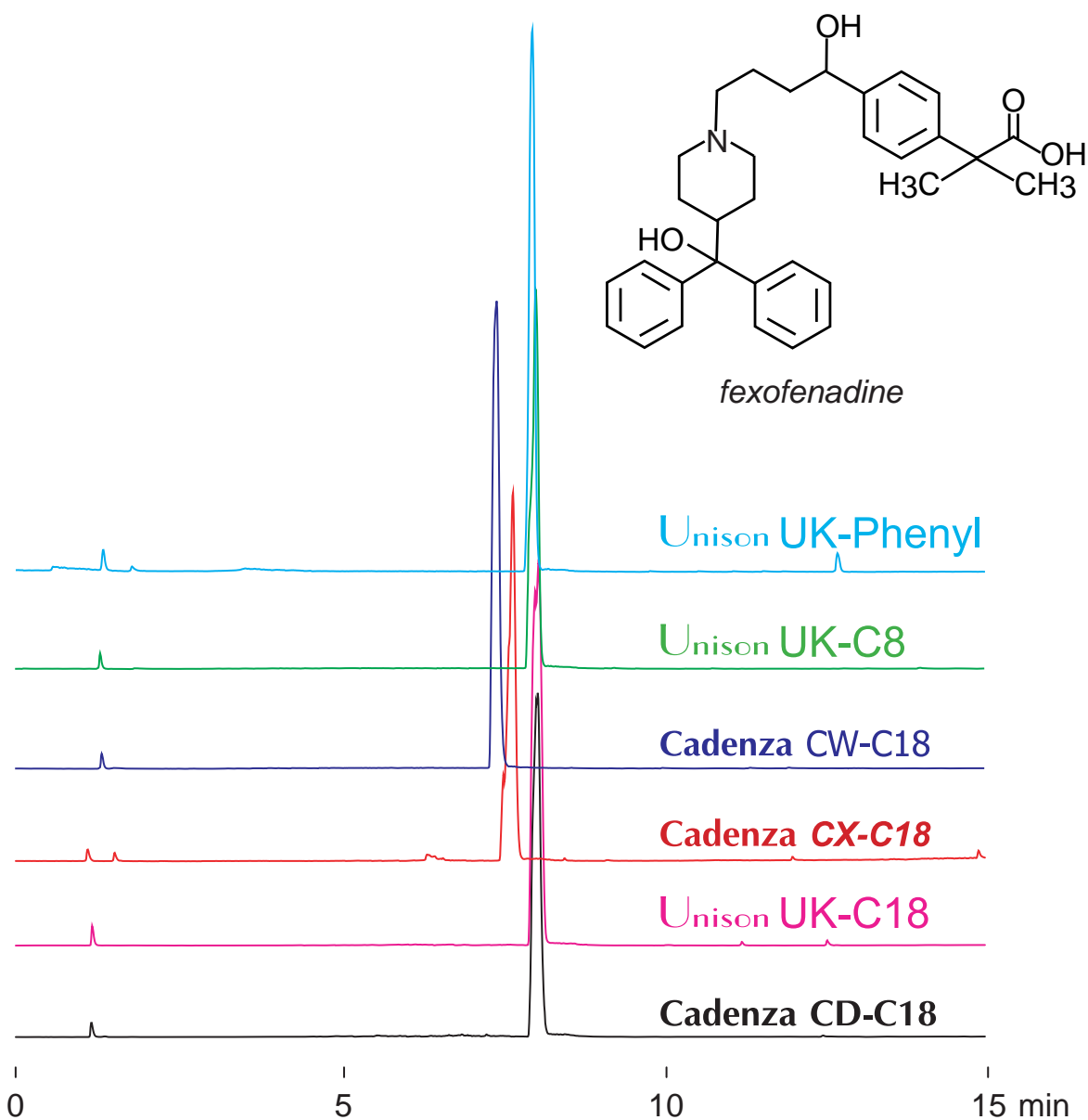

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

0-100%B (0-15min)

1.0mL/min (8MPa), 37°C

ELSD (spray chamber 30°C, drift chamber 60°C)

10µL (0.25mg/mL, acetonitrile / formic acid = 100 / 1)