


Scherzo SS-C18
 Scherzo SM-C18
 Scherzo SW-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Fexofenadine (Allergic rhinitis, Opioid receptor agonist)

フェキソフェナジン (アレルギー鼻炎薬 , ヒスタミン受容体拮抗薬)

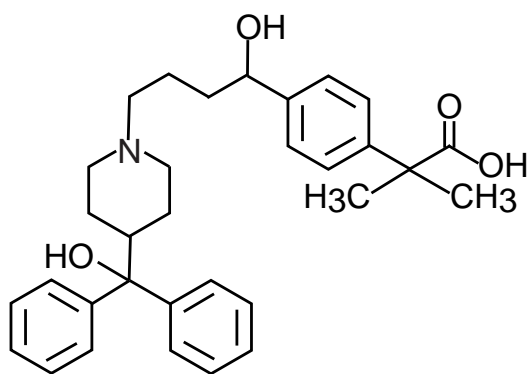

fexofenadine

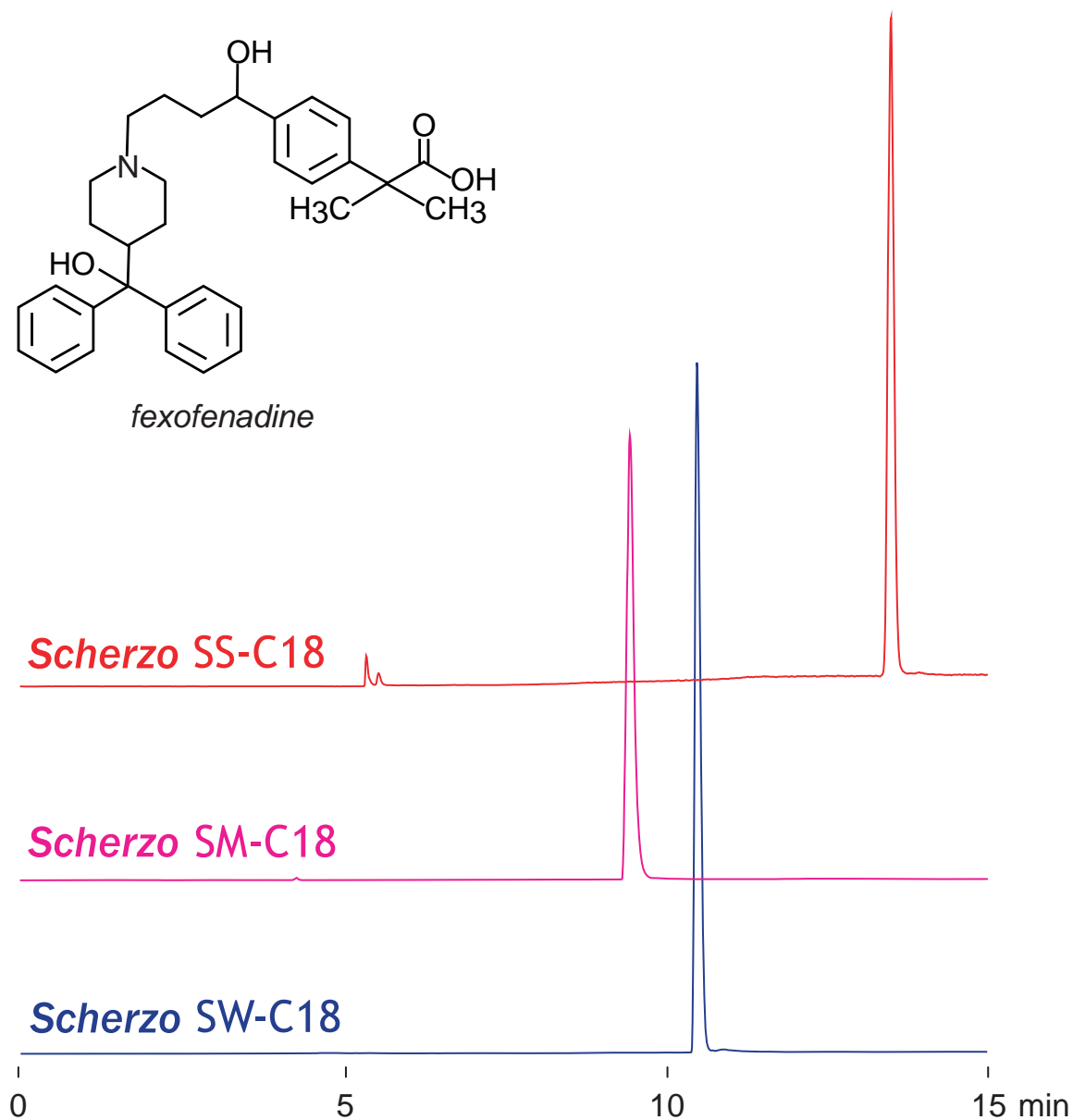

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / 100mM ammonium formate = 50 / 50

0-100%B (0-10min), 100%B(10-15min)

1.0mL/min (10MPa), 37°C

ELSD (spray chamber 30°C, drift chamber 90°C)

10µL (0.25mg/mL, methanol)