


DACAPO DX-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Fursultiamine ( Vitamin B1 derivative )

フルスルチアミン ( ビタミン B1 誘導体 )

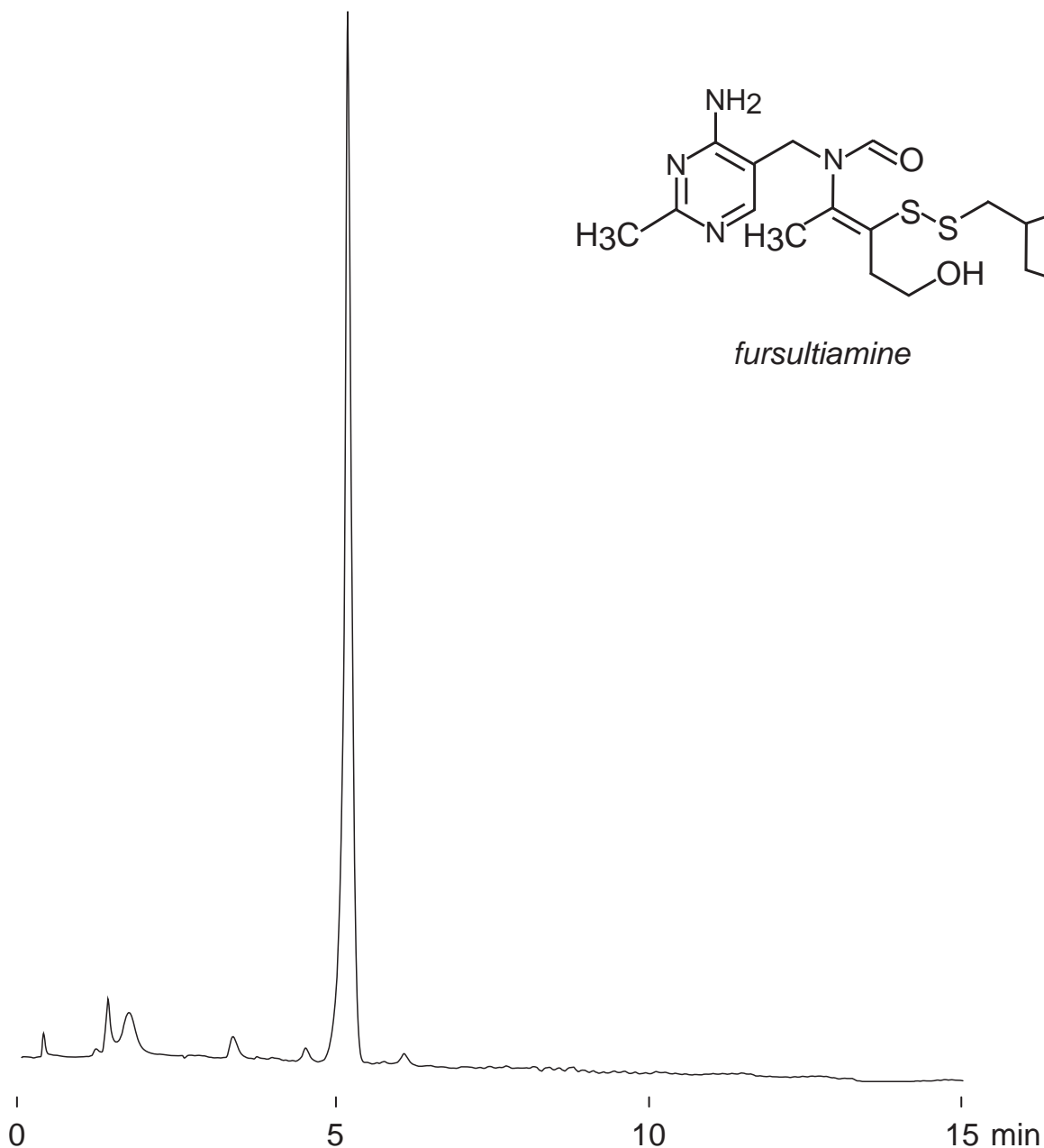

DACAPO DX-C18, 100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

0-100%B (0-15min)

1.0 mL/min (14MPa), 37°C, 260nm

10µL (0.5mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)