

Cadenza CW-C18

Cadenza CX-C18

Unison UK-C18

Cadenza CD-C18

Unison UK-Phenyl

Unison UK-C8



DACAPO DX-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Basic Violet 1 (Dye)

ベーシックバイオレット1 (染料)

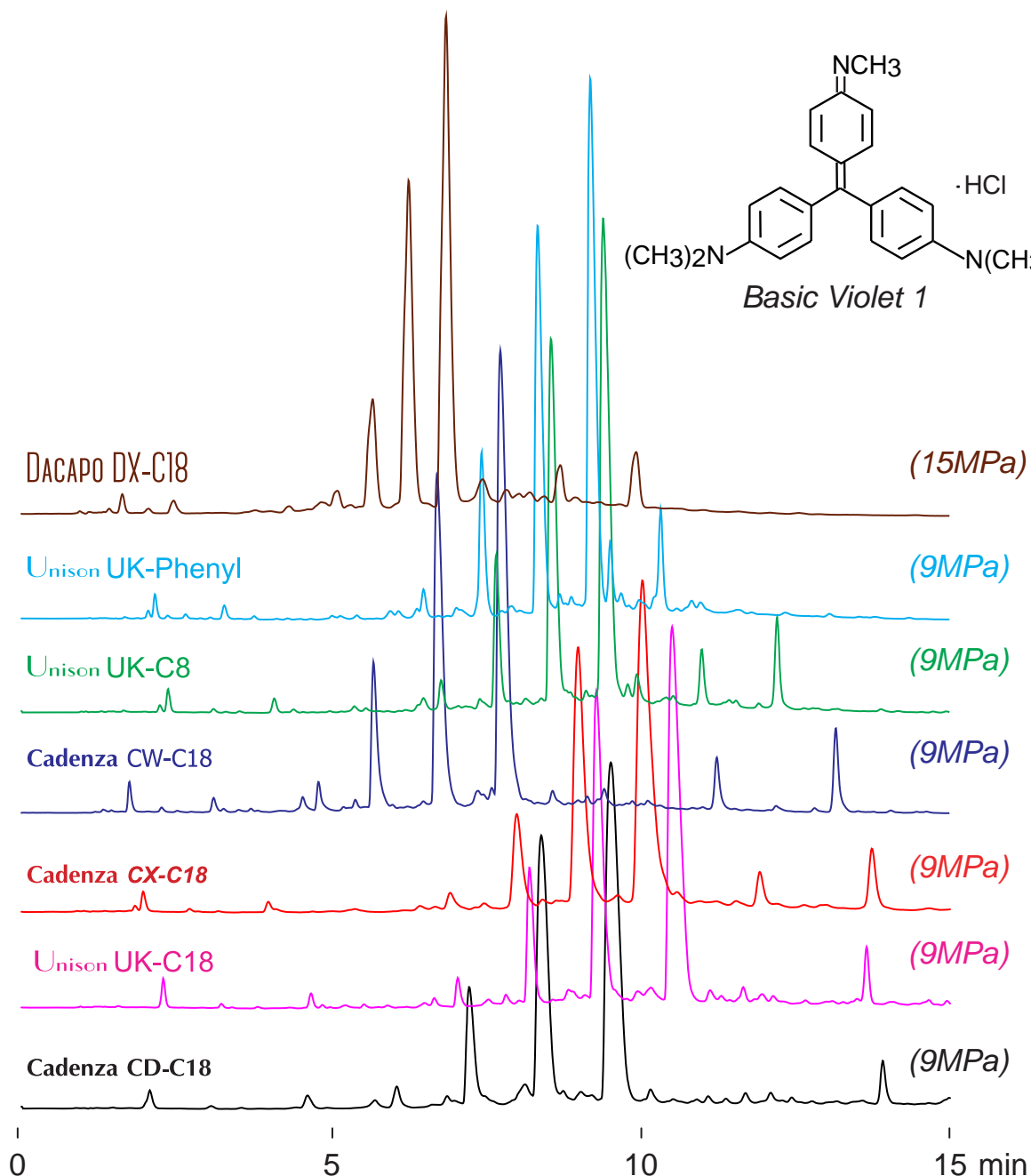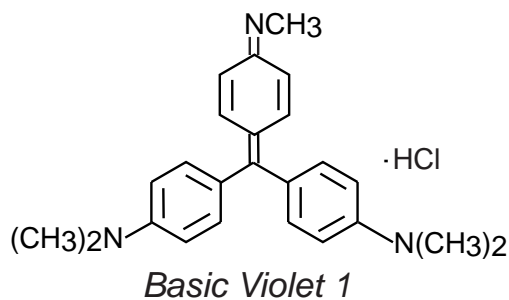

100x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

35-100%B (0-15min)

1 mL/min, 37°C, 260nm

10 μ L (1mg/mL, ACN / water / formic acid = 50 / 50 / 1)