


Intrada Organic Acid

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: Tyrosinemia type 1 related compounds

チロシン血症 型関連化合物の LC-MS 分析

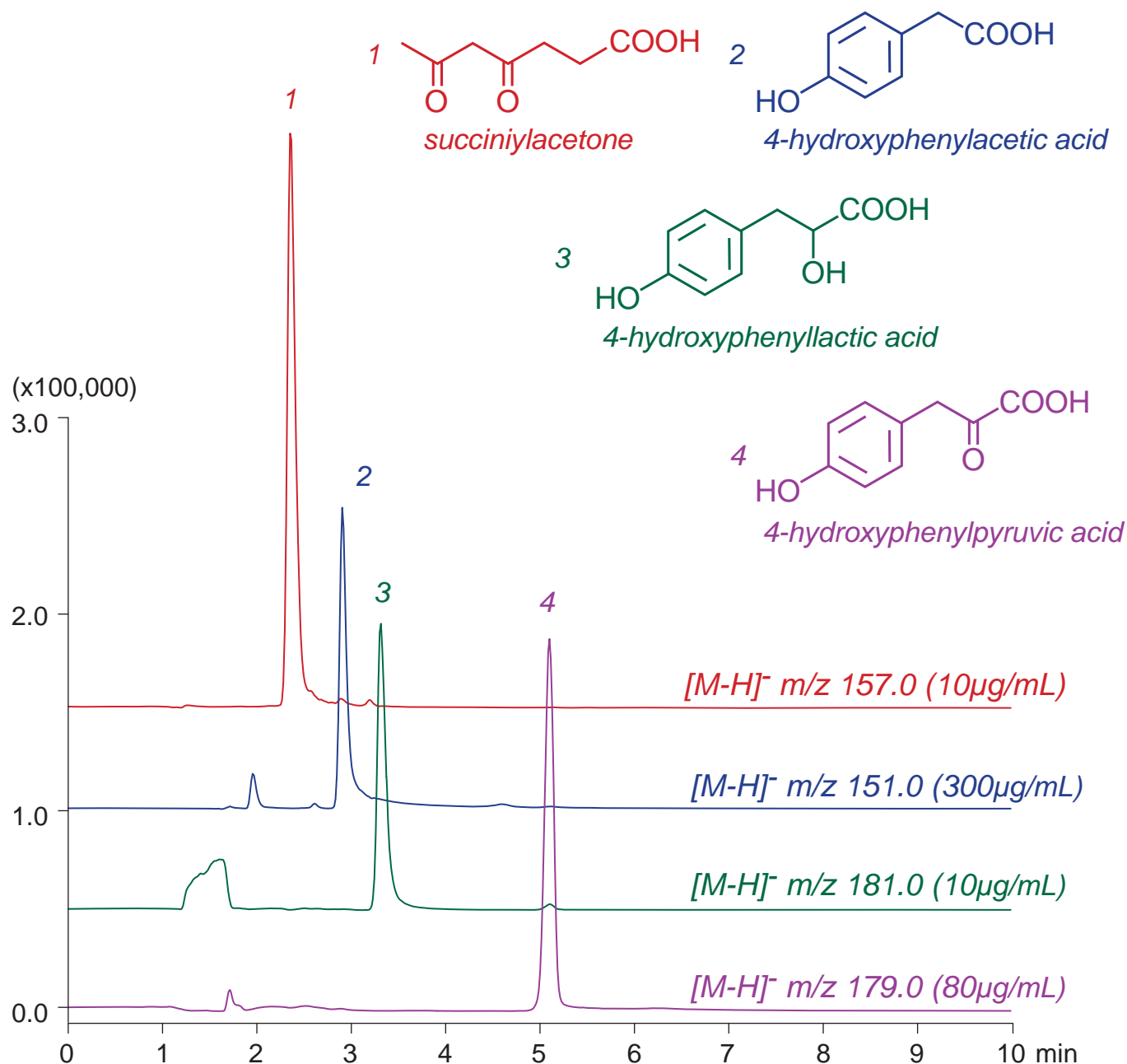

Intrada Organic Acid, 100 x 2 mm

A: acetonitrile / water / formic acid = 10 / 90 / 0.1

B: acetonitrile / 100mM ammonium formate = 10 / 90

0%B (0-1min), 0-100%B (1-7min), 100%B (7-10min)

0.2mL/min (12MPa), 37°C, MS (ESI, SIM, Negative)

5µL (10-300µg/mL, 50mM ammonium formate)