

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8
Cadenza CX-C18



Unison UK-C18
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Indigo Carmine (Food additive)

青色 2号 (食品添加物)

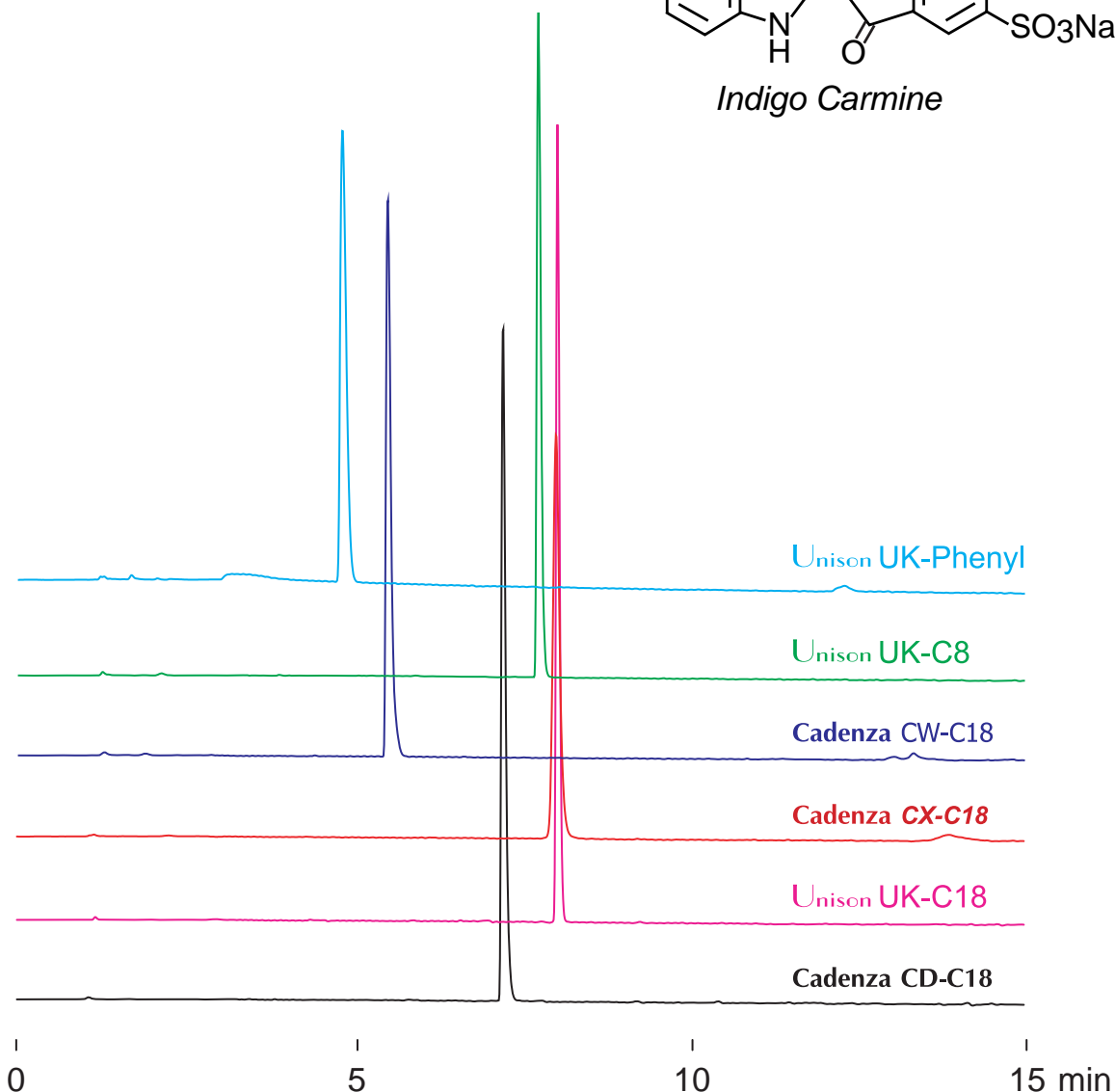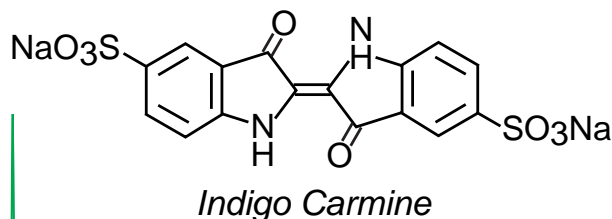

100x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-30%B (0-15min)

1 mL/min (9MPa), 37°C, 260nm

10 µL (1mg/mL, 100mM ammonium formate)