

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8  
Cadenza CX-C18


Unison UK-C18  
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Brilliant Blue ( Food additive )

青色 1号 ( 食品添加物 )

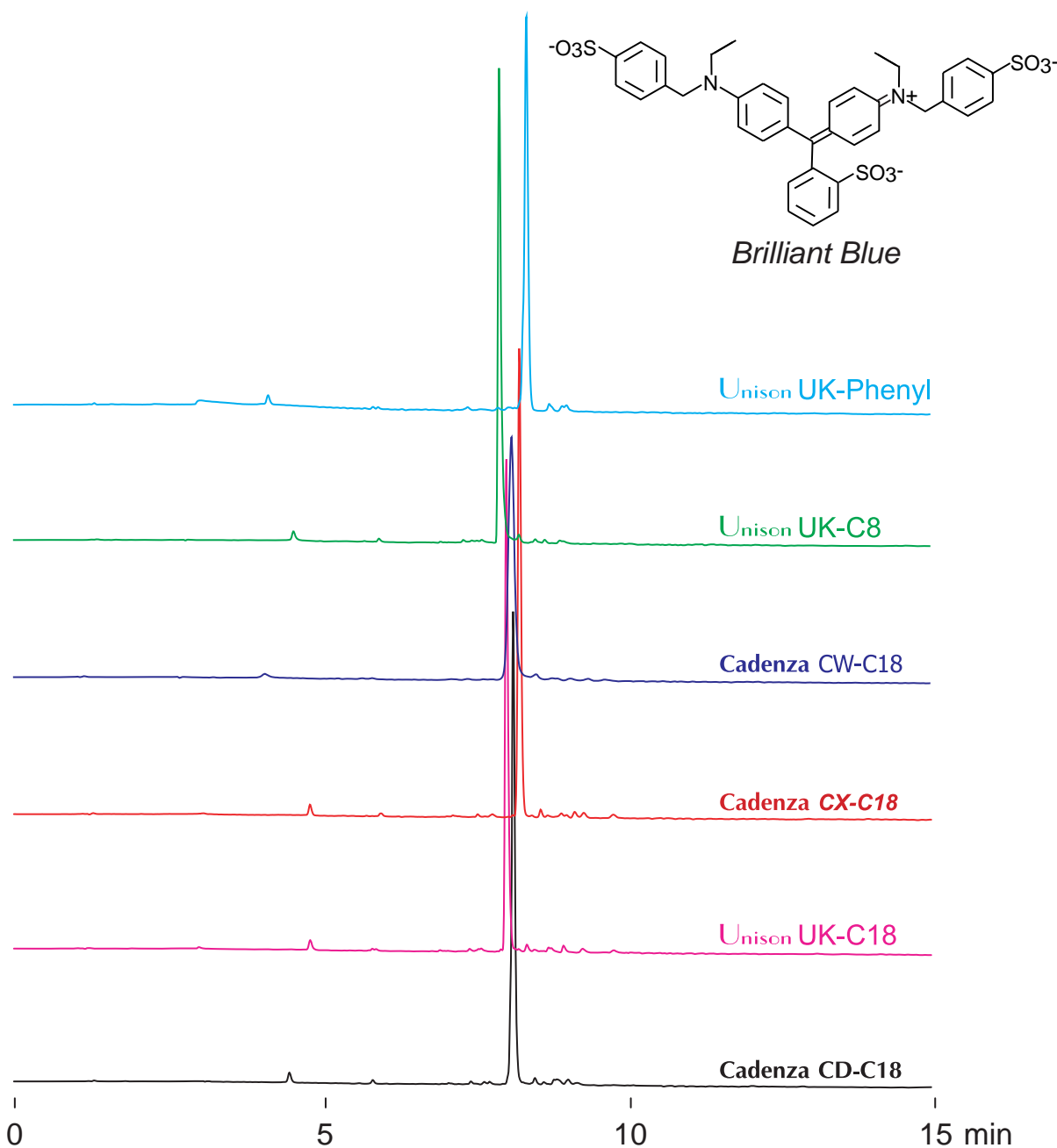

100x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-75%B (0-15min)

1 mL/min (9MPa), 37°C, 260nm

10 µL (0.5mg/mL, 100mM ammonium formate)