

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Acid Blue 3 ( Dye )

アシッドブルー 3 ( 染料 )

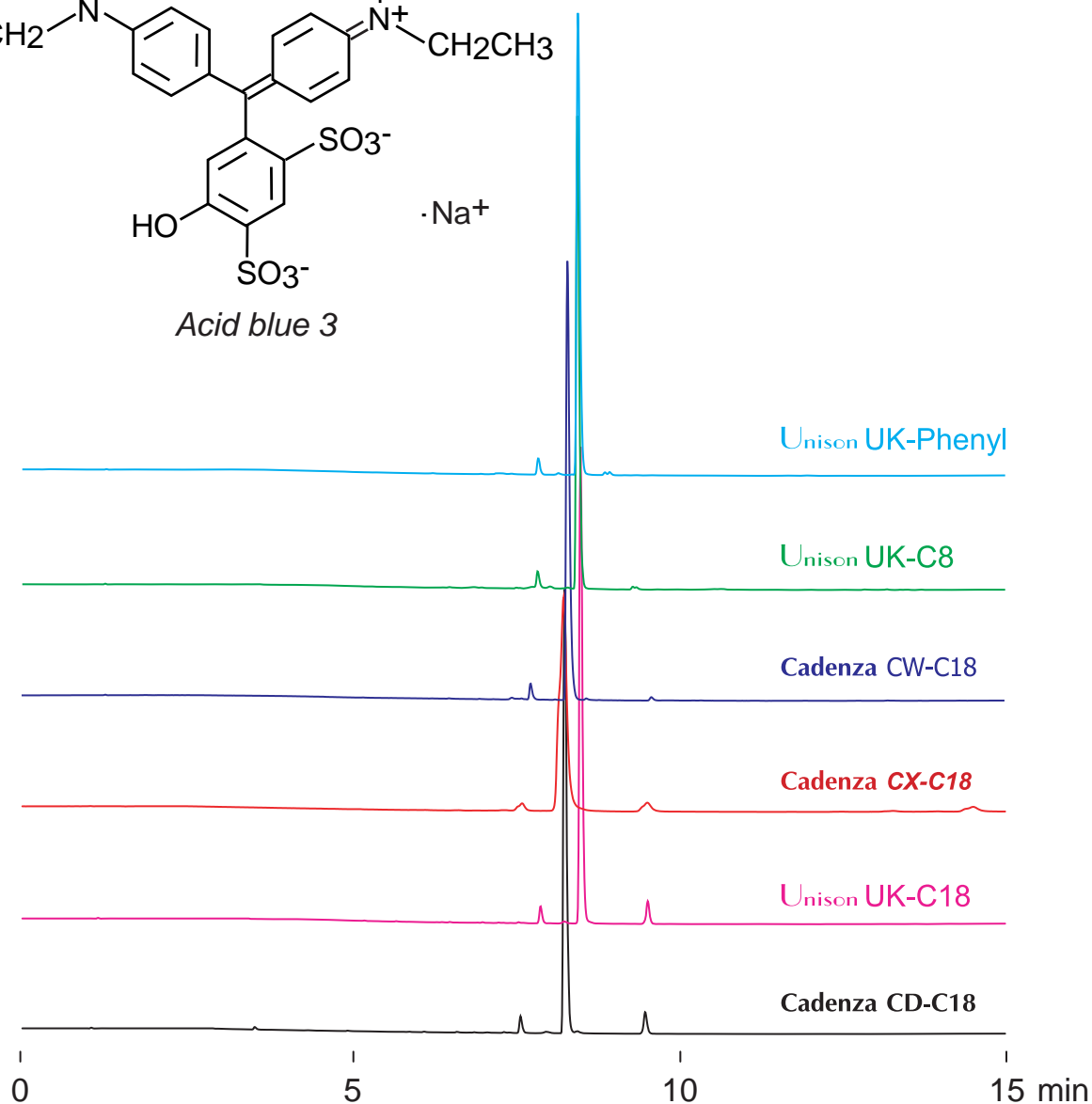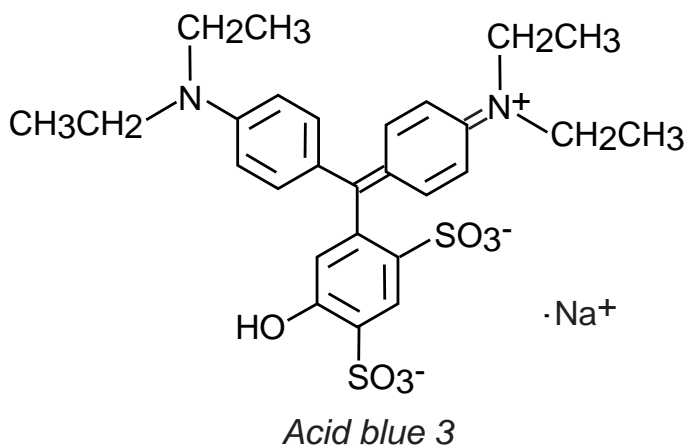

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium formate= 90 / 10

0-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, water / formic acid = 100 / 1)