

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8
Cadenza CX-C18



Unison UK-C18
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Tartrazine (Food additive)

黄色 4号 (食品添加物)

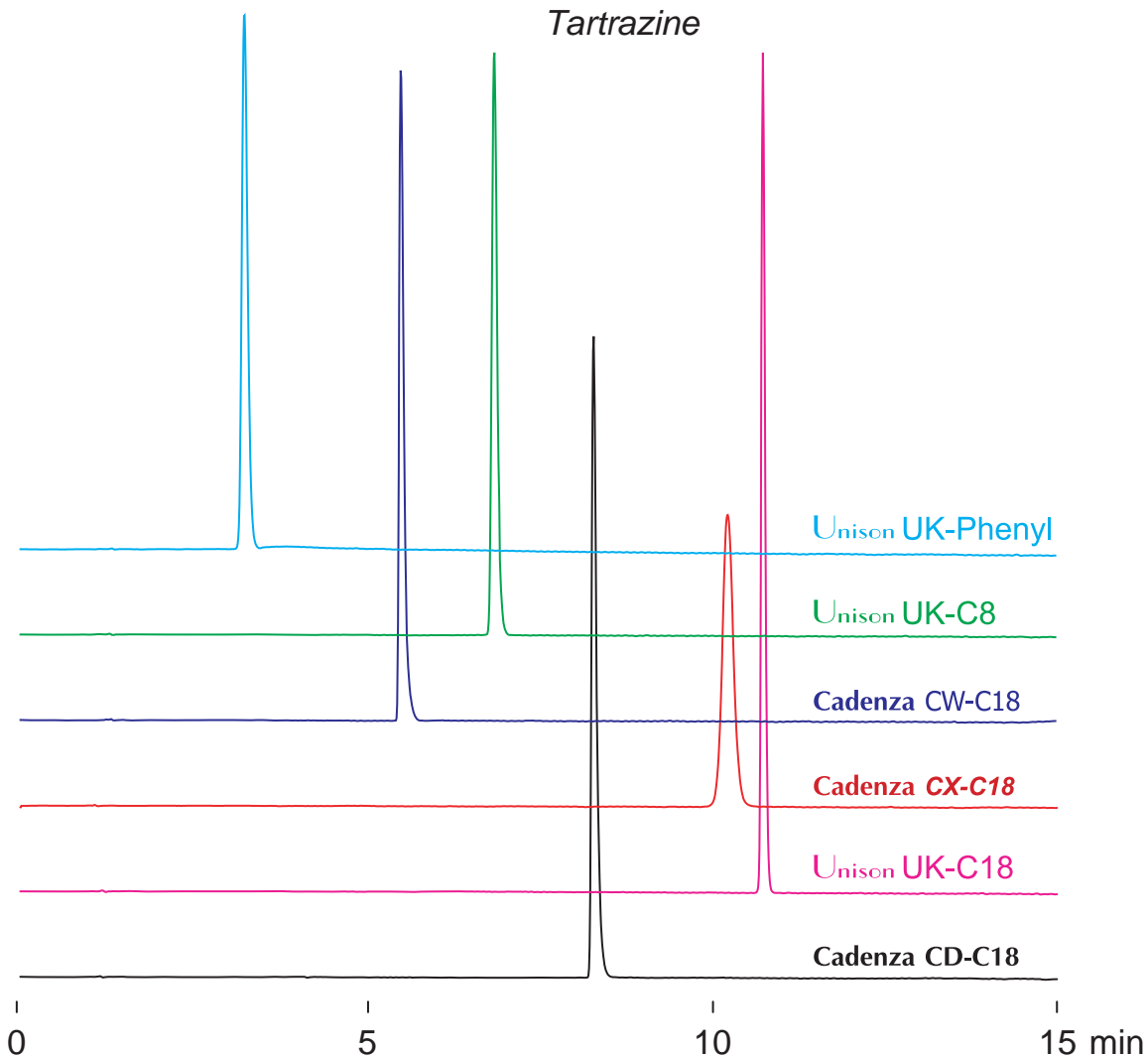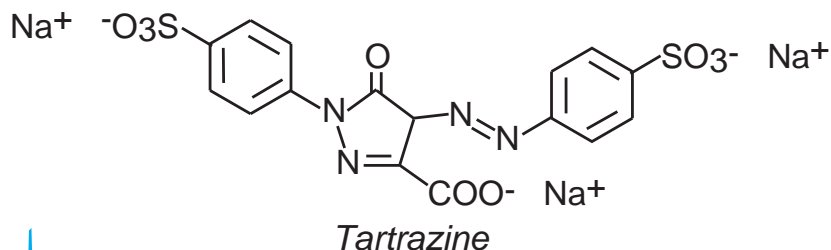

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-15%B (0-15min)

1mL/min (9MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate)