

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

Ponceau 3R (Artificial coloring)

ポンソー 3R (人工着色料)

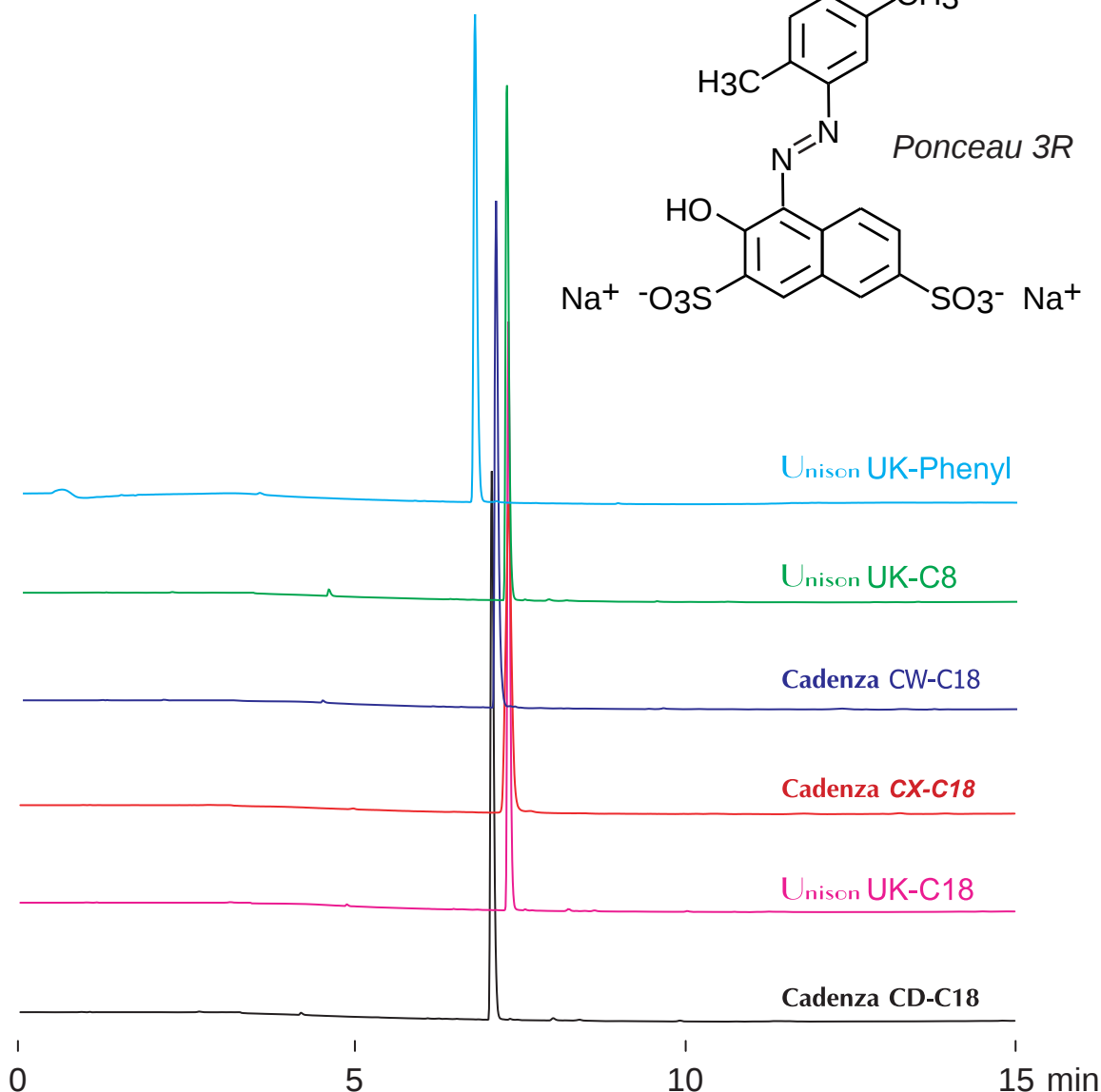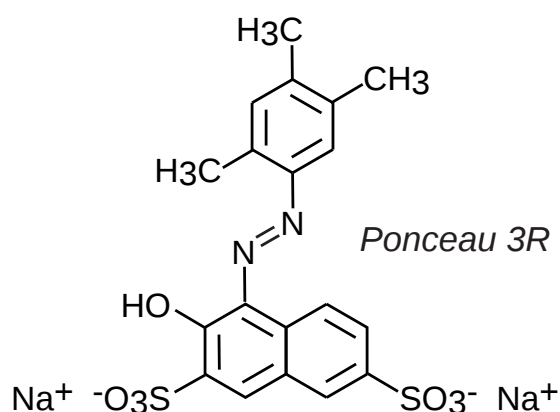

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

0-100%B (0-15min)

1mL/min (9MPa), 37°C, 260nm

10μL (0.25mg/mL, 100mM ammonium formate)