

Unison UK-Phenyl

Unison UK-C8

Unison UK-C18


Cadenza CW-C18

Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Phloxine B (Food additive)

赤色 104号 (食品添加物)

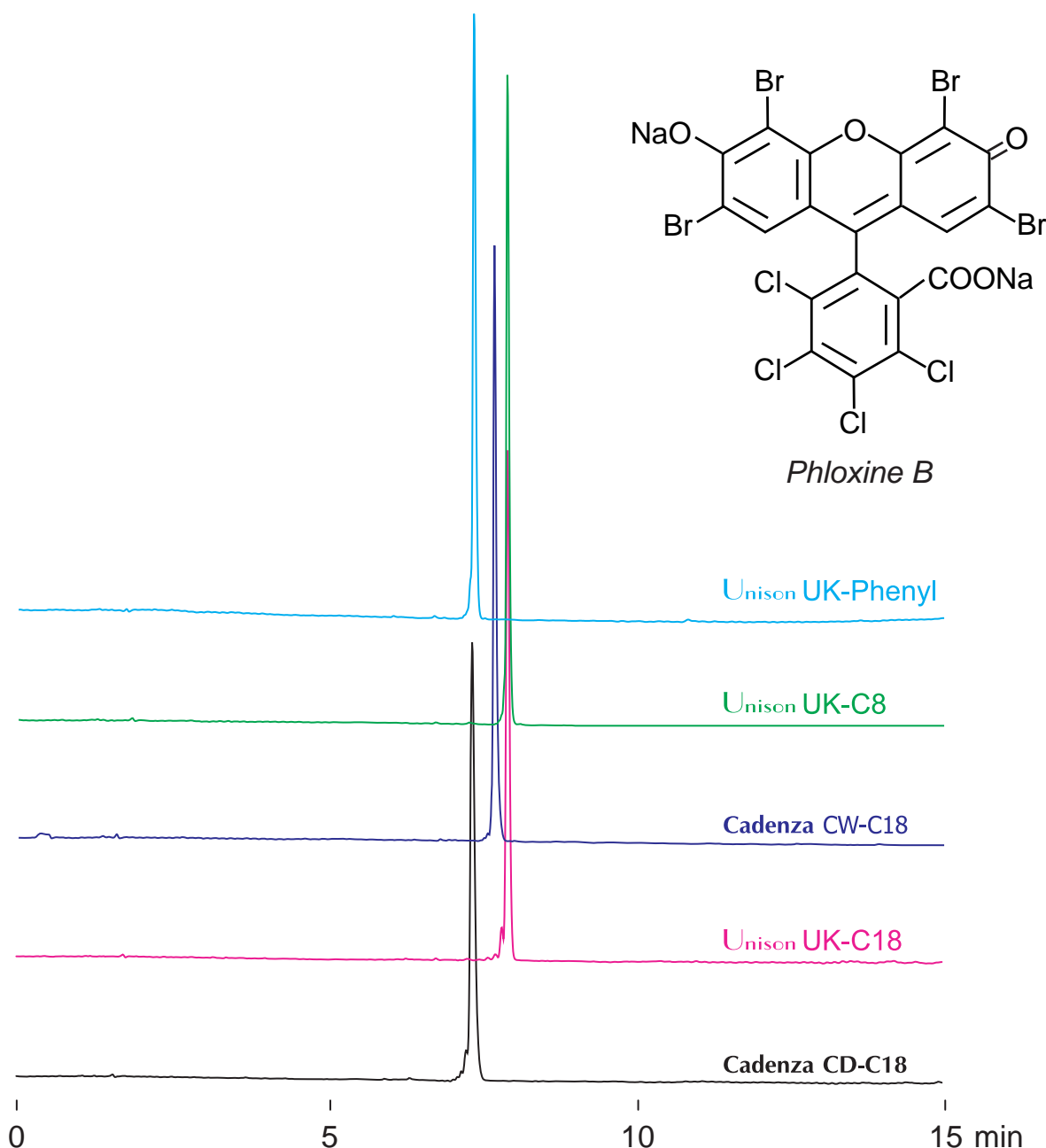

100 x 4.6 mm

A: 100mM ammonium acetate

B: acetonitrile / 100mM ammonium acetate = 90 / 10

15-100%B (0-15min)

1mL/min (8MPa), 37°C, 260nm

10µL (0.1mg/mL, 100mM ammonium formate)