


Unison UK-Phenyl Cadenza CW-C18  
 Unison UK-C8 Cadenza CX-C18  
 Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Light Green SF Yellowish ( Artificial coloring )

ライトグリーンSF黄口 (人工着色料)

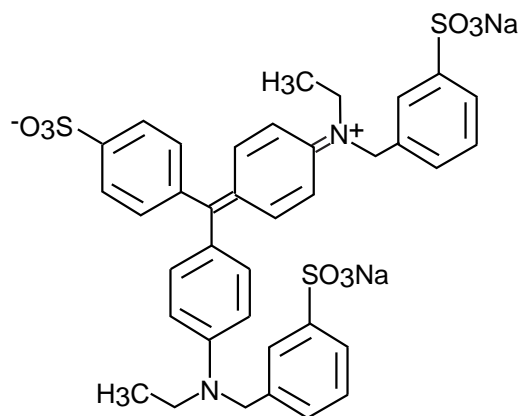

Light Green SF Yellowish

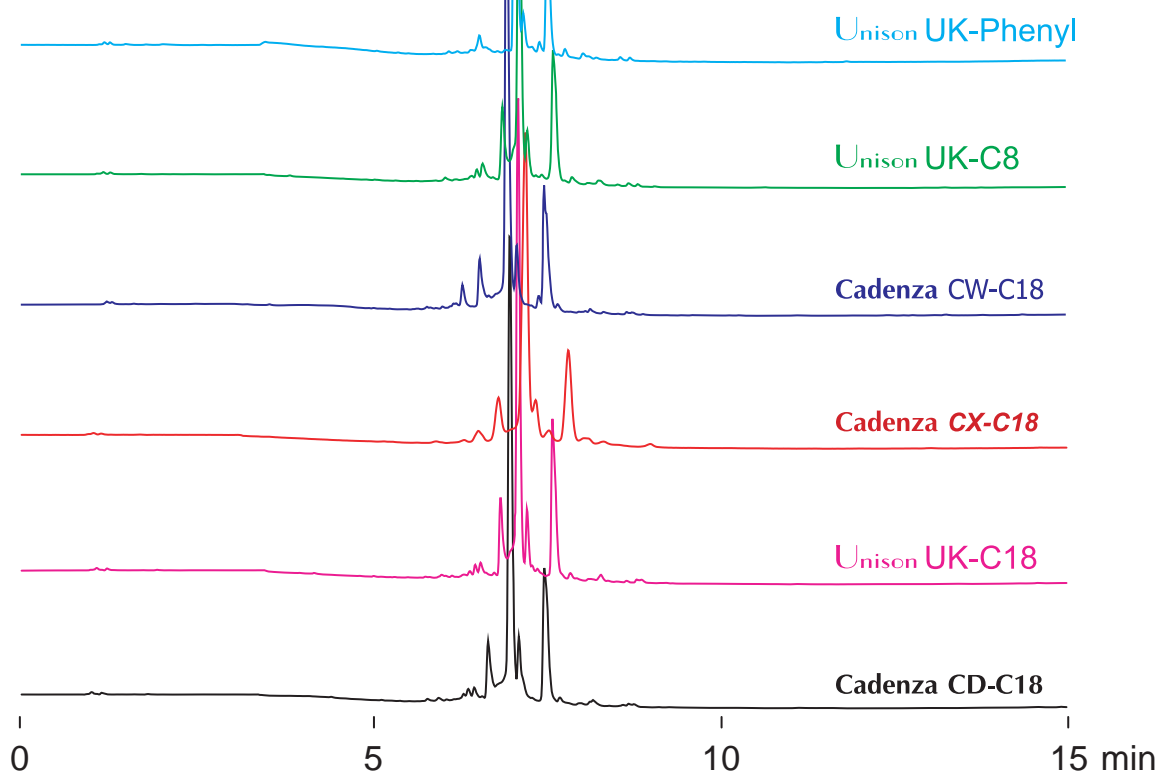

100 x 4.6 mm  
 A: 100mM ammonium acetate  
 B: acetonitrile / 100mM ammonium acetate = 90 / 10  
 0-100%B (0-15min)  
 1mL/min (9MPa), 37°C, 260nm  
 10µL (0.5mg/mL, 100mM ammonium formate)