
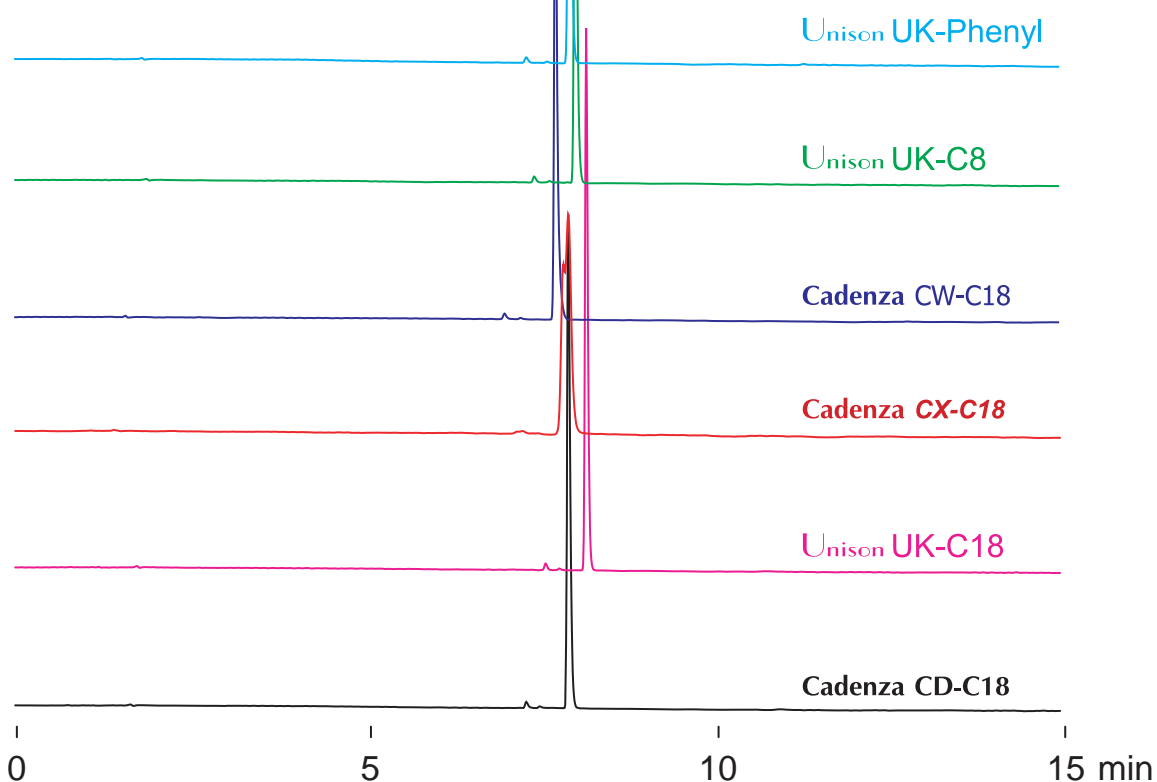


Unison UK-Phenyl Cadenza CW-C18  
 Unison UK-C8 Cadenza CX-C18  
 Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

### Acid Red ( Food additive )

赤色 106号 (食品添加物 )

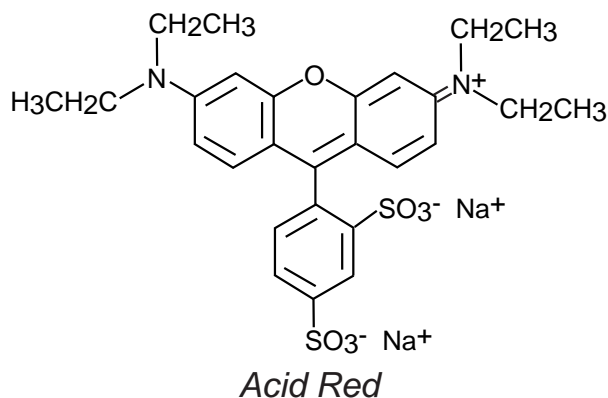

100 x 4.6 mm  
 A: 100mM ammonium acetate  
 B: acetonitrile / 100mM ammonium acetate = 90 / 10  
 10-100%B (0-15min)  
 1mL/min (8MPa), 37°C, 260nm  
 10µL (0.08mg/mL, 100mM ammonium formate)