



Scherzo SS-C18 MF

100 x 2 mm

Application

## LC-MS: N-nitrosoamines

N-ニトロアミンのLC MS分析

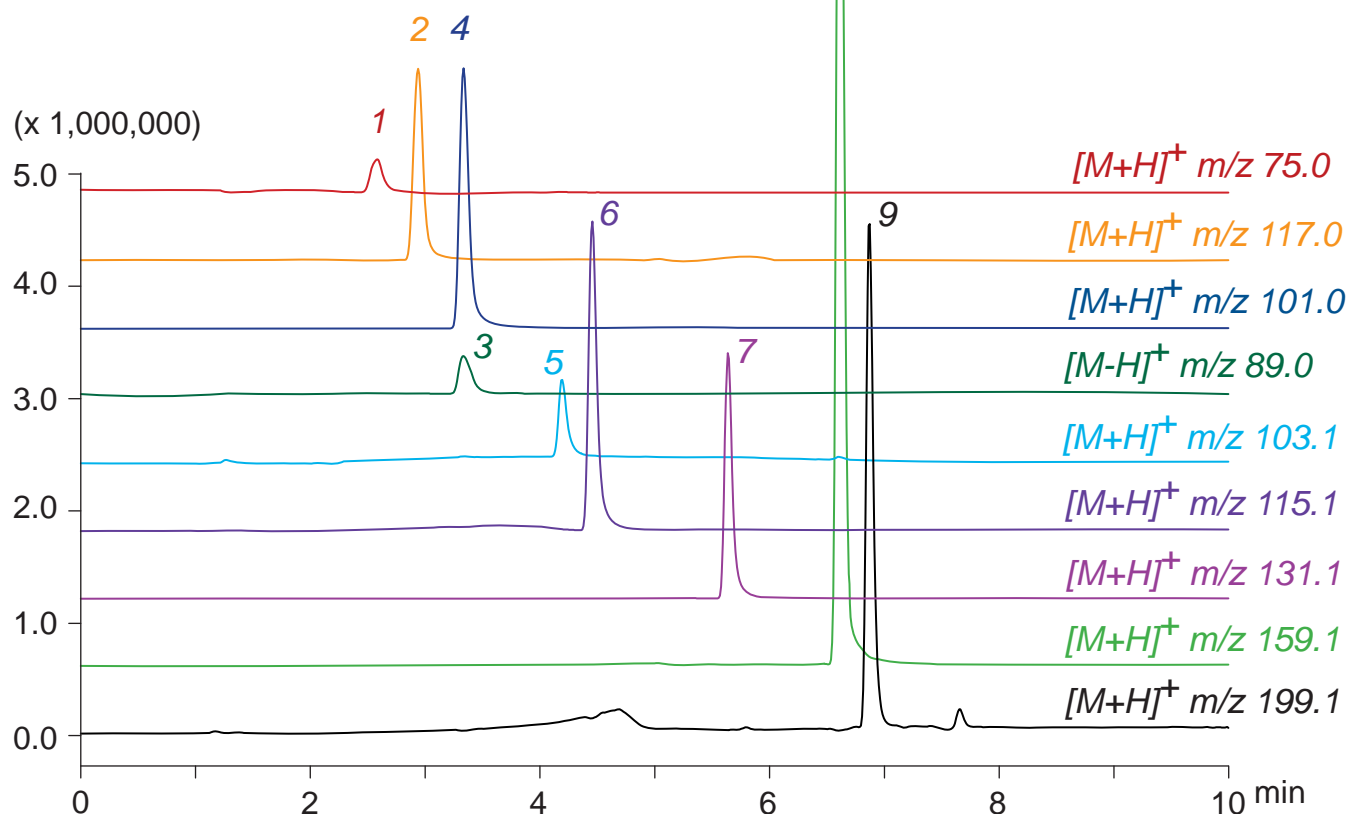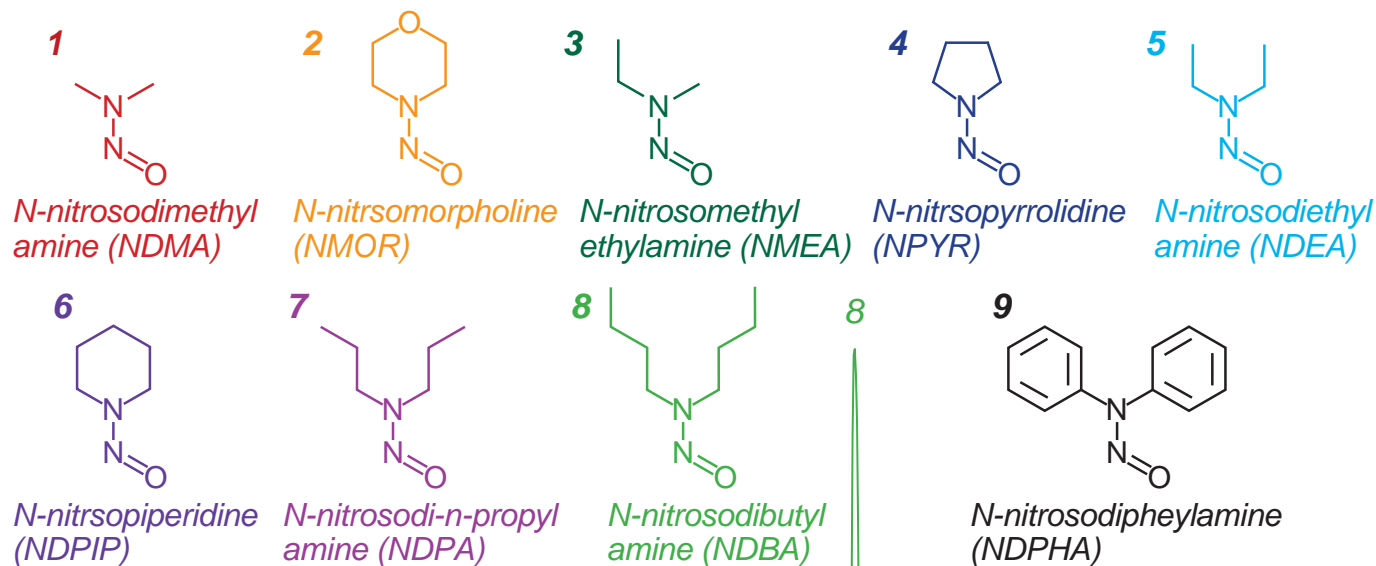

Scherzo SS-C18 MF, 100 x 2 mm  
 A: water / formic acid = 100 / 0.1  
 B: methanol / formic acid = 100 / 0.1  
 0-90%B (0-6min), 90%B (6-10min)  
 0.2mL/min (7MPa), 37°C, MS (ESI, SIM, positive)  
 5µL (4µg/mL, MeOH/water/HCOOH = 50/50/0.1)