
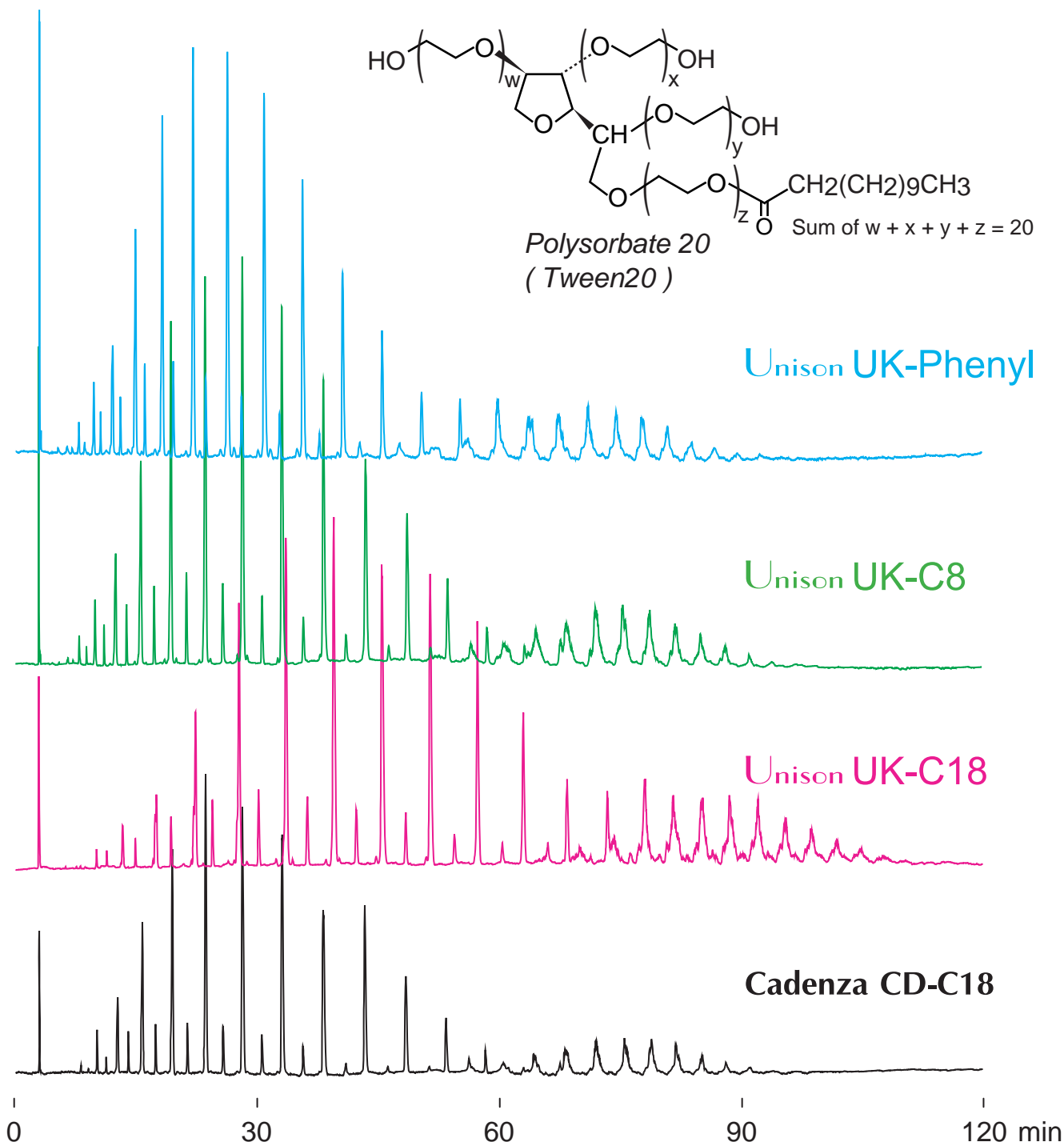


Unison UK-C18 Cadenza CD-C18 Unison UK-Pheny Unison UK-C8 250 x 4.6 mm

HPLC COLUMN APPLICATION

Polysorbate 20, (Surfactant)

ポリソルベート 20 (界面活性剤)

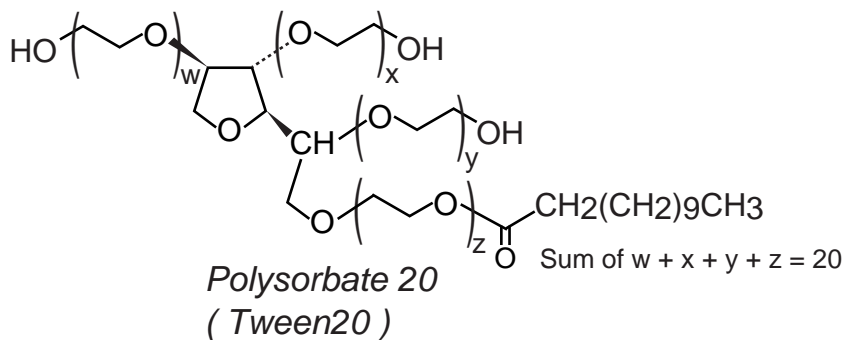

250 x 4.6 mm
 A: water / acetic acid = 100 / 0.1
 B: acetonitrile / acetic acid = 100 / 0.1
 13-30%B(0-120min), 1.0mL/min(14MPa,)
 60°C, ELSD(spray chamber 70°C, drift chamber 120°C)
 20µL (50mg/mL, 100mM ammonium formate)