

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8

Cadenza CX-C18

Unison UK-C18



Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Clopidogrel ( Anti platelet agents )

クロピドグレル (抗血小板薬)

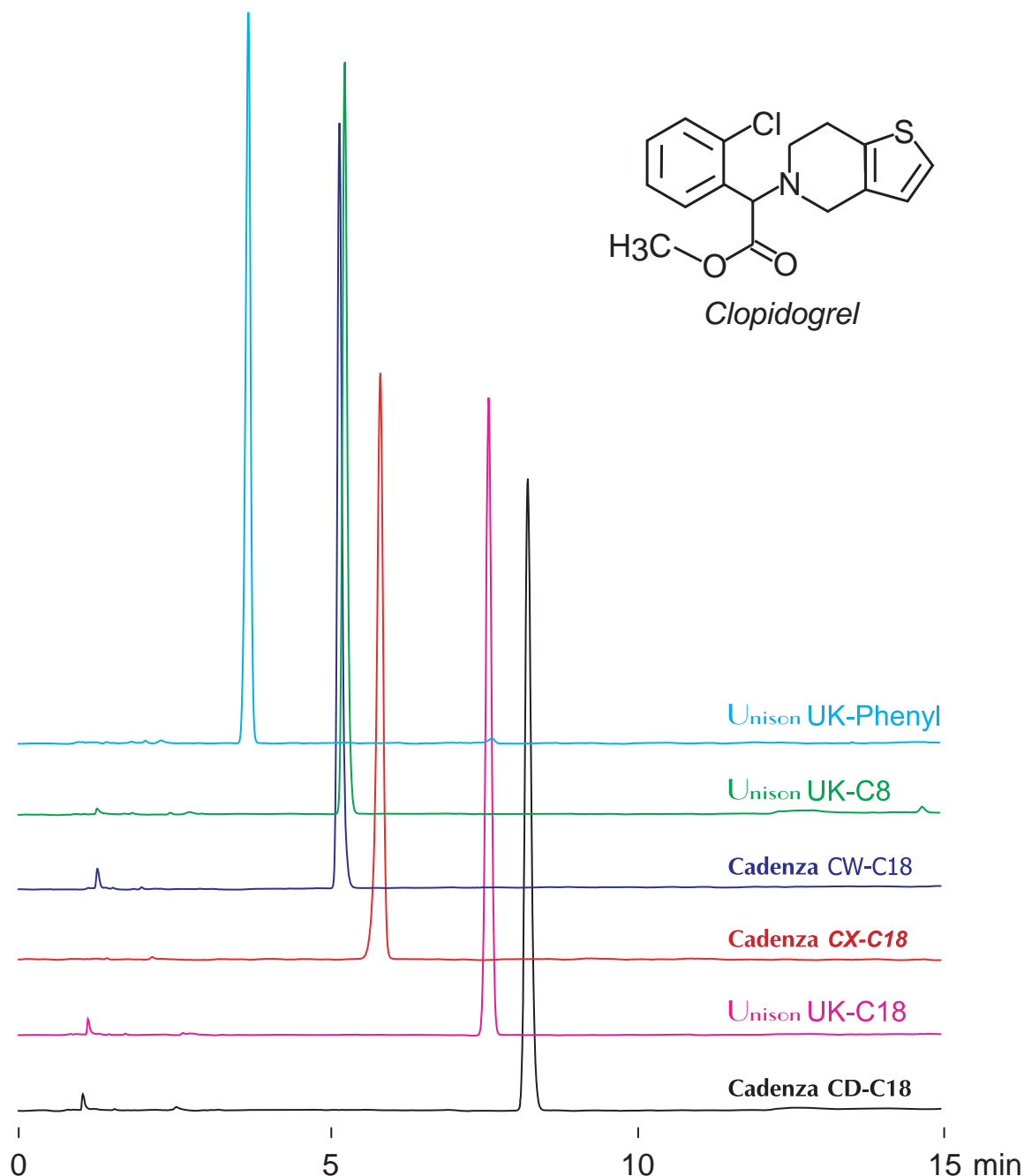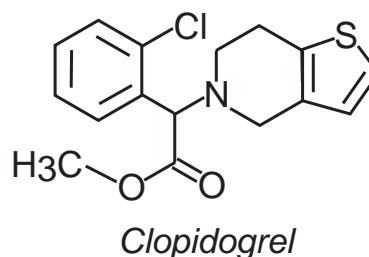

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate= 90 / 10

60-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (1mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)