




**Cadenza CD-C18**  
**Unison UK-C18**

250 x 4.6 mm

**HPLC COLUMN APPLICATION**

**Polysorbate 20 ( Nonionic surfactant )**

ポリソルベート20 (非イオン性界面活性剤 )

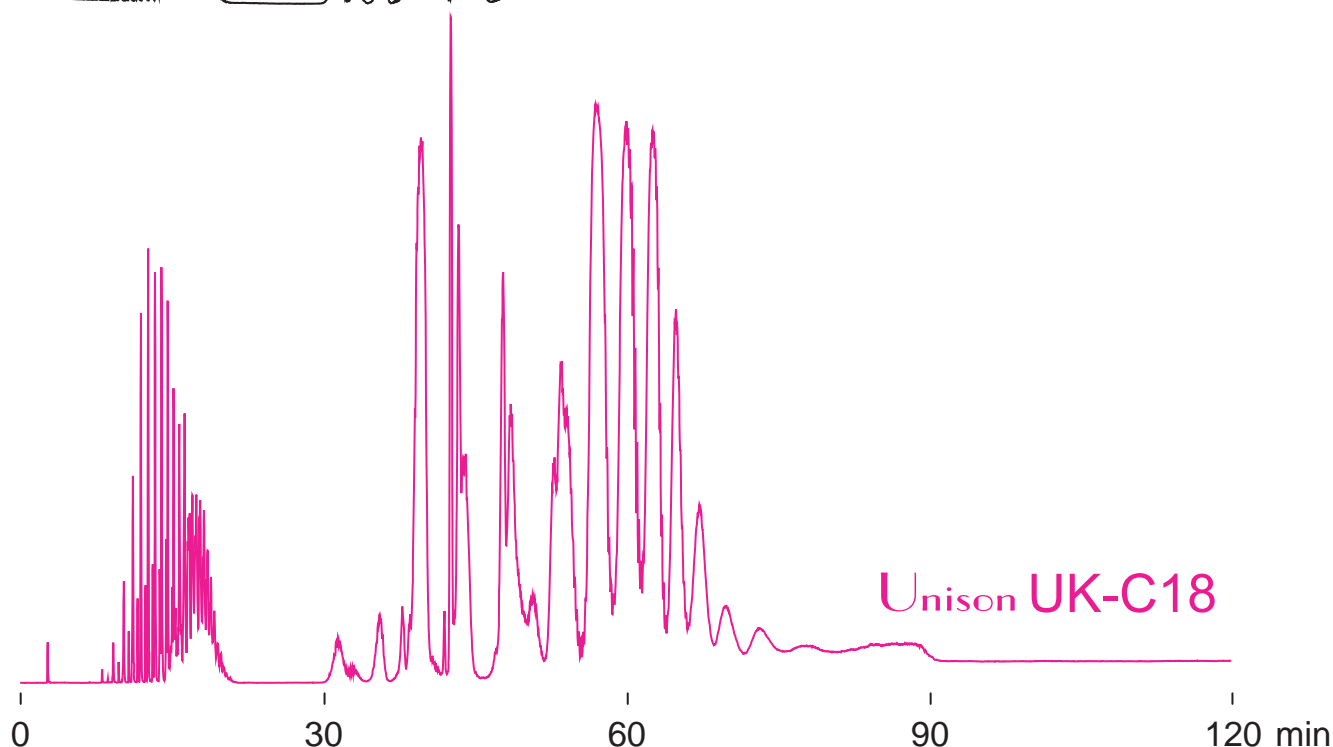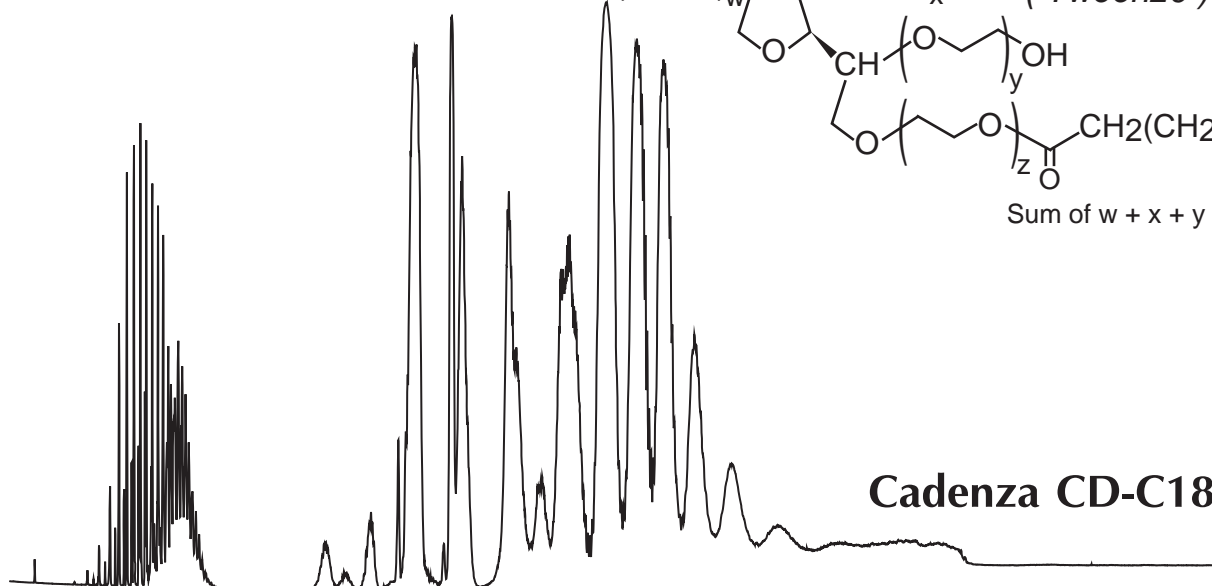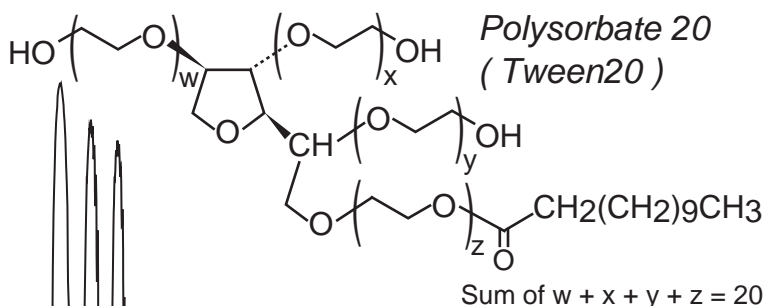

250 x 4.6 mm

A: water / acetic acid = 100 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

13-100%B(0-60min), 100%B(60-120min),

1.0mL/min(14MPa.), 60°C, ELSD(spray chamber 70°C, drift chamber 120°C)

20µL (50mg/mL, 100mM ammonium formate)