



Unison UK-C1
Unison UK-C18

250 x 4.6 mm


HPLC COLUMN APPLICATION

Polysorbate 20 (Nonionic surfactant)

ポリソルベート20 (非イオン性界面活性剤)

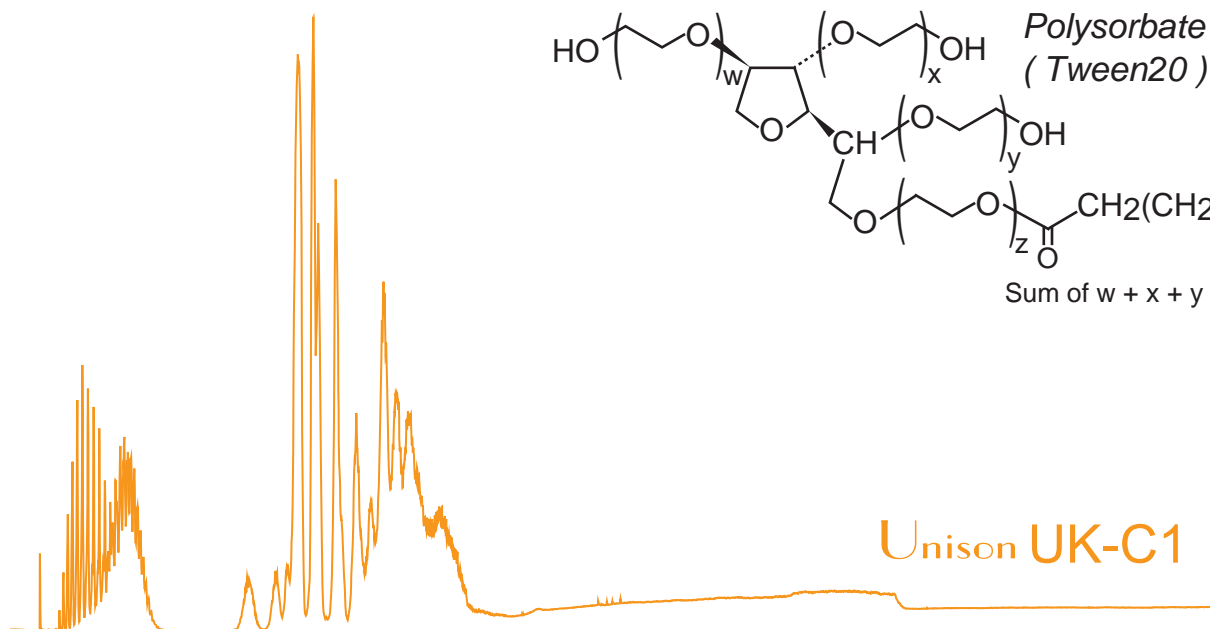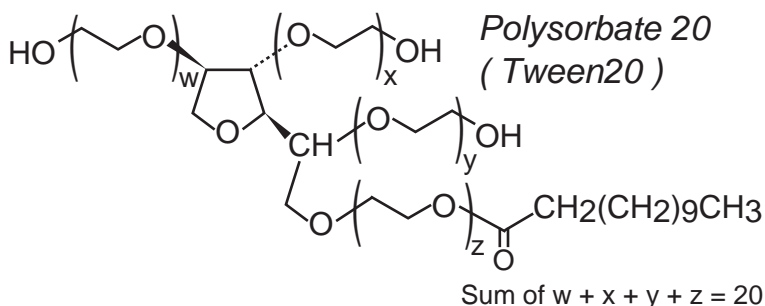

Unison UK-C1

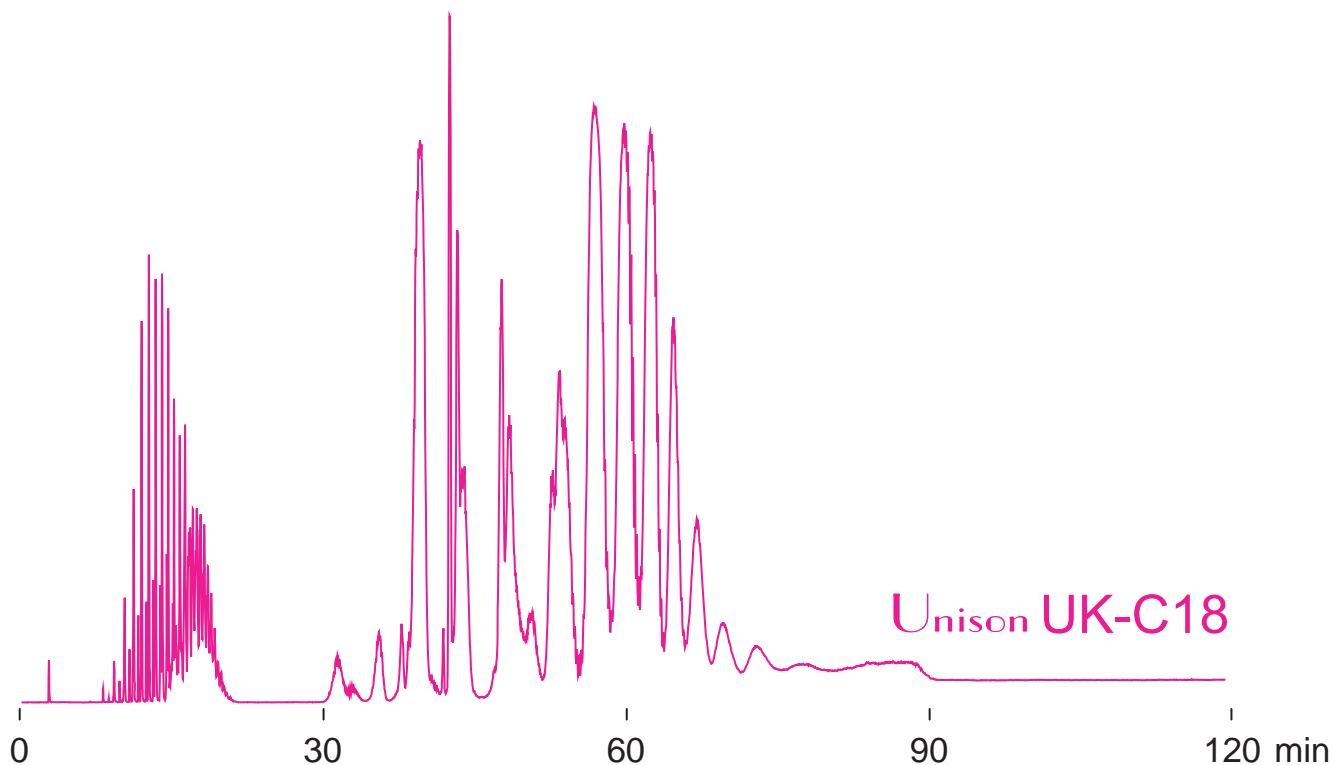

Unison UK-C18

250 x 4.6 mm

A: water / acetic acid = 100 / 0.1

B: acetonitrile / acetic acid = 100 / 0.1

13-100%B (0-60min), 100%B (60-120min),

1.0mL/min (14MPa.), 60°C, ELSD (spray chamber 70°C, drift tube 120°C)

20µL (50mg/mL, 100mM ammonium formate)