




Intrada Organic Acid

150 x 2 mm

HPLC COLUMN APPLICATION

LC-MS: foseetyl, Ethephon

ホセチル, エテホンの LC-MS 分析

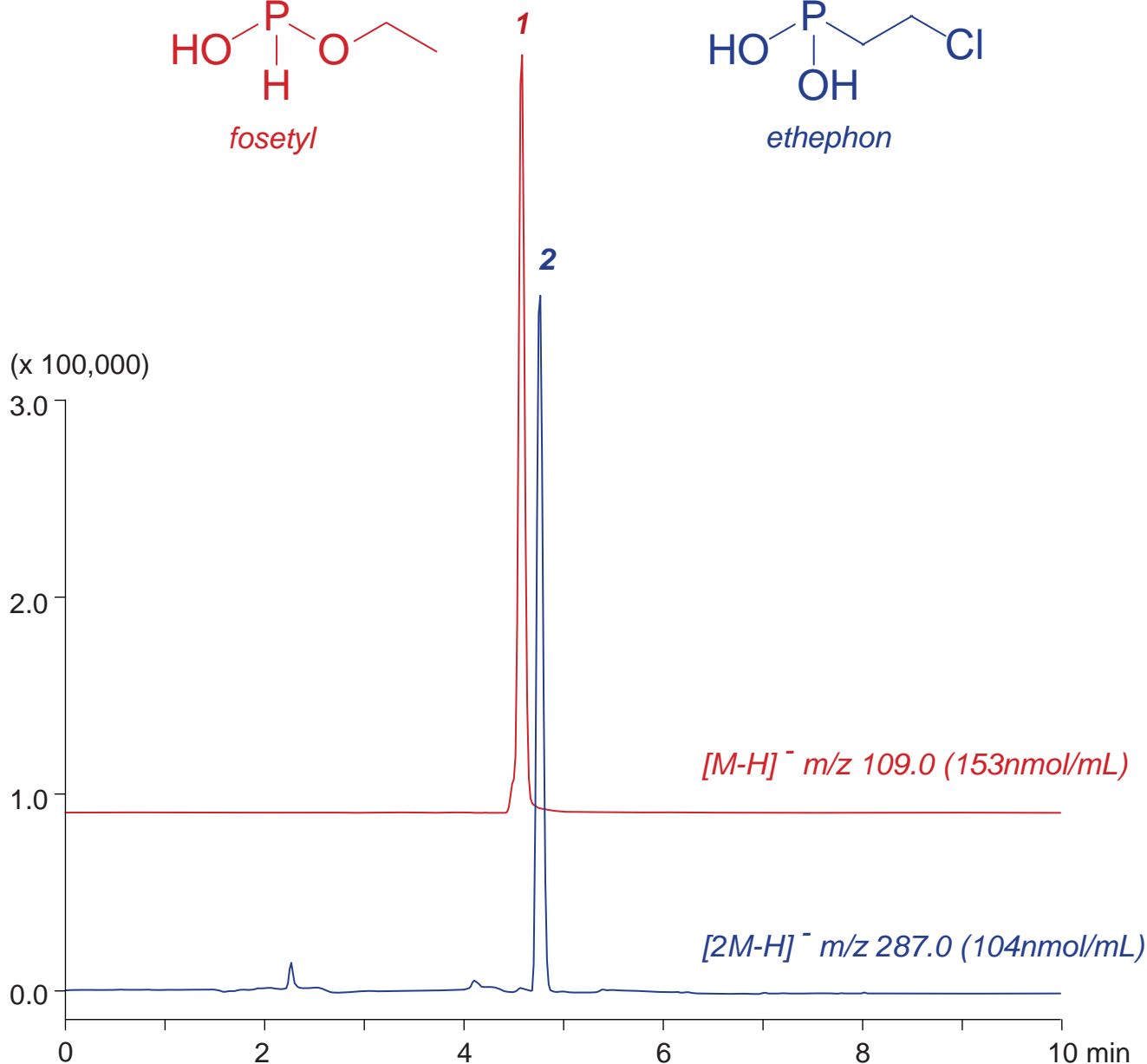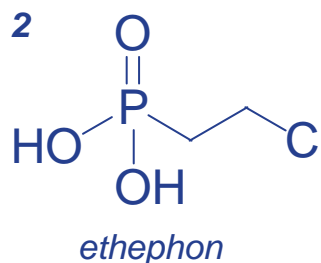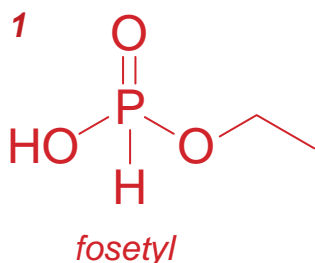

Intrada Organic Acid, 150 x 2 mm
 A: acetonitrile / 100mM ammonium formate = 95 / 5
 B: acetonitrile / 100mM ammonium formate = 10 / 90
 0%B (0-1min), 0-100%B (1-5min), 100%B (5-10min)
 0.2mL/min (5MPa), 37°C, Single Quad. MS (ESI)
 5µL (104-153nmol/mL, water / HCOOH = 100 / 2)