



Scherzo SM-C18

150 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Triazine herbicides

トリアジン系除草剤の LC-MS 分析

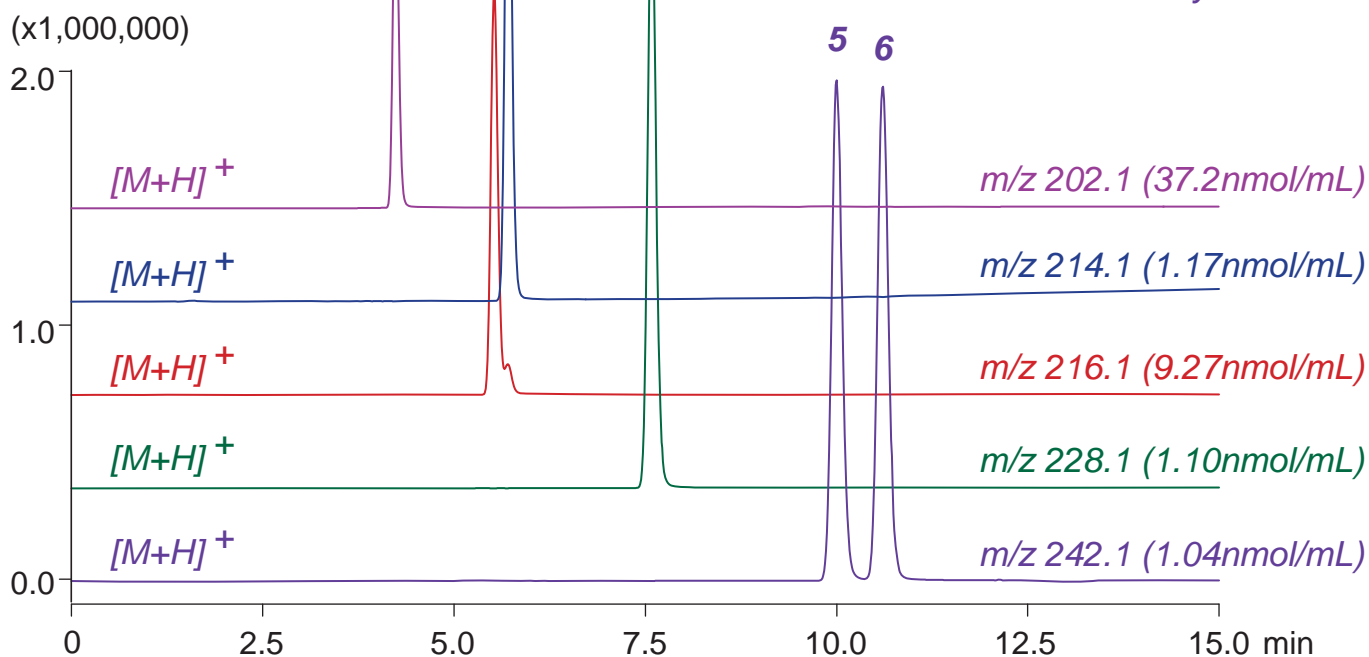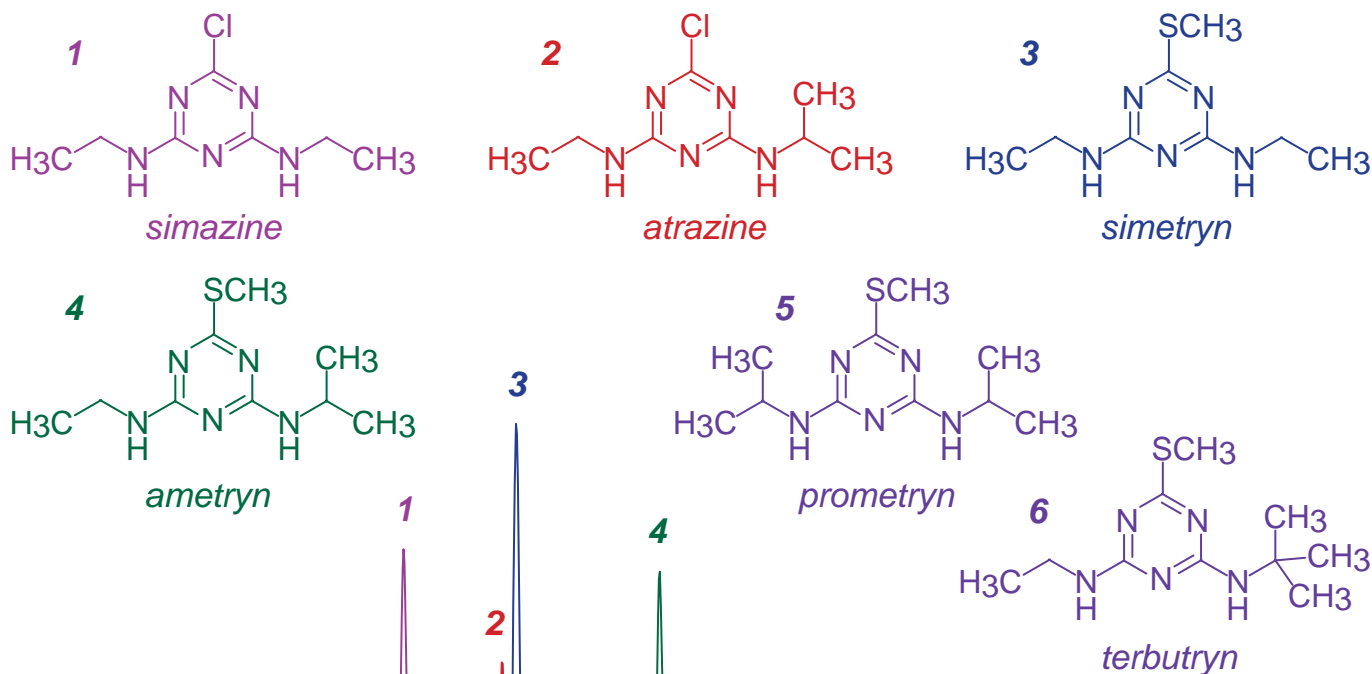

Scherzo SM-C18, 150 x 3 mm

A: acetonitrile / 20mM ammonium formate = 50 / 50

B: acetonitrile / 100mM ammonium formate = 70 / 30

0-70%B (0-15min)

0.4mL/min (8MPa), 37°C, Single Quad. MS (ESI)

5µL (1.04-37.2nmol/mL, MeOH / 100mM HCOONH<sub>4</sub>=60/40)