

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8  
Cadenza CX-C18



Unison UK-C18  
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Warfarin ( Anticoagulant )

ワルファリン ( 抗凝固薬 )

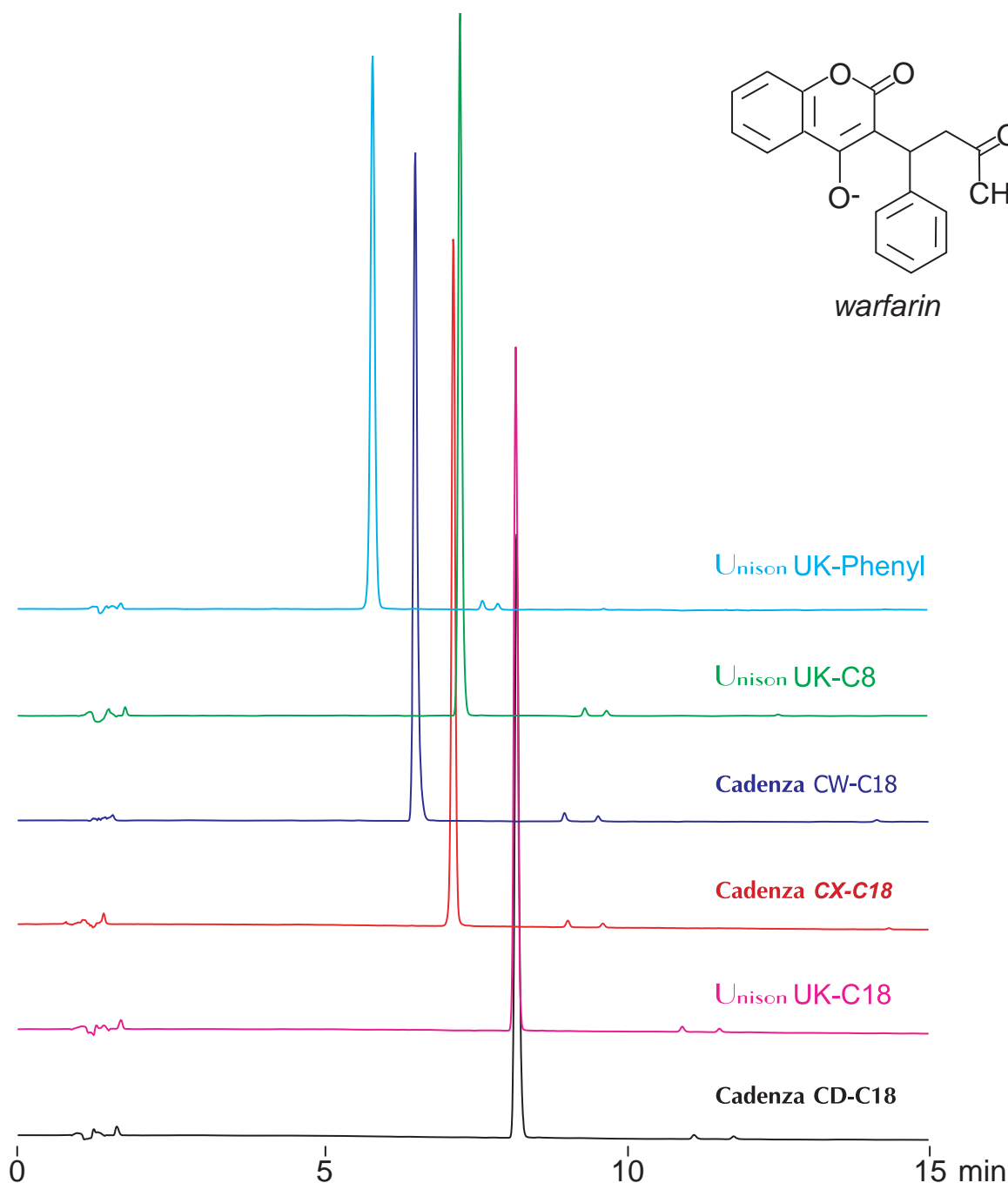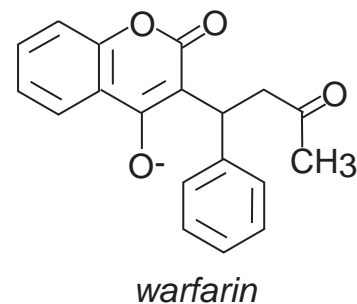

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

35-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.25mg/mL, methanol / water / HCOOH = 70 / 30 / 1)