

Unison UK-Phenyl Cadenza CW-C18

Unison UK-C8  
Cadenza CX-C18



Unison UK-C18  
Cadenza CD-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Warfarin ( Anticoagulant )

ワルファリン ( 抗凝固薬 )

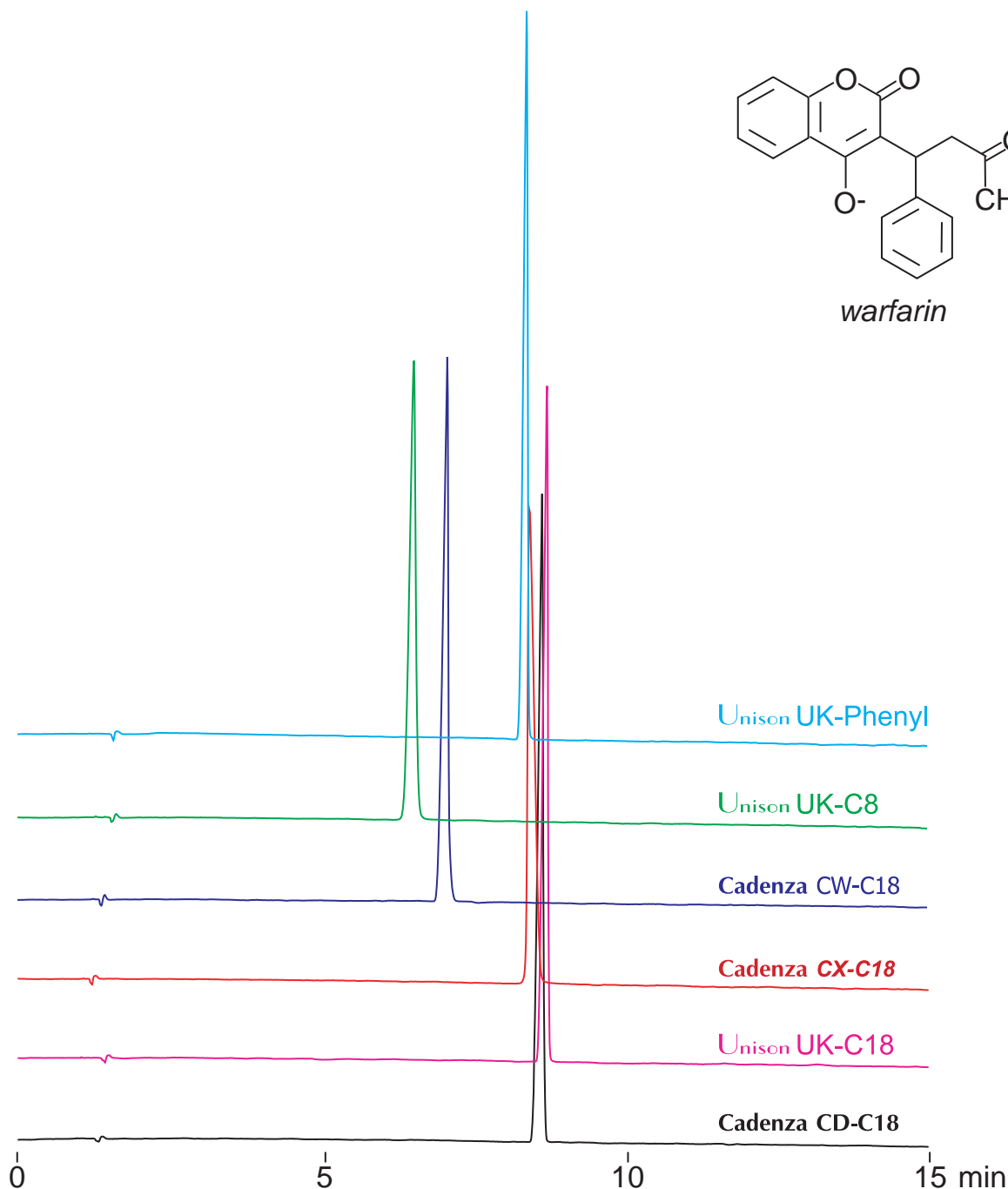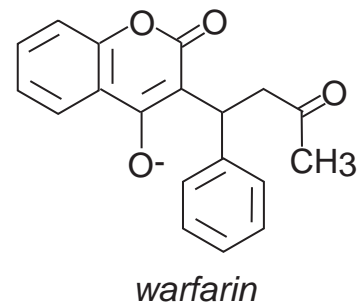

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

20-100%B (0-15min)

1.0mL/min (12MPa), 37°C, 260nm

10µL (0.25mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)