


**Cadenza CD-C18**

150 x 3 mm

**HPLC COLUMN APPLICATION**

## LC-MS: Organophosphorus pesticides ( 3 )

有機リン系殺虫剤 ( 3 ) の LC-MS 分析

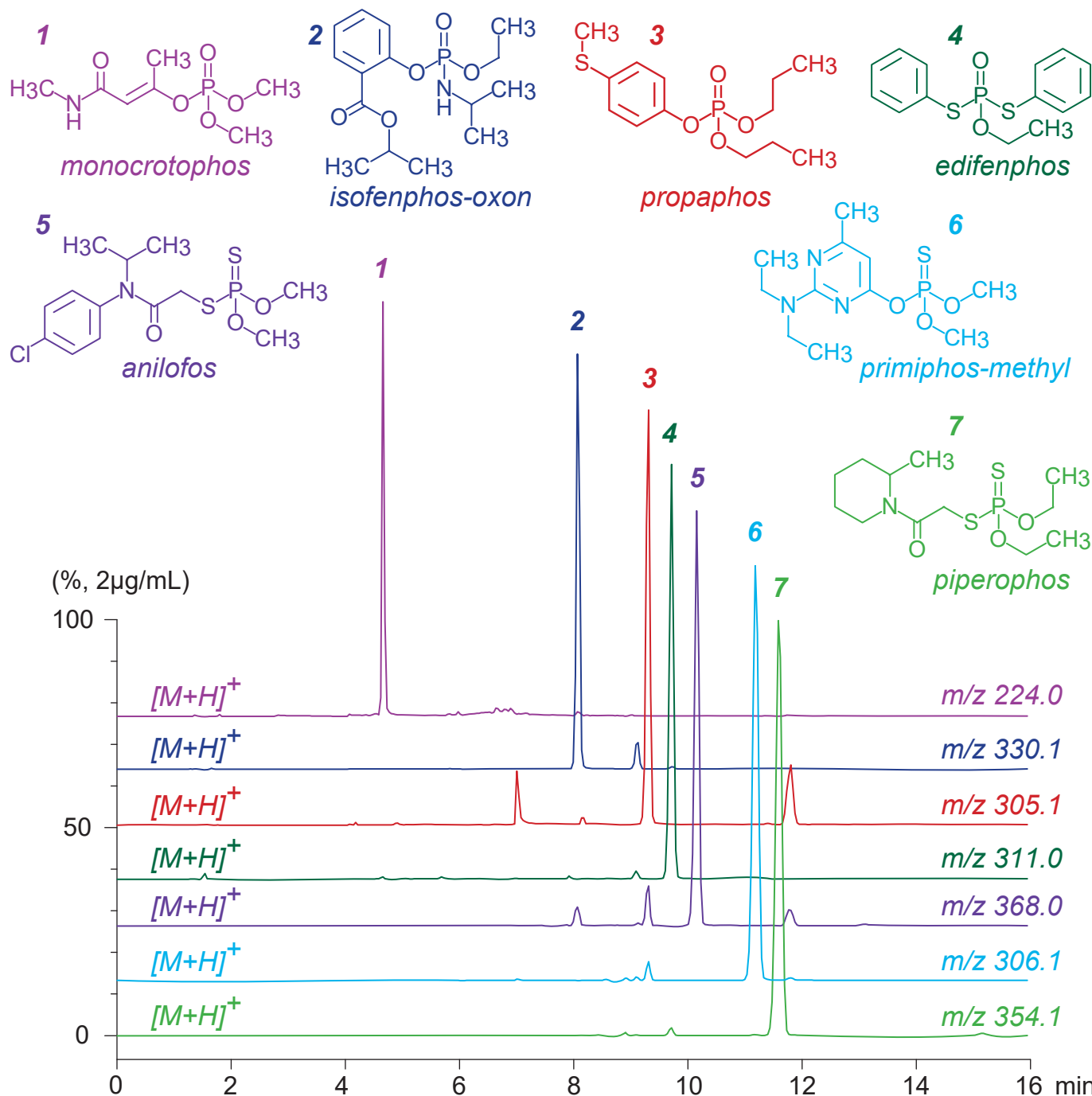

Cadenza CD-C18, 150 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

0-75%B (0-5min), 75%B (5-16min)

0.4mL/min (9MPa), 37°C, Single Quad. MS (ESI)

5µL (2µg/mL, ACN / water / HCOOH = 50 / 50 / 0.1)