

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

Coumatetralyl (Coumarin anticoagulant)

クマテトラリル (クマリン系抗血液凝固薬)

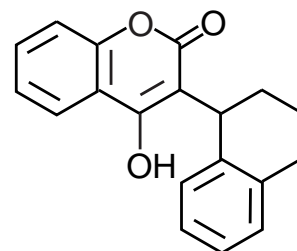

coumatetralyl

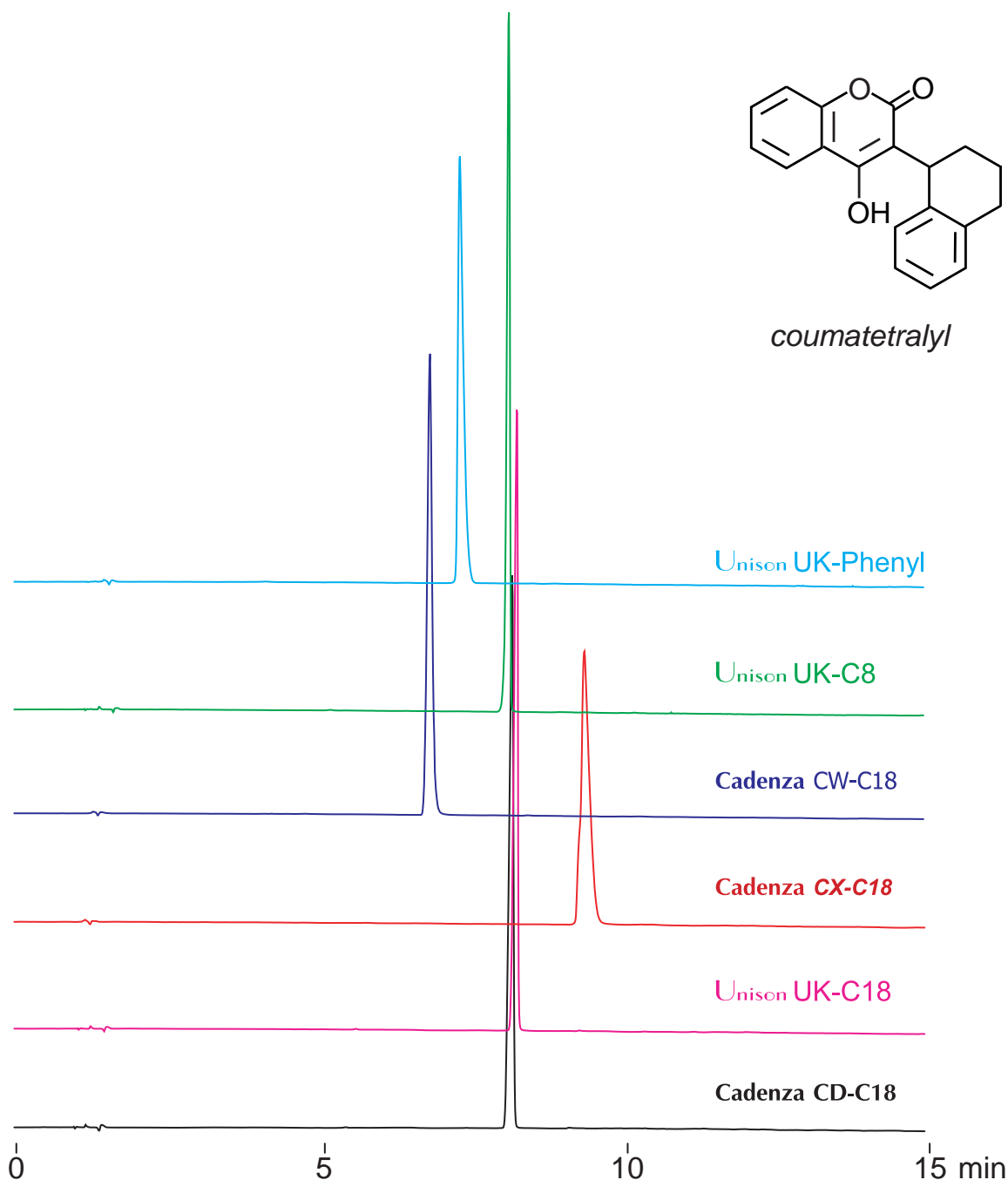

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

10-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 280nm

10µL (0.2mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)