


Scherzo SM-C18

100 x 4.6 mm


HPLC COLUMN APPLICATION

Trimetazidine (Vasodilator)


トリメタジジン (血管拡張薬)

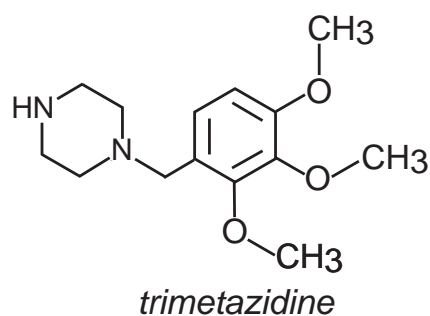

A: 10mM ammonium formate
 B: acetonitrile / water / formic acid
 = 70 / 30 / 0.1

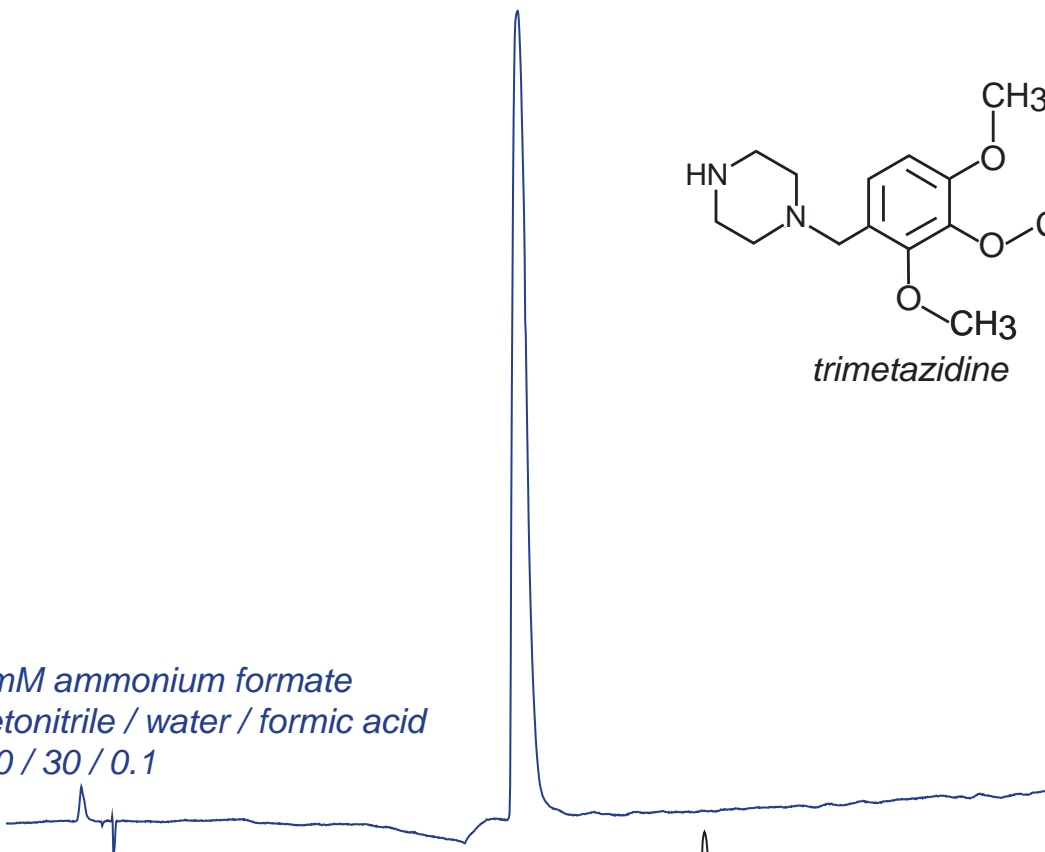

A: 5mM ammonium formate
 B: 100mM ammonium formate / acetonitrile
 = 30 / 70

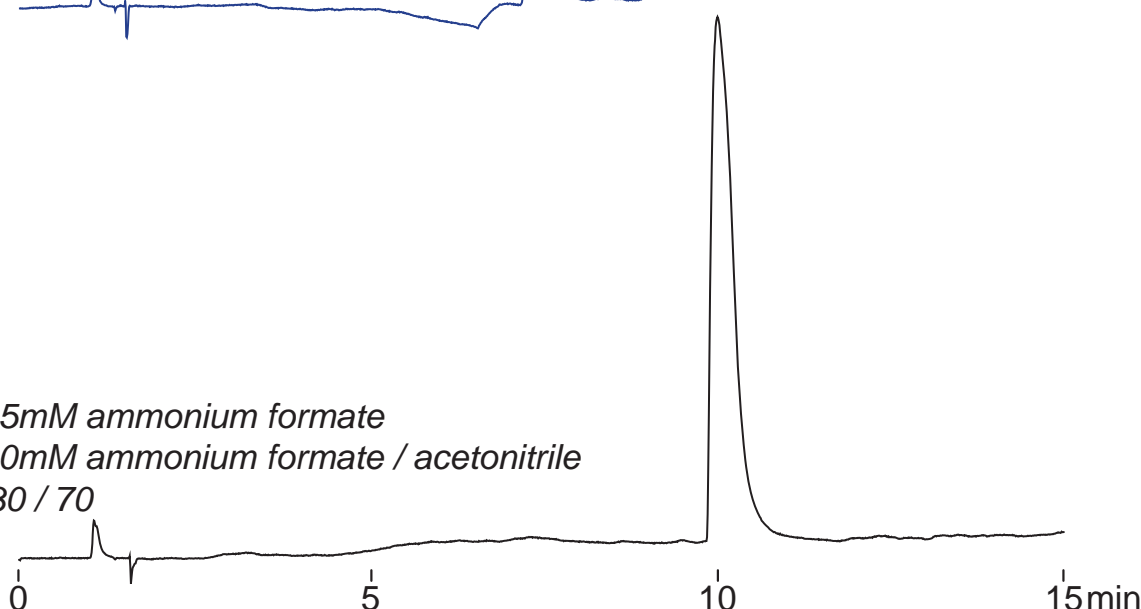

Scherzo SM-C18, 100 x 4.6 mm

0-100%B (0-15min)

1.0mL/min (10MPa), 37°C, 270nm

10µL (1.0mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)