

Unison UK-Phenyl Cadenza CW-C18


Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm


HPLC COLUMN APPLICATION

## Olmesartan medoxomil ( Antihypertensive )

オルメサルタン メドキシミル ( 血圧降下薬 )

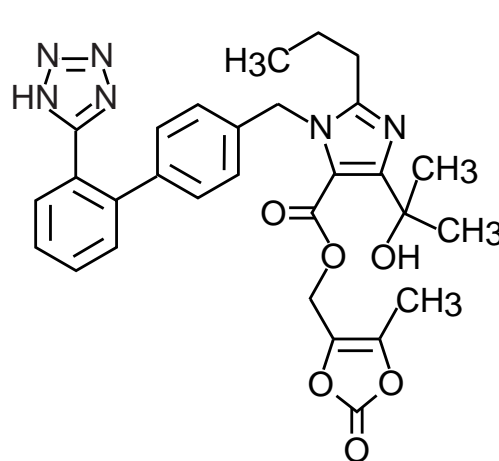

olmesartan medoxomil

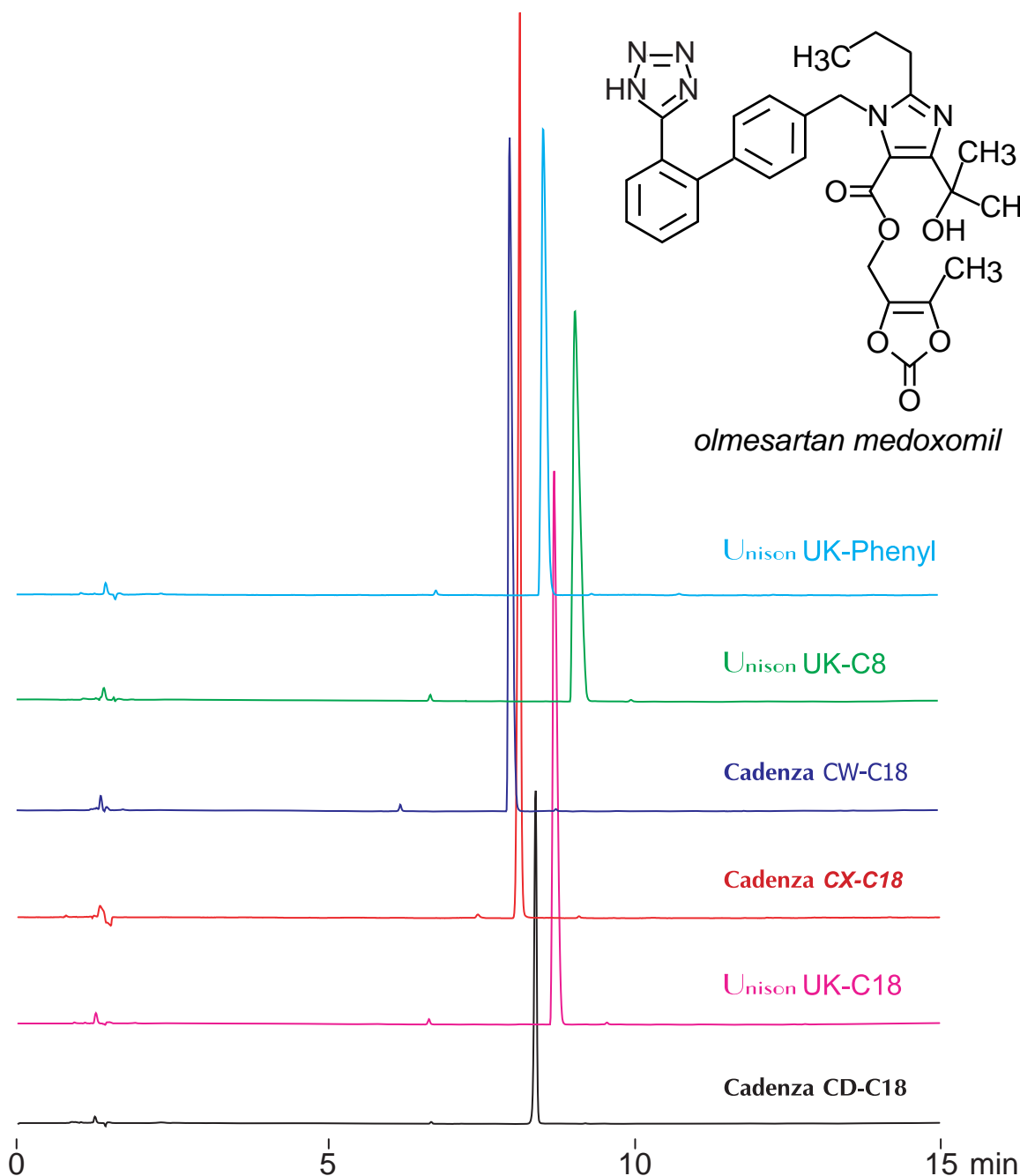

100 x 4.6 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

10-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 260nm

10µL (0.2mg/mL, acetonitrile / water / formic acid = 50 / 50 / 1)