

Unison UK-Phenyl Cadenza CW-C18



Unison UK-C8 Cadenza CX-C18

Unison UK-C18 Cadenza CD-C18 100 x 4.6 mm

HPLC COLUMN APPLICATION

Pemafibrate (Antihyperlipidemic)

ペマフィブラート (脂質異常治療薬)

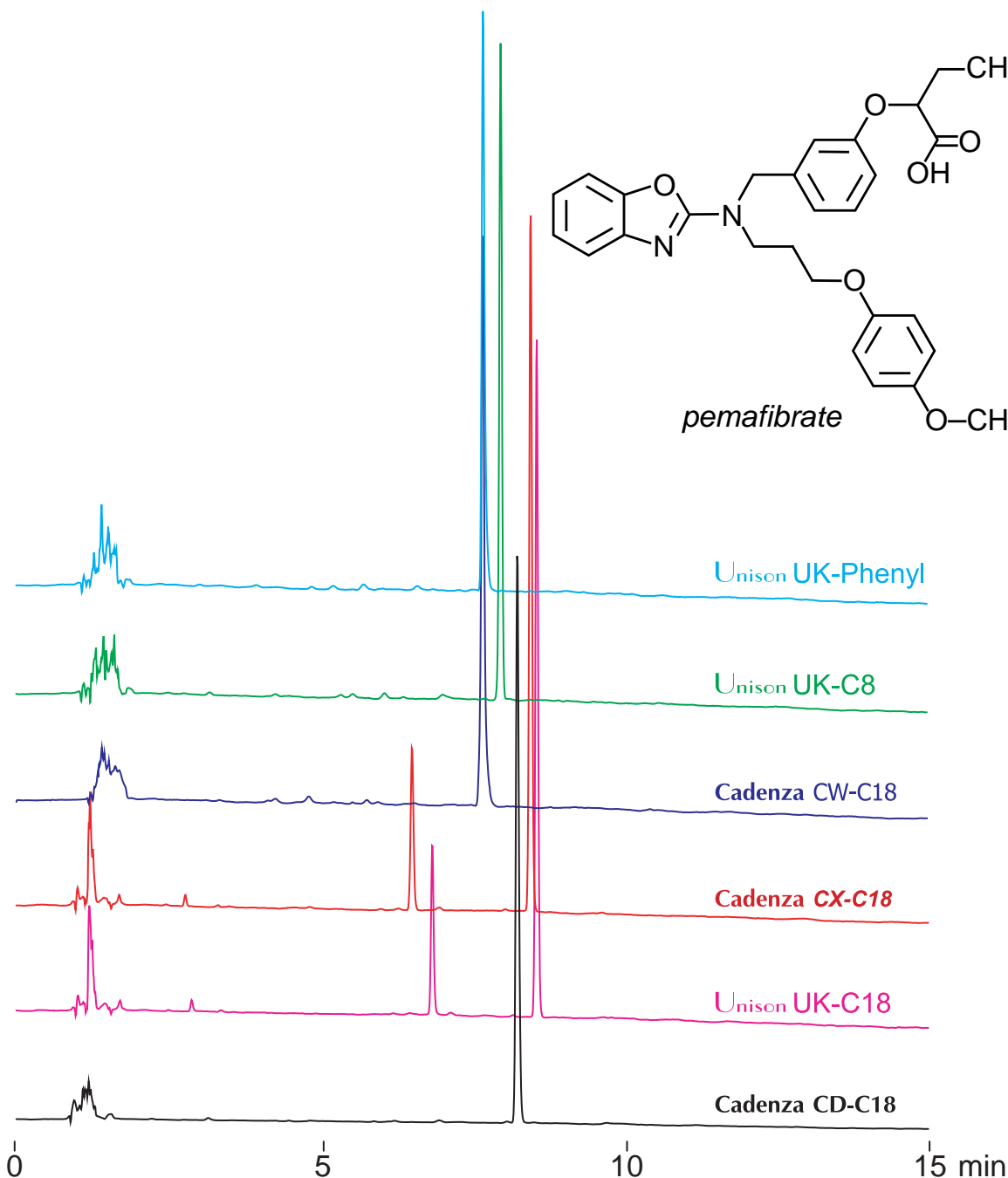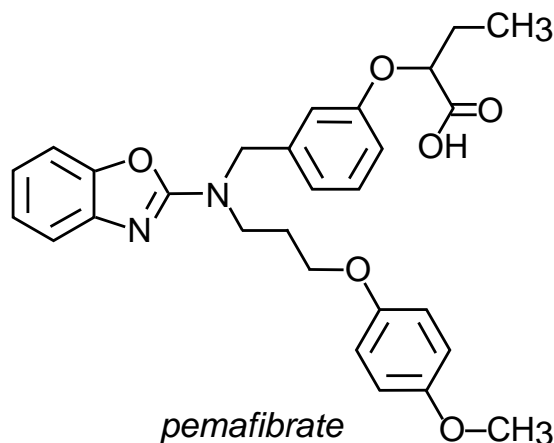

100 x 4.6 mm

A: 100mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 90 / 10

25-100%B (0-15min)

1.0mL/min (9MPa), 37°C, 280nm

10µL (1.0mg/mL, dimethyl sulfoxide)