



**Cadenza CD-C18**

100 x 3 mm

**HPLC COLUMN APPLICATION**

**LC-MS: Pemafibrate ( Antihyperlipidemic )**

ペマフィブラート ( 脂質異常治療薬 ) の LC-MS 分析

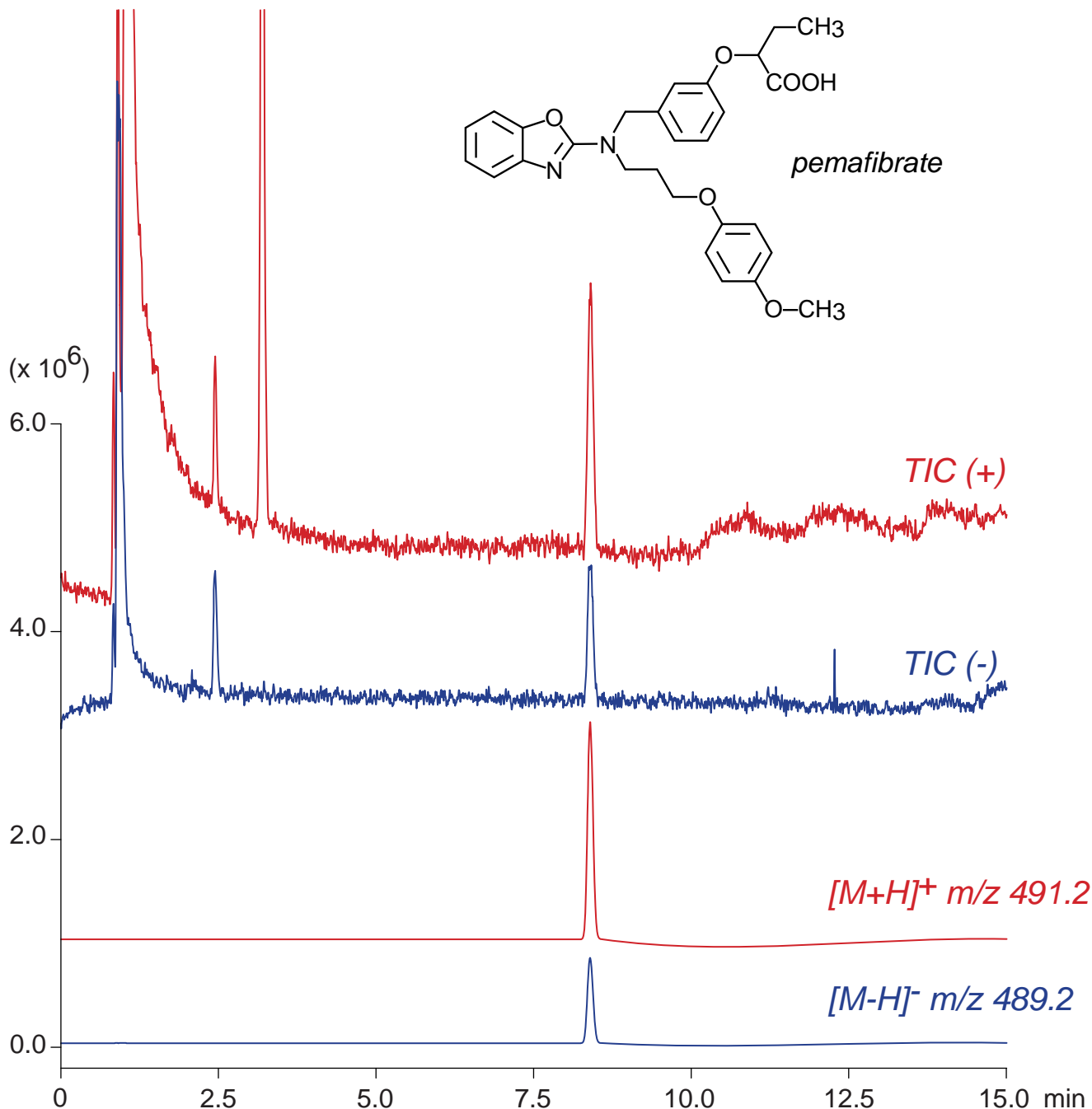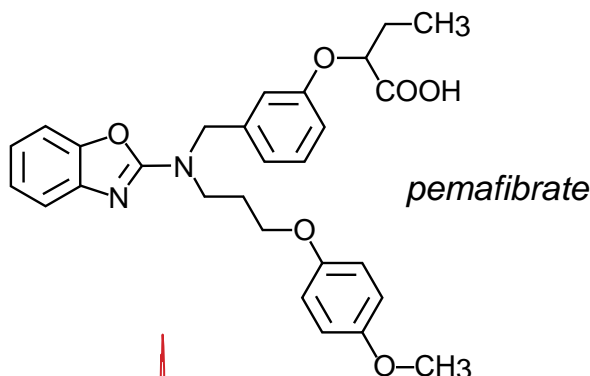

Cadenza CD-C18, 100 x 3 mm

A: water / formic acid = 100 / 0.1

B: acetonitrile / formic acid = 100 / 0.1

45-100%B (0-15min)

0.4mL/min (7MPa), 37°C, Single Quad. MS (ESI)

2 $\mu$ L (20 $\mu$ g/mL, dimethyl sulfoxide)