





Scherzo SM-C18

100 x 3 mm

HPLC COLUMN APPLICATION

LC-MS: Anticouglants

抗凝固薬の LC-MS 分析

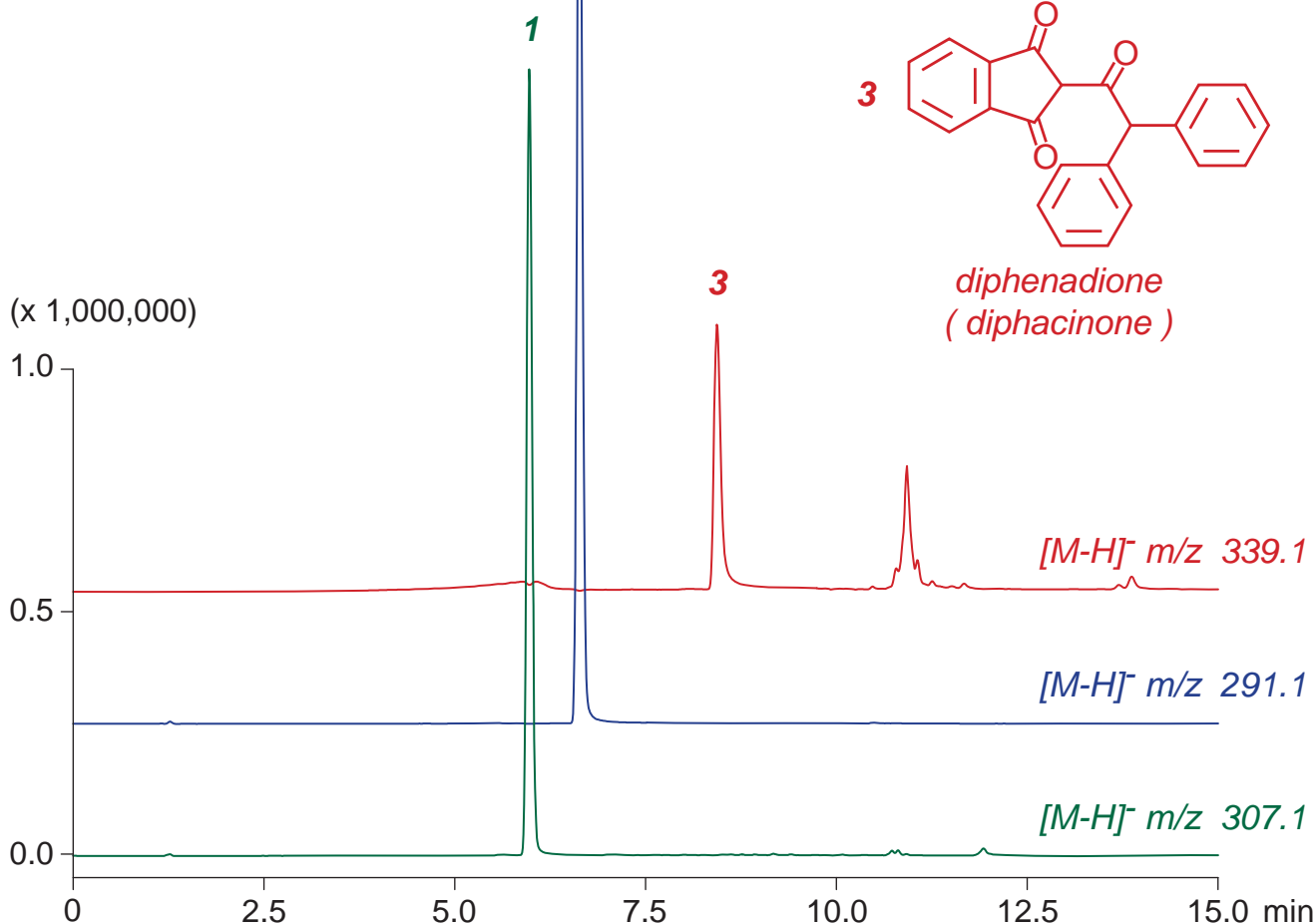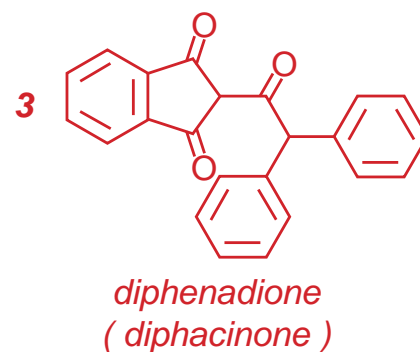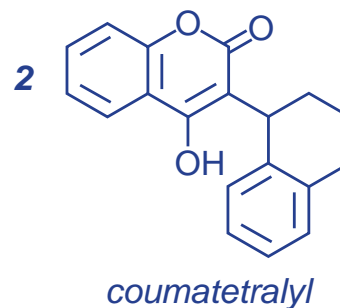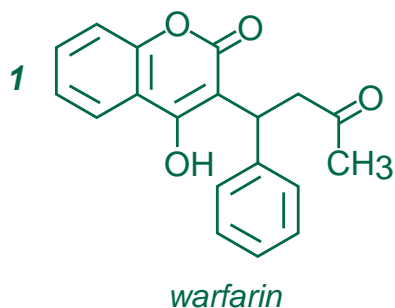

Scherzo SM-C18, 100 x 3 mm

A: methanol / 100mM ammonium acetate = 10 / 90

B: methanol / 100mM ammonium acetate = 95 / 5

10-100%B (0-9min), 100%B (9-15min)

0.4mL/min (8MPa), 60°C, Single Quad. MS (ESI)

5µL (2.5µg/mL, 100mM AcONH₄ / ACN = 50 / 50)