


Scherzo SW-C18  
 Scherzo SM-C18  
 Scherzo SS-C18

100 x 4.6 mm

HPLC COLUMN APPLICATION

## Edoxaban ( Anticouglants )

エドキサバン ( 抗凝固薬 )

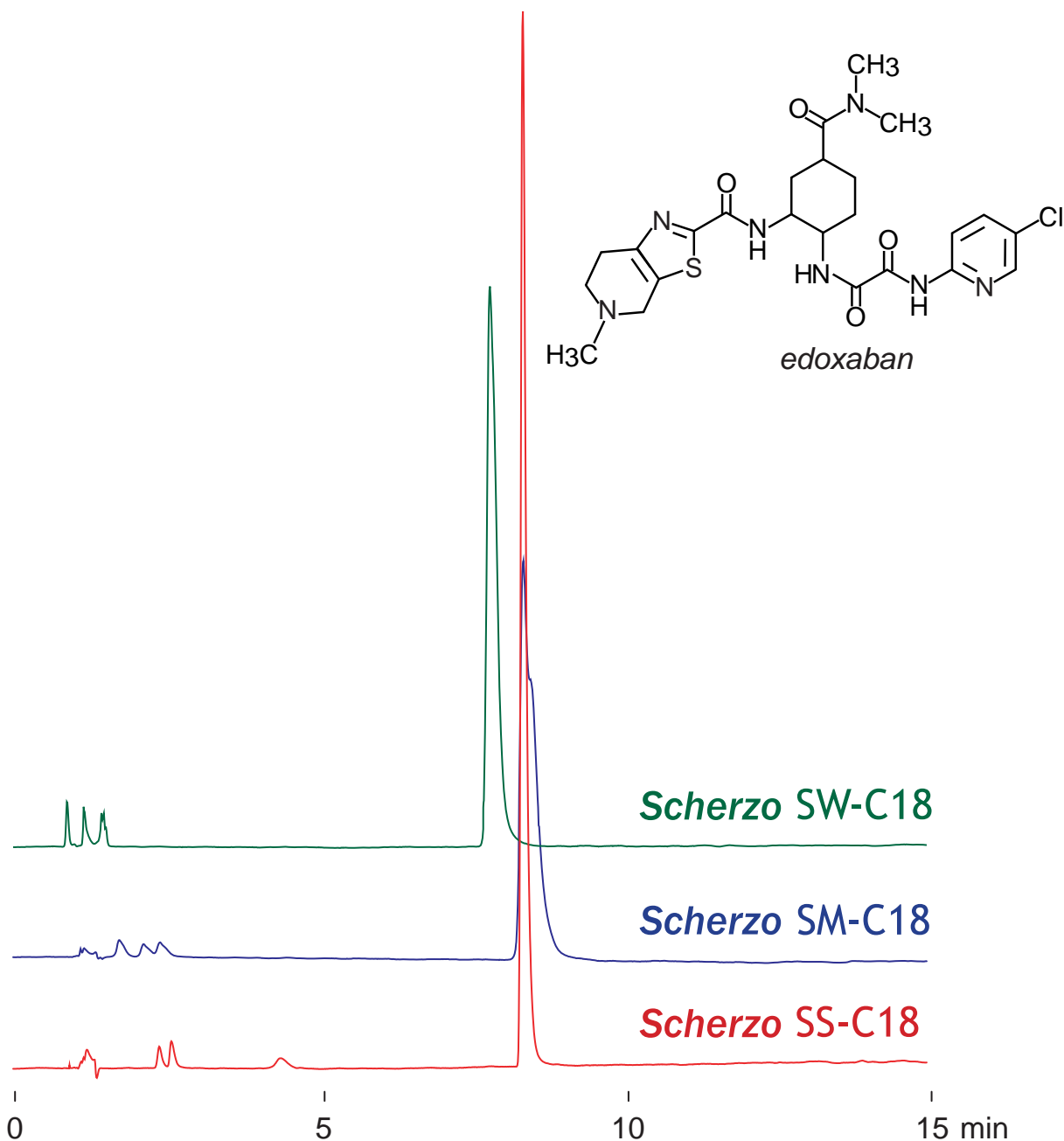

100 x 4.6 mm

A: 20mM ammonium formate

B: acetonitrile / 100mM ammonium formate = 70 / 30

30-100%B (0-12min), 100%B (12-15min)

1.0mL/min (11MPa), 37°C, 260nm

10µL (0.25mg/mL, acetonitrile / 100mM ammonium formate = 50 / 50)